# Sampl'Air Lite

bioMérieux SA  
69280 MARCY L'ÉTOILE / FRANCE  
RCS LYON 673 620 399  
Tel: + 33 (0)4 78 87 20 00  
Fax: + 33 (0)4 78 87 20 90

[www.biomerieux.com](http://www.biomerieux.com)

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---

Notice originale.

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## Document history

---

Translation of the original document.

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# TABLE OF CONTENTS

<b>1. PREFACE</b> .....	<b>35</b>
<b>2. PRESENTATION</b> .....	<b>36</b>
<b>3. TECHNICAL DESCRIPTION</b> .....	<b>37</b>
3.1 Front view .....	37
3.2 Rear view .....	37
3.3 User Interface.....	38
3.4 Remote control.....	40
<b>4. SAFETY</b> .....	<b>41</b>
<b>5. USE</b> .....	<b>41</b>
5.1 SWITCH ON /OFF .....	41
5.1.1 Switching on.....	41
5.1.2 Switching off .....	41
5.2 Battery loading.....	41
<b>6. SETTINGS</b> .....	<b>42</b>
6.1 Installation .....	42
6.1.1 Choice of the language .....	42
6.1.2 Setting the backlight.....	42
6.2 Setting the delay .....	43
<b>7 ACHIEVING A QUALITY SAMPLING</b> .....	<b>43</b>
7.1 Sampling programs .....	43
7.1.1 Creating a sampling program .....	44
7.1.2 Deleting a sampling program .....	45
7.1.3 Starting a sampling program.....	45
7.2 Positioning the Petri dish and starting the sampling .....	46
7.3 End of collection.....	47
7.4 Achieving quality sampling .....	47
7.5 Collection Method .....	47
7.5.1 Precautions for a good collection .....	47
7.5.2 Collection procedures.....	47
7.5.3 Cleaning procedure .....	48
7.6 Counting and calculation.....	48
7.7 Elaborating a sampling plan and interpreting results .....	49
7.7.1 Mapping of the location and determination of risks levels .....	49
7.7.2 Elaboration of a control plan .....	50
<b>8. CARE AND MAINTENANCE</b> .....	<b>51</b>
8.1 Cleaning.....	51
8.2 Calibration of the tachometer .....	51
8.3 Battery charging.....	51
8.4 Autotest .....	52
8.5 Troubleshooting .....	53
8.6 Reset procedure .....	54
<b>9. OPTIONS</b> .....	<b>55</b>
9.1 Tripod – ref : AESOC0018 .....	55
9.2 Telescopic pole – ref : AESSA0610 .....	55
9.3 Carrying case for Petri dishes– ref : AESSA0609 .....	55
<b>10. TECHNICAL FEATURES</b> .....	<b>56</b>
10.1 Correspondence table : sampling time/volume sampled .....	56
10.2 Environmental conditions for good operational practice :.....	56
<b>11. WARRANTY CONDITIONS</b> .....	<b>57</b>
<b>12. SPARE PARTS LIST</b> .....	<b>57</b>

**Check that the device has not been damaged during transport.  
Ensure that all of the accessories are included.**

#### **RECYCLING OF NiMH BATTERIES:**

In certain countries, it is forbidden to dispose of NiMH batteries with domestic or office waste. It is necessary to get proper information about local regulations. Do not use NiMH batteries that are damaged or are leaking.

#### **WARNING**

To prevent fire or shock hazard, do not expose this device to rain or humidity.

To avoid electrical shock, do not open the cabinet of this device. Servicing should be carried out by qualified personal only.

#### **FCC Compliance:**

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generate, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

#### **Caution:**

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate this device.

#### **Industry Canada (IC) Compliance**

This Class B digital apparatus complies with Canadian ICES-003.

**AES Blue line® is a registered trademark of BIOMÉRIEUX S.A.**

#### **Information Contained in this Manual**

Due to the informative nature of the information and diagrams contained in this manual, bioMérieux reserves the rights to make changes to the document without prior notice

**Users are required to read all the accompanying documents, including the statutory information, before using the equipment.**

#### **Meaning of the pictograms used in this manual:**

An alarm tone sounds

**Informative message**

**Important message**

## 1. PREFACE

Thank you for buying the Sampl'air Lite air sampler.

This one was designed to offer you a high level of performance and to be user friendly.

Your Sampl'air Lite has been developed to meet all your needs in microbial air analysis:

- The quality and efficiency of analysis are guaranteed by a **Quality Control Certificate** delivered with each machine.

- The sampling procedure based on impaction is recommended by the **ISO NF EN ISO 14698-1 standard**.

In addition to the technical potential of the instrument, the design features aim at providing ultimate usability for the customer.

A **light** and **autonomous** instrument combined with a **user-friendly programming panel** is a strong asset for routine analysis. **Sturdy**, the Sampl'air Lite is a bioimpactor that will be used for years.

## 2. PRESENTATION


Many equipment and methodologies (sedimentation) are suggested for the measurement of air contamination. Quality sampling requires performing equipment ensuring reproducible operation. Sampl'air Lite meets this requirement.

Manufactured by BioMérieux S.A, the Sampl'air Lite is a bioimpactor whose aim is to measure the air contamination level, using Petri dishes culture.

### **The impaction method:**

Collection of micro organisms is by air aspiration, through a grid which is situated a few millimetres above a Petri dish. The viable organisms in the air are impacted on the agar.

The following diagram describes this technique:

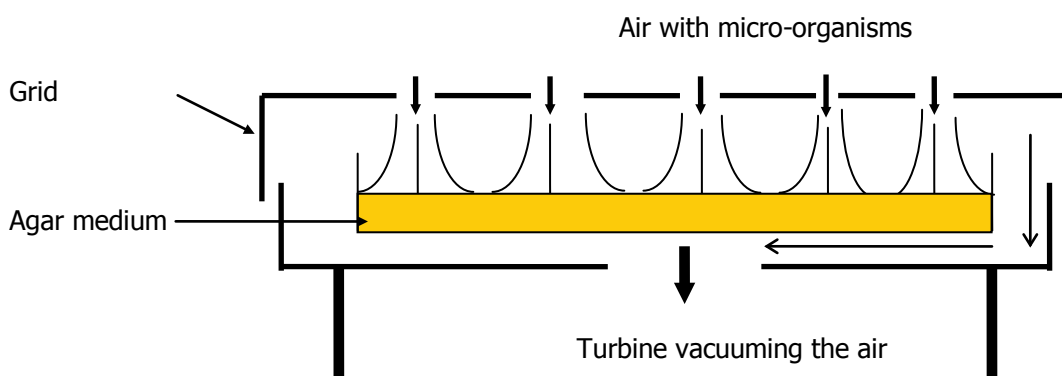The flow rate is known (100L/min). The impacted volume is proportional to the impaction duration. The number of microorganisms in the sampled air is easy to estimate by counting the CFU (Colony Forming Unit) on the Petri dish. The most probable number is assumed, using the table of correspondence (see appendix Number 2, in this manual).

The impaction technique is the most commonly used because it is reliable, easy to execute, economic (standards Petri dish) and safe (sterilisation with a disinfectant and/or autoclave of the Sampling head).


Handy and user-friendly

Monobloc and portable the Sampl'air Lite is handy and easy to use. An infrared control supplied with the device allows a remote controlled sampling.


The user interface consists of an intuitive numerical keyboard and a LCD screen.

### 3. TECHNICAL DESCRIPTION

#### 3.1 Front view

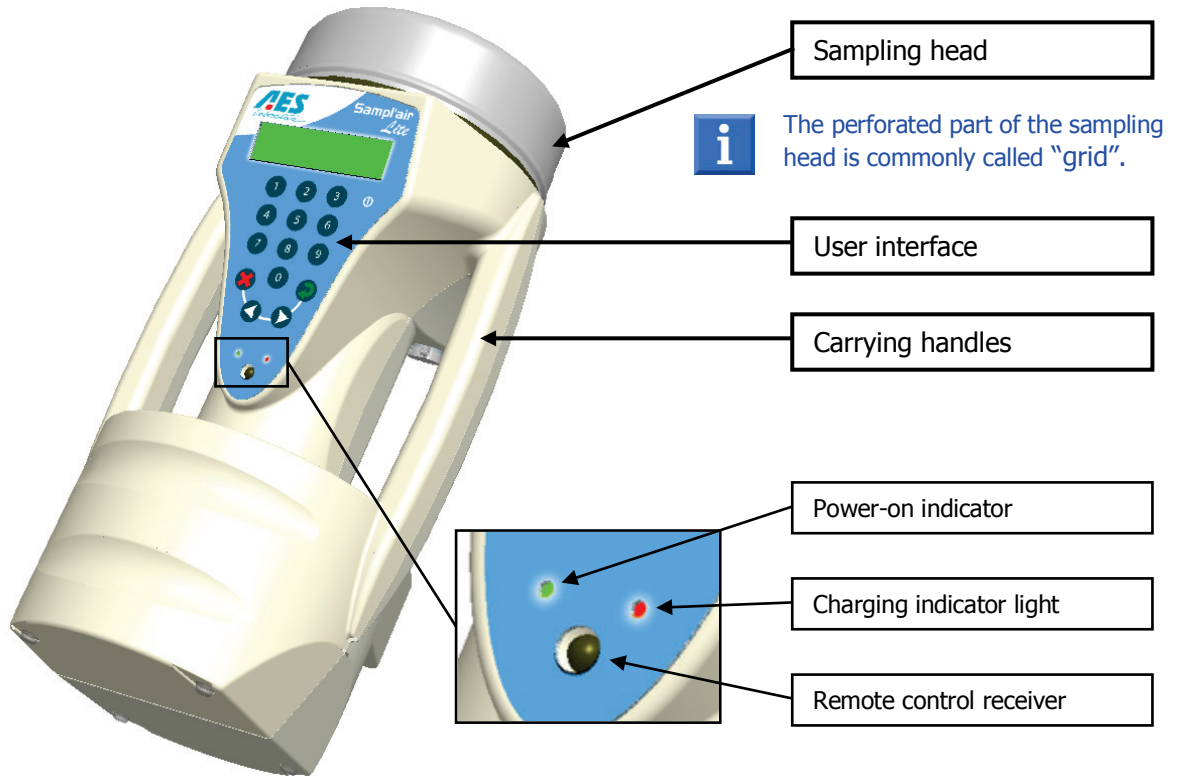
#### 3.2 Rear view
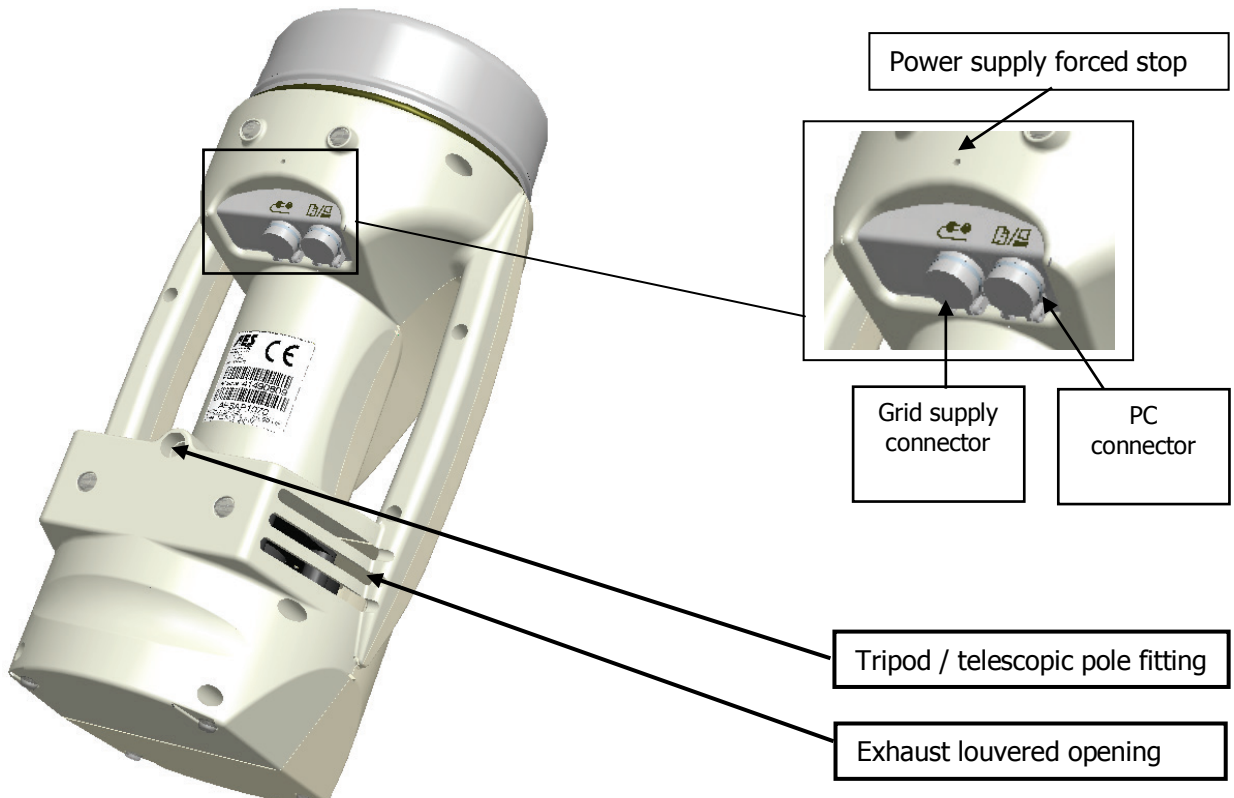
### 3.3 User Interface

The lower part of the screen displays the menus.

This is a numerical keypad.  
To dial figures, you simply need to press the appropriate key.

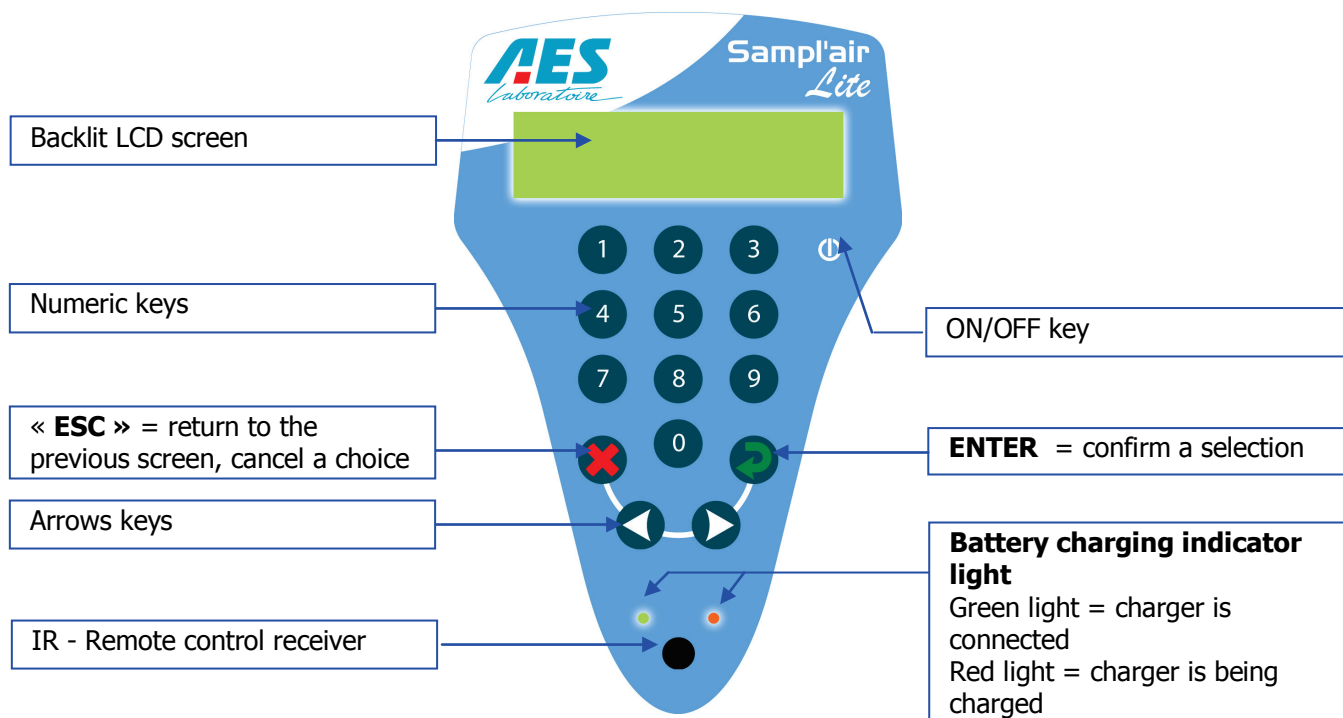The 'Enter' key is used to confirm a selection or a numerical value.

The 'Escape' key is used to move backward, to cancel a selection or an action.


These arrows are used to select an option in a menu (always displayed on the lower line of the screen).  
It may also be used to insert a blank space or to delete a figure.


**When a cursor ( \_ ) is displayed, you may enter data.**

A default value can be displayed.


**Sampl'air Lite synoptic**

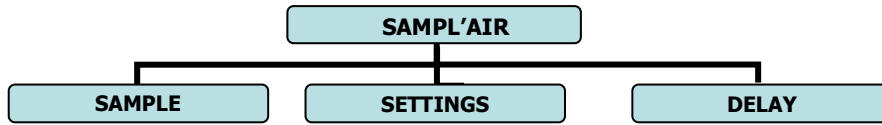**3 MAIN MENUS**

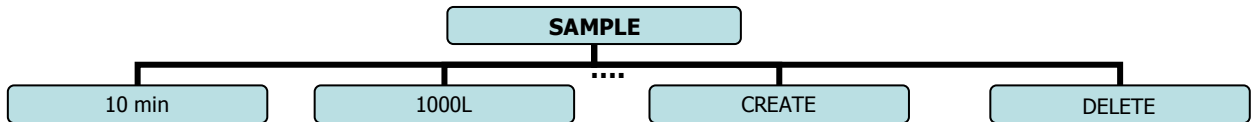→ In this menu, the user can select, create or delete a sampling program. (see § 7).

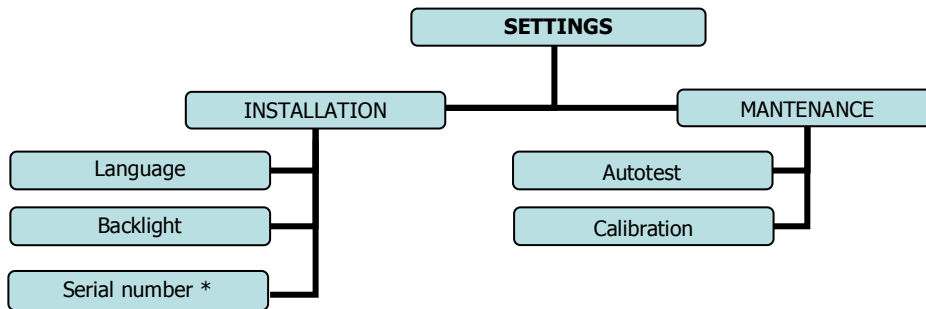→ SETTING menu enables the administrator to set the operating parameters for the Sampl'air Lite (§ 6), and the technicians of BioMérieux S.A to reach maintenance functionalities and advanced parameter setting.

\* Only BioMérieux S.A technicians may access this menu.

→ DELAY menu enables setting of a postponed start of the sampling cycle (§ 6.3).

### 3.4 Remote control

The remote control allows access to three functions accessible with 2 keys: « **Start** » and « **Stop** ».

Here are the functions:

- **Start a cycle:** short pressure the Start key (green arrow) once the program has been selected (the screen displays: "Start")
- **Interrupt a cycle:** during a cycle, short pressure the Stop key (red cross)
- **Switch the Sampl'air off:** long press the « Stop » key (red cross)

## 4. SAFETY

**RECYCLING NiMH BATTERIES:**  
 In certain countries, it is forbidden to dispose of NiMH batteries with domestic or office waste. It is necessary to get proper information about local regulations. Do not use NiMH batteries that are damaged or are leaking.

**WARNING**  
 To prevent fire or shock hazard, do not expose this device to rain or humidity.  
 To avoid electrical shock, do not open the cabinet of this device.  
 Only qualified personnel should carry out servicing


## 5. USE

### 5.1 SWITCH ON /OFF

#### 5.1.1 Switching on

Press and hold the **ON/OFF key**  for one second.

**Starting sequence:**
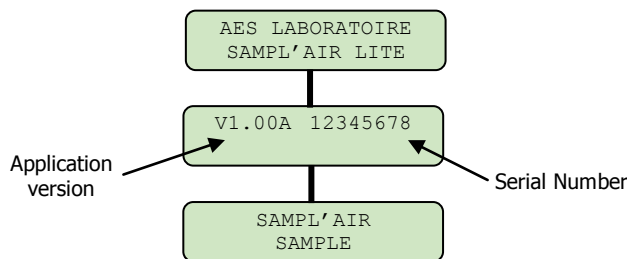
#### 5.1.2 Switching off

Press and hold the **ON/OFF key**  for at least one second, to switch the unit off, or press the "Stop" key on the remote control.

### 5.2 Battery loading

When receiving the Sampl'air it is advised to load the battery (refer to § 8.3 – battery charging).


## 6. SETTINGS

### 6.1 Installation

#### 6.1.1 Choice of the language

At first use, English is the default language. Language options are English, French, German and Spanish.

To change the language:

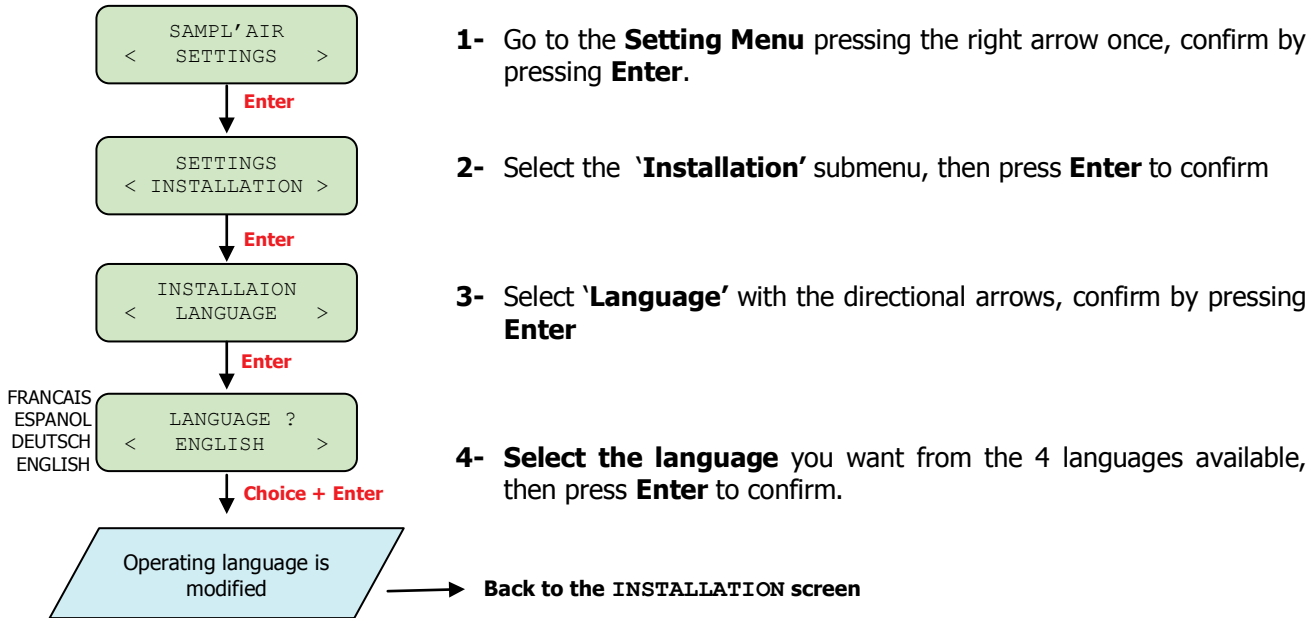
#### 6.1.2 Setting the backlight

The user interface backlight may be set in three ways:

- **Disabled**
- **Permanent** ( backlight always on when the Sampl'air is on)
- **Temporised** (After 20 seconds, the backlight turns off until the next key is pressed on)


To select the backlight:
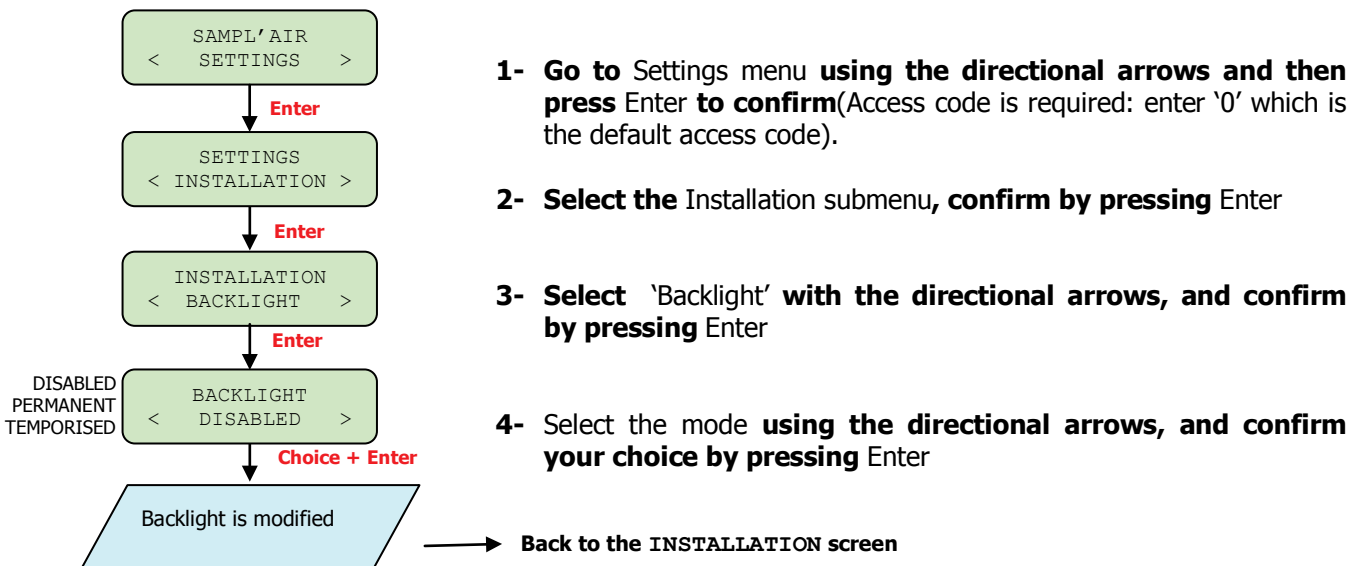
## 6.2 Setting the delay

When using the Sampl'air Lite for the first time, no is set (the delay is 00:00)

If you want a delay time before starting sampling:

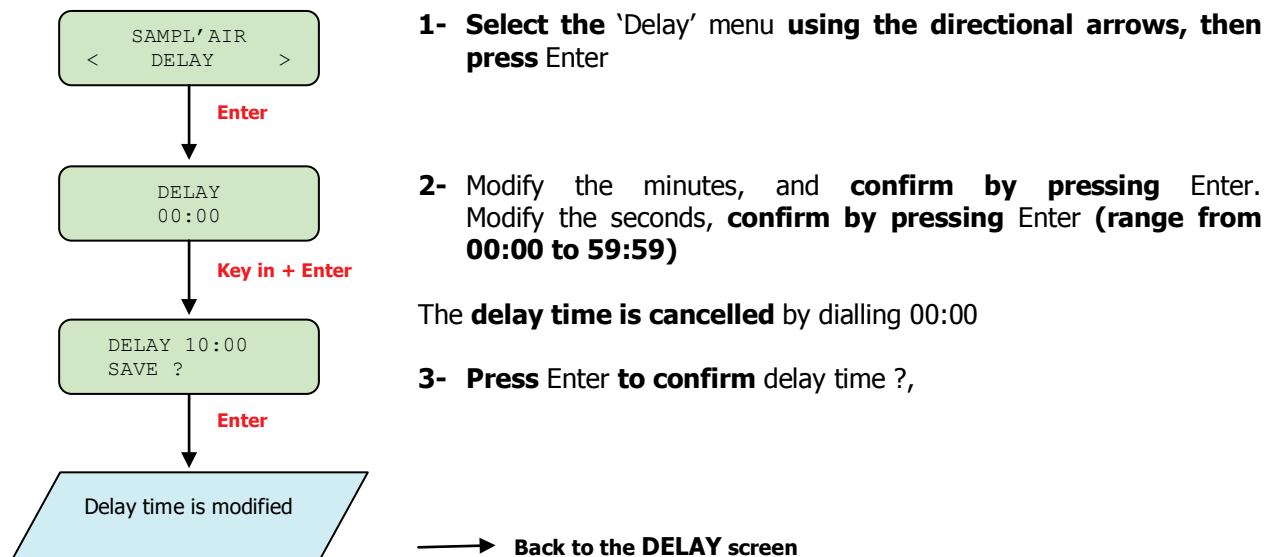The countdown timer allows you to enter a time value that will be displayed on the start confirmation screen. Count down starts as soon as you have pressed the Enter key. When the time is up, the collection starts.

## 7 ACHIEVING A QUALITY SAMPLING

A sampling cycle is started from the Sample menu.

### **SAMPLING FLOW RATE**

**The Sampl'air Lite flow rate is set for 100 litres / minute, +/- 5 minutes (in normal operating conditions). With 100l/min, the impaction velocity onto the agar is 16.8m/s.**

**The grid (258 holes diam. 0.7mm) has been developed in association with the best specialists of aerodynamic modelling and underwent cycles of technical and operational validation in real conditions.**

### 7.1 Sampling programs

When using the Sampl'air Lite for the first time, two programs are displayed

- **1000 L**
- **10 MIN**


These programs may be deleted and replaced by programs created by the operator (up to 10 sampling programs.)

To select the sampling program you want to start, use the directional arrows and confirm your choice by pressing **Enter**.

7.1.1 Creating a sampling program

Up to **10 programs** may be created. Sampling will be achieved according to a **time or volume entered by the operator.** (See table of correspondence § 10.1)

To create a program, proceed as follows:

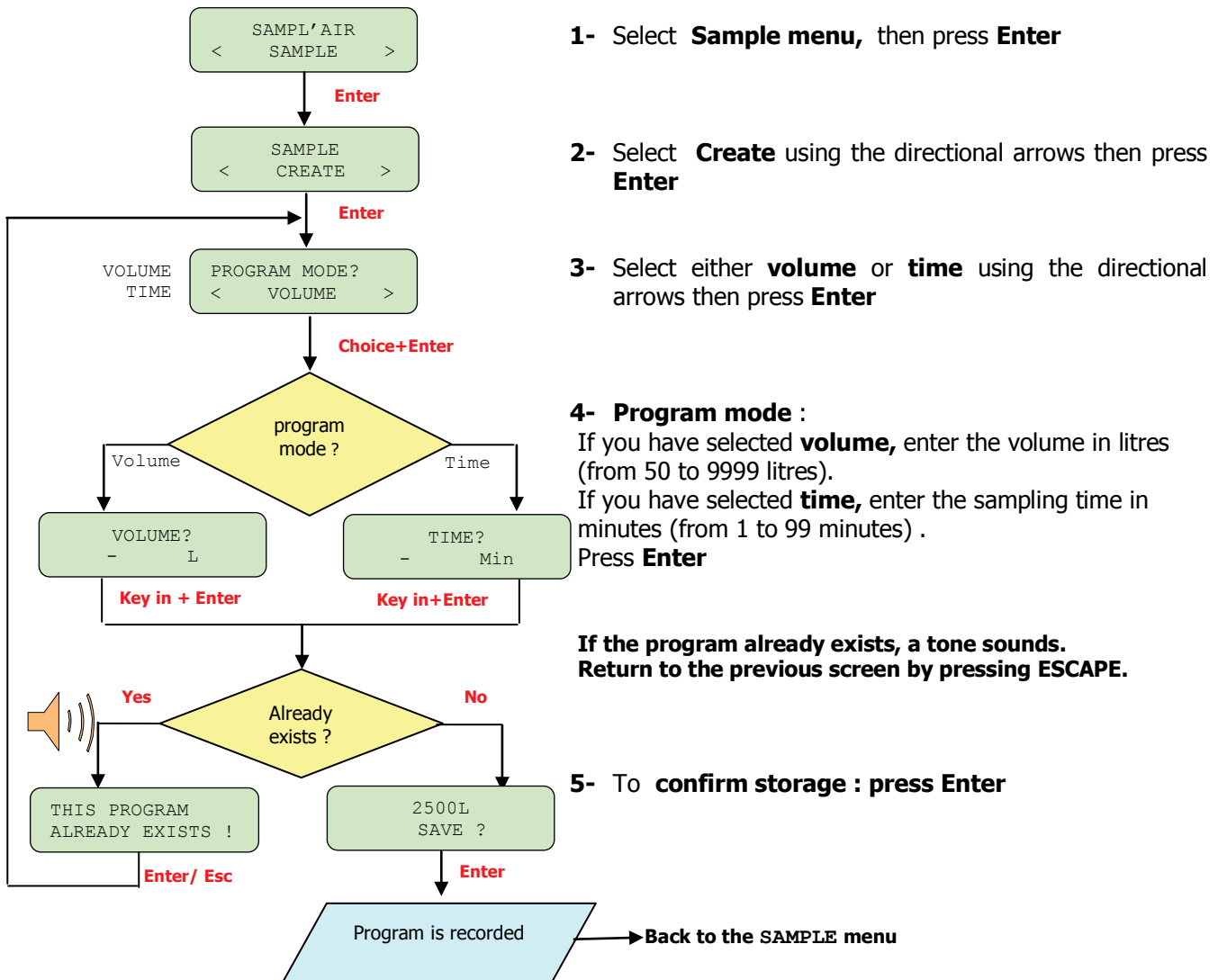1- Select **Sample menu**, then press **Enter**

2- Select **Create** using the directional arrows then press **Enter**

3- Select either **volume** or **time** using the directional arrows then press **Enter**

4- **Program mode :**  
 If you have selected **volume**, enter the volume in litres (from 50 to 9999 litres).  
 If you have selected **time**, enter the sampling time in minutes (from 1 to 99 minutes).  
 Press **Enter**

**If the program already exists, a tone sounds.**  
**Return to the previous screen by pressing ESCAPE.**

5- To **confirm storage** : press **Enter**


 **For the most usually used media, a cycle longer than 15 minutes may lead to the culture medium drying, and affect its fertility.**

### 7.1.2 Deleting a sampling program


You may delete programs if there is at least **one** program recorded.

To delete a program, proceed as follows:
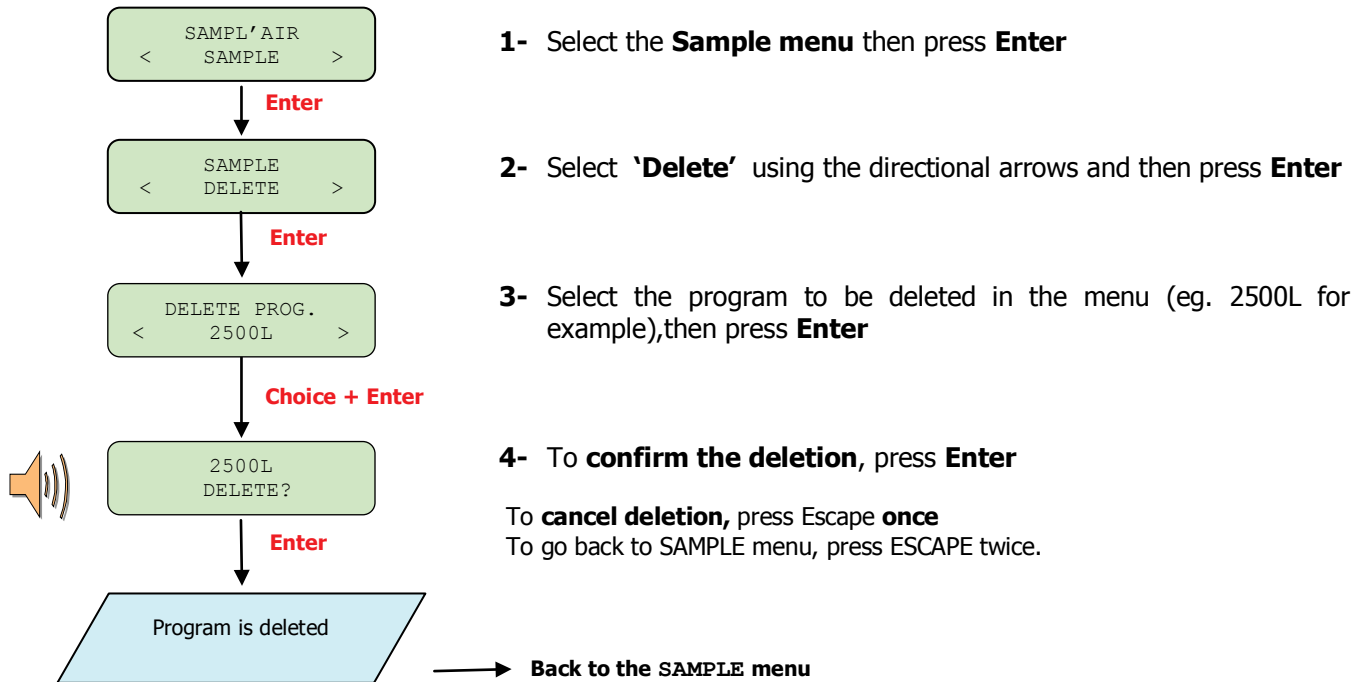
### 7.1.3 Starting a sampling program

To start a sampling cycle, proceed as follows:
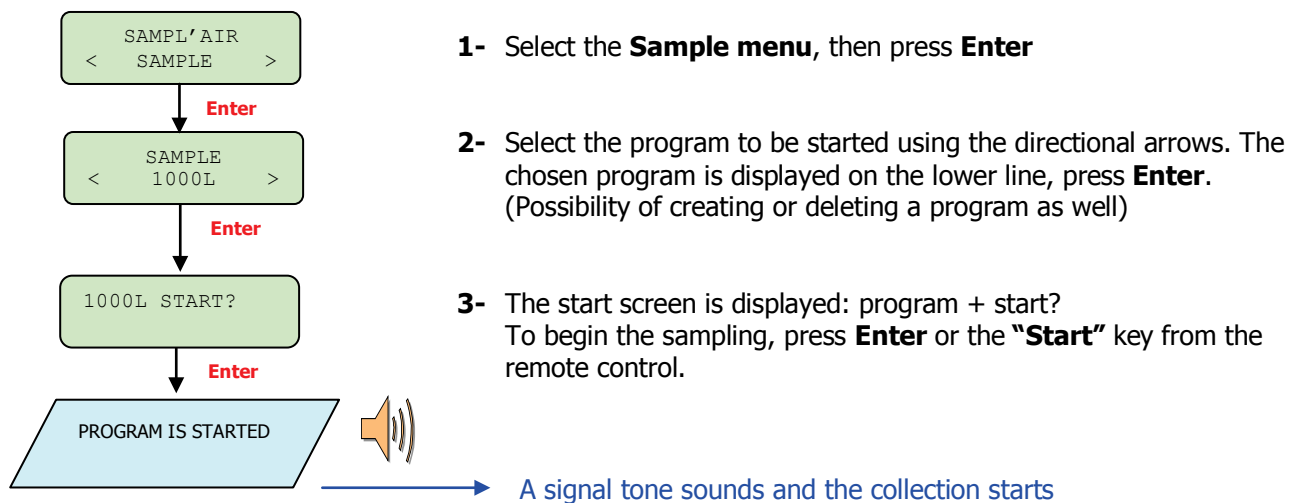
## 7.2 Positioning the Petri dish and starting the sampling

**Step 1:** Choose the program to be started

**Step 2:**

Place the sampling head vertically. Pull up the stainless steel grid to remove it from the unit; place the grid on a clean surface.

The sampling head's grid must not be touched.

**Step 3:**

Fit the Petri dish on the clips and make sure that the Petri dish is placed horizontally. The clips help with adjusting the Petri dishes in the center above the turbine. Make sure that the dish is properly inserted.

**Step 4:**


Remove the Petri dish's lid and put it on a clean surface. Place the Sampl'air Lite grid back. The Sampler is ready to start a cycle.

**Step 5:**

Start the collection by pressing the **Enter** key or using the remote control\*.

**\*Important:** the field of the remote control might be decreased because of nearby neon lighting.

The following screen appears :

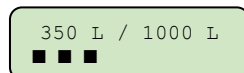If a delay has been set, after the start signal, the time is counted down and displayed before starting the cycle.

You may stop a collection while in progress by pressing the **ESCAPE** key or pressing the **STOP** key of the remote control.

### 7.3 End of collection

At the end of the sampling period a signal tone sounds and the message "End of sampling" is displayed.

Press Enter to stop the signal and go back to the SAMPLE menu.

Remove the impaction grid and close the Petri dish. Place the Petri dish into an incubator.

In case of problems occurring during the cycle, the error codes are successively displayed at the end of the cycle. Press ENTER to acknowledge each Error message. Some errors may stop. In this case, the message "Sample interrupted" is displayed.

### 7.4 Achieving quality sampling

To get good quality sampling a standard methodology is required. Indeed large ranges of factors have an impact on the harvest of micro organisms and their growth (air movements, atmospheric pressure, impaction speed, physiological state of the micro organisms).

#### **Quality medium sampling**

The users are free to choose the media which are the most relevant to grow the organisms they wish to isolate. In any case, they must ensure that the following precautions are taken:

- The depth of the agar in the Petri dish should be reproducible and measure from 2,5 to 3,5 mm
- The agar surface should be flat and levelled.
- The storage methods should be appropriate to the medium
- The agar should not present any dryness

**If the instructions above are not followed, the measurement parameters might be difficult to reproduce and the results might be false. Besides, some media, which have been developed for clinical bacteria, are not recommended for environment microorganisms research.**

- The agar should not present water droplets when used. Any excessive humidity on the agar may affect results

Some adapted media, suggested by BioMérieux S.A suit perfectly this kind of analysis.

### 7.5 Collection Method

Different procedures need to be followed in order to prevent fluctuations due to exterior phenomena or users. Those procedures should be strictly followed when collecting in Zones 3 and 4.(see § 7.7.1).

#### 7.5.1 Precautions for a good collection

- Hands and nails should be washed before using the Sampl'air Lite (using products such as PROSEPTOL® is recommended)
- Protection clothing should be worn in accordance with the collection environment.
- Avoid breathing and coughing in direction of the Petri dish when open.
- Complete the work quickly and efficiently and avoid any delays.

#### 7.5.2 Collection procedures

- Collect from low contaminated areas moving up to higher contaminated areas.
- It is recommended to achieve sampling in the middle of the room or in different points especially when the size of the room is significant.
- Duplicate sampling (at the same place) in order to eliminate false results and to allow a rate of contamination calculation.

### 7.5.3 Cleaning procedure

After cleaning the device (See "Care and maintenance" chapter), vacuum or dry the device with warm air in order to dry the sampling head.

**Follow carefully the operation instructions and sampling plan described by the team that analyses the results. File the test report carefully after every analysis and identify the Petri dish.**

### 7.6 Counting and calculation

After incubation, count the number "n" of CFU (Colony Forming Units) which have grown.

In order to take into account the probability of having coinciding impacts of 2 or more micro organisms (implying only one countable colony), an algorithm is applied according to Peto and Powel work <sup>(1)</sup> :

$$E_n = \left\{ \frac{1}{T} + \frac{1}{T-1} + \frac{1}{T-2} + \dots + \frac{1}{T-n+1} \right\}$$

T = number of holes in the sampling head

n = number of numerated colonies

$E_n$  is the estimated number of viable micro organisms giving n positive impacts, when using a grid whose number of holes = T.

This formula suggests that the micro organisms flow stops when a micro organism reaches the hole n.

This flow being irregular, the estimated number of particles, when n impacts are observed, has to be equal or superior to  $E_n$ , but inferior to  $E_{(n+1)}$ .

The adjusted number N is obtained using:

$$N = \frac{1}{2} \{ E_n + (E_{n+1} - 1) \}$$

Consult the table in the appendices section (appendix number 2). A table indicates the n value (number of colonies actually numerated) corresponding to N (probable number of viable impacted micro organisms - either estimated or adjusted).

Once N is identified, the controlled volume ( $V = \text{sampling duration} \times 0.1 \text{ m}^3 / \text{min}$ ) and the airborne contamination level can be obtained using  $\frac{N}{V}$  Colony Forming Units /  $\text{m}^3$ .

#### **Example**

After sampling and incubation, 120 colonies are counted on Petri dish. Then  $n=120$ .

$$N = 161 \quad V = 5 \text{ min} \times 0.1 \text{ m}^3 / \text{min} = 0.5 \text{ m}^3$$

$$\text{Level of airborne contamination: } \frac{161}{0.5} = 322 \text{ CFU} / \text{m}^3.$$

(1) **Peto S. and Powel E. O. (1970). The assessment of aerosol concentration by means of the Andersen sampler. J. App. Bacteriol. 33,582-598.**

## 7.7 Elaborating a sampling plan and interpreting results

The sampling plan establishment and the aerobiocontamination evaluation are part of a HACCP plan (Hazard Analysis Critical Point). It is a voluntary, autonomous process based on prevention. The user has to evaluate the contamination level to be met in order to insure that the critical point is reached (place, product process) and controlled.

The aerobiocontamination management requires a multidisciplinary team whose aim is to define the areas at risk, the preventive actions, the sampling points, their frequency, the control specifications (target, alert and action levels), and the corrective measures to be taken.


Prior to results interpretation, control specifications must be defined. The control specifications are based on existing referentials (GMP, APSEC, usual Standards...), or referentials created thanks to the inherent data. The international ISO/14698-1 standard about air quality, recommends a recording of contamination levels, at least for several weeks, in order to create a reference and assess target, alert, and action levels.

### 7.7.1 Determination of control specifications

#### Referential for clinical environment:

According to the ASPEC recommendations (Association for the prevention and study of contamination), it is possible to determine 4 areas at risk. The risk levels of these areas are linked to the patients’ fragility, and to the risk due to the type and duration of the medical act.

Risk level of a medical act in clinical environment:

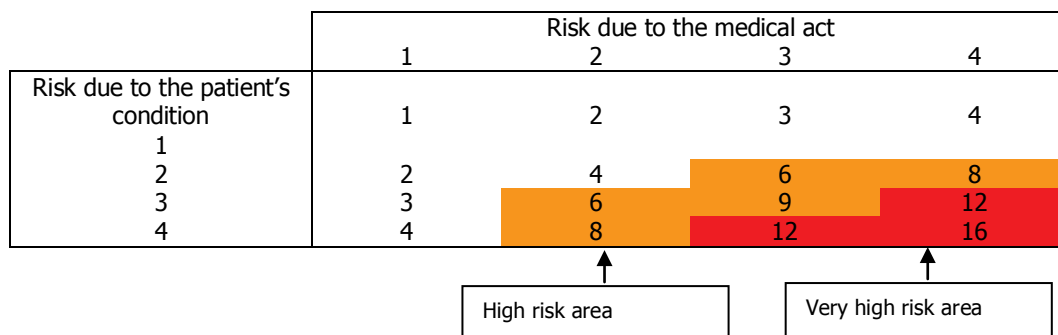Control specifications for aerobiocontamination in very high risk and high-risk areas:

	High risk area		Very high-risk area	
	Micro-organisms CFU*/ m <sup>3</sup>	Molds CFU*/ m <sup>3</sup>	Micro-organisms CFU*/ m <sup>3</sup>	Molds CFU*/ m <sup>3</sup>
Target level	500	1	10	1
Alert level	100	1	5	1
Action level	10	<1	<1	<1

#### Referential for pharmaceutical environment:

The table below gives a summary of the aerobiocontamination rates according to clean rooms class, as defined in the "Good Manufacturing Practice" in Pharmaceutics:

GMP Classes / Categories	A	B	C	D
Contamination (recommended limits) UFC* / m <sup>3</sup>	<1	10	100	200

\* CFU: Colony Forming Unit

The sampling places should be determined according to critical points identified in a precise area. With a mapping of the area to be controlled, critical points may be identified.

7.7.2 Determination of target, alert, and action levels

Once the critical points are defined in each area to be controlled, target, alert, and actions levels must be determined. Appropriate corrective actions must be taken as soon as the action level is reached.

The method for the target level determination should be defined in a quality insurance procedure, and be mentioned in the Quality Process Certificate.

*Example for the identification of risk levels:*

Reference area	Target level		Alert level		Action level	
	Total Flora	Yeast and molds	Total Flora	Yeast and molds	Total Flora	Yeast and molds
Area 1						
Area 2						
Area 3						

It is then the aim of the person in charge through his experience and knowledge, to assess the acceptable level of contamination in various areas. He/ she will also have to determine the corrective actions (removing filters, new waste management, and new organisation of services ...) to reduce the impact of the contamination levels, and improve the environment management.

Within a quality insurance process, preventive monitoring and maintenance must prevail against curative actions.

7.7.3 Control plan establishment

Once the reference plan is determined, a control plan should be described in order to validate the hygiene level. The users depending on their aims can develop different sampling plans.

Sampling plan may be:

- ✓ A routine sampling plan for a global prevention.
- ✓ A specific sampling (Areas at risk...)
- ✓ A sampling to measure the incidence of an unforeseen phenomenon (Impact of rebuilding work, removing of filtration filters...)
- ✓ or any other process if needed.

Whatever the aims are, protocol, sampling locations, and duration of sampling should be clearly defined to compare and analyse the results.

Finally a quality control sheet (such as the one given in appendices of this manual) should be filled in to gather any useful information.

It is important to specify on this sheet any parameter that might have an impact on the result (Level of personal activity, rebuilding works, air recycling system, sampling area ...)

Air contamination control can be obtained only if concurrently an effort for staff training is made.

## 8. CARE AND MAINTENANCE

### 8.1 Cleaning

The Sampl'air Lite should be cleaned after each sampling campaign because a wide range of stains can contaminate the sampling head and then reduce the reliability of the results.

The part that has to be cleaned the more carefully is the sampling grid, as it is between the air and the medium.

The unit surfaces are not deteriorated by usual cleaning products (Ethanol 90 °, bleach, «Dégragerme», «Proseptol»).

#### Between each sampling

External surfaces and the sampling grid should be washed with a disinfecting spray or a lotion between each sampling. Some cleaning and disinfectant wipes may be used (AESANI205228).

#### Between each sampling campaign

The stainless steel impaction grid may be autoclaved (at 134°C for 18 minutes or at 121°C for 15 minutes)

**ONLY THE STAINLESS STEEL GRID MAY BE AUTOCLAVED.**

This equipment is not to be used in water. Consequently it must not be swilled out.

### 8.2 Calibration of the tachometer

The tachometer defines the rotation speed of the turbine allowing the nominal airflow of 100 litres per min. To ensure a good quality of the samplings, a control and calibration of should be done **every year**. The calibration of the tachometer is **strictly reserved** for the qualified personnel only (BioMérieux S.A staff, or persons who have been trained by BioMérieux S.A staff).

### 8.3 Battery charging

- **Before using Sampl'air Lite for the first time, it is advised to charge it for at least 12 hours.**
- **If left unused for a certain amount of time, a fully charged battery will discharge over time. It is recommended to recharge the battery at least once a month in order to preserve its lifetime.**

- The charger's temperature may be higher during charging. This is normal, and does not indicate a problem.
- Charging may take longer if the battery is not new or the temperature increases during charging.
- A battery that is new or has not been used for a long time may get discharged more quickly. This is normal. A battery stays loaded longer once it has been used and then charged several times.
- If a battery that is not completely discharged is put in charge, it will discharge quickly when used afterwards. Its normal performances may be retrieved when it is totally discharged then charged again.
- A battery that stays charged only a short time, even after being totally charged, is faulty and needs to be replaced.
- Do not try to charge a fully charged battery; this may reduce its performances.
- The battery becomes performant more quickly if you charge it at room temperature, and when the low battery warning occurs.

**LOW BATTERY WARNING**

A 'low battery' indicator light appears on the screen (symbol  ) and indicates that your battery has enough autonomy to run up to 4 cycles of 10 minutes before switching off. A few minutes before switching off, a tone is sounded and on the message 'Low battery' is displayed.

Sampl'air Lite may be connected to the charger to prevent the cycle interruption.

**Charging procedure:**

Connect the battery charger to the device as shown below:

Connect the charger to the mains power supply (voltage 100 to 240 VAC+/- 10%, 50-60 Hz, 700mA). The charger supplies a voltage of 15 V cc, 2A.

**Only the charger delivered by BIOMÉRIEUX S.A must be used (FRIWO type MPP30 FW7362).**

The green light indicates that the charger is connected. When no charger is connected, no light switches on.

The red light switches on showing that the device is charging. The red light switches off when the battery is loaded.

 **The total charging time depends on the energy level of the battery. It is recommended to let the Sampl'air Lite loading at least 5 hours after a sampling campaign.**

**RECYCLING OF NiMH BATTERIES:**  
 In certain countries, it is forbidden to dispose of NiMH batteries with domestic or office waste. Please, be sure to be informed about any local regulation. Do not use NiMH batteries that are damaged or have leaked.

**8.4 Autotest**

Control and diagnosis function are present in the Sampl'air Lite. These functions are reserved for BioMérieux S.A qualified personnel only.

## 8.5 Troubleshooting

<b>Problem</b>	<b>Possible cause</b>	<b>Corrective action</b>
<b>No display</b>	1 <sup>st</sup> case : the battery is empty or the charger is not connected	Charge the Sampl'air Lite (§ 8.3)
	2 <sup>nd</sup> case: The battery is charged and the charger is unplugged. Electronic problem	Contact the after-sales service
<b>Battery charger is connected, the red warning light does not light up</b>	1 <sup>st</sup> case : the charger of the battery is not correctly connected	Check the connection of the charger
	2 <sup>nd</sup> case : the voltage supplied by the charger is not correct	Use a battery charger adapted: 15Vcc, 2A minimum
	3 <sup>rd</sup> case : electronic problem	Contact the after-sales service
<b>The charger is connected, the green charging led lights up, the red charging led does not light up.</b>	1 <sup>st</sup> case : the battery is already charged	the Sampl'air Lite can be used
	2 <sup>nd</sup> case : an error message is displayed	See below for the meaning of the error message
	3 <sup>rd</sup> case : the battery is not charged, any error message is displayed	Contact the after-sales service
<b>Screen display : EMPTY BATTERY</b>	Battery is empty	put the Sampl'air Lite on charge (§ 8.3)
<b>Screen display : LOW BATTERY</b>	Charge with low battery	At the end of the sampling cycle, put the Sampl'air Lite on charge (§ 8.3)
<b>Screen display : unsuitable Power supply</b>	The voltage supplied by the battery charger is <13V or >17V.	Disconnect immediately the battery charger, use a battery charger delivering 15Vcc, 2A mini. (§ 8.3)
<b>Le Sampl'air Lite is switched on, but some keys do not work.</b>	1 <sup>st</sup> case : Program is jammed	Press the "forced stop" button (§3.2)
	2 <sup>nd</sup> case : keyboard is out of order or electronic problem	Contact the After Sales Service
<b>Screen display : ERREUR 11</b>	Electronic problem : battery circuit breaking	Contact the After Sales Service
<b>Screen display : ERROR 62</b>	Turbine's speed variates abnormally	Make sure that the air circulates properly, that the dish is placed correctly. If everything seems OK, contact the After Sales Service (electronic problem).
<b>Screen display : ERROR 63</b>	Electronic problem : Tachymètre turbine	Contact the After Sales Service

## 8.6 Reset procedure

Reset is not an usual operation, but it may be necessary to do so after the update of the software or if the device is not operating correctly.

To reset your Sampl'air Lite, press simultaneously the **Enter and Escape key** during 5 seconds.

The following data are reset to default values:

- Program list (10 min / 10 L)
- Configuration options
  - language : english
  - backlight : off

**NOTE :** After resetting, it is necessary to resume the complete programming ( list of the programs, choice of the language, backlight and administrator code).

Are not reset:

- serial number of unit
- calibration data
- operational counters

In case of jamming of the device, stopping the sampler is also possible when inserting a stick in the « **forced stop** » hole that is situated at the bottom of the device. Note that the data you have keyed in and saved in the memory are not deleted when proceeding like this, it resets only the main PCB.

## 9. OPTIONS


### 9.1 Tripod – ref : AESOC0018

The tripod allows supporting of the device

### 9.2 Telescopic pole – ref : AESEA0610

The telescopic pole allows the Sampl'air Lite to be used in places hardly accessible.
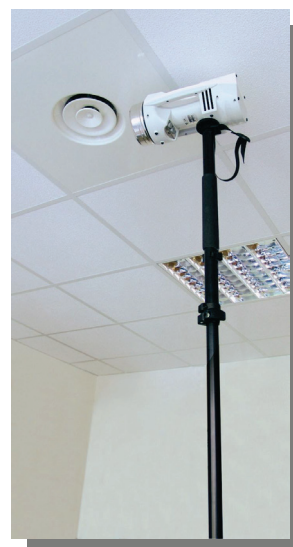
### 9.3 Carrying case for Petri dishes– ref : AESEA0609

Lightweight and washable case equipped with a strap. It includes a foam rubber where 40 Petri dishes may be placed (or 20 Petri dishes and 50 disinfectant wipes).

## 10. TECHNICAL FEATURES

<b>Dimensions :</b>	130 x 170 x 300 mm
<b>Weight :</b>	1,8 kg
<b>Noise :</b>	50 dBA at 1 meter
<b>Sampling head :</b>	Stainless steel AISI 316L, autoclavable, 258 holes diameter 0,7 mm
<b>Airflow :</b>	100 litres / mn, $\pm$ 5%
<b>Sampling mode:</b>	time (1 to 99 minutes) or volume (50 to 9999 litres)
<b>Battery :</b>	NiMH - 9,6V - 3,8 Ah, autonomy : 25 cycles of 10 min (Ambient temperature higher or equal to 20°C, backlight deactivated)

<b>Battery charger :</b>	FRIWO MPP30 FW7362/15 <b>Necessarily ordered from BIOMÉRIEUX S.A</b> Input : 100 – 240 V~ / 50 – 60 Hz / 700 mA Output : 15 V= / 2 A ; Low Voltage (SELV) Interchangeable mains plugs : - EN 60320 type, France / Europe double pole, outlet plugging - IEC double pole type, use with a cord
--------------------------	---

**NOTE: qualified personnel only should carry out replacing batteries**

### 10.1 Correspondence table : sampling time/volume sampled

TIME	VOLUME
30s	50 l
1 mn	100 l
2 mn	200 l
4 mn	400 l
8 mn	800 l
10 mn	1m <sup>3</sup>

Average speed : - 16,8 m/s for 100 l/min

### 10.2 Environmental conditions for good operational practice :

The running of Sampl'air Lite is ensured under the following conditions:

- When operating : Temperature within 15°C and 35°C
- When charging battery :Temperature within 15°C and 35°C
- Hygrometry within 30% and 85% of relative humidity non condensing
- Atmospheric pressure within 700 and 1040 HPa
- According to the type of lighting used (strong or indirect natural light, light bulbs, neon tubes), the action of the infrared remote control may be reduced.
- Connectors are protected from corrosive gas or liquid with folding caps.

**After unpacking or after storage at low temperature, it is strongly recommended waiting at least 2 hours before using the device (to let it acclimatize to room temperature)  
The company is not responsible for any damage due to failure to operate the device without following the instructions.**

## 11. WARRANTY CONDITIONS

The Sampl'air Lite is guaranteed against manufacturing or components problem (except batteries) for twelve months from the delivery date. BioMérieux S.A will repair broken parts free of charge in its premises during warranty period. Any instrument with a defective working problem should be sent back carriage paid by BIOMÉRIEUX S.A agent or customer when necessary. The warranty will not be valid for damages due to a transport problem, accident or internal interferences caused by a bad electricity network.

Any reparation is carried out by a technician from our After Sale service department or by a person who has received appropriate training.

A contract to perform periodic calibration of the Sampl'air Lite can be offered. For more information, contact us. The calibration is carried out in our factory.

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