



Peak of Perfection

Made for a faster and easier daily routine:
Eppendorf Liquid Handling Consumables

Our Benchmark: Your Requirements

How liquid handling consumables can influence your assay results

The reliability and consistency of your research results are our priority when developing consumables.

Chemicals like slip agents, plasticizers or biocides, used as manufacturing additives, can leach out of the plastic into the sample and substantially inhibit enzymatic assays and binding studies. As described in recent publications* these chemicals, such as slip agents (e.g. oleamide), can slow down evaporation, skew absorbance readings and lead to erroneous DNA quantification. Some of these slip agents have also been shown to negatively affect the outcome of biological tests like enzyme activity or receptor-binding assays.

*Clear information and certification on additives and trace metals are available.

Providing excellence for scientific experts ensures unaffected and reproducible results.

To achieve highest reliability and consistency for your experiments, Eppendorf has optimized materials and processes to minimize the risk of interference.

- > Eppendorf liquid handling consumables are made of highest quality, virgin polypropylene (PP) free of biocides, plasticizers and latex
- > Optimized, highly polished molds produce without the use of slip agents like oleamide, erucamide, stearamide
- > Used dyes are free of organic additives and heavy metals



»We have seen substantial inhibition of our enzyme assays by chemicals leaching from disposable plastic consumables. To obtain the best possible reproducibility we use consumables from manufacturers that can confirm the absence of critical manufacturing additives.«

Dr. Andrew Holt
*Department of Pharmacology,
University of Alberta, Canada*



»We need to avoid that contaminants from the plastic material enter the sample and inhibit bacterial growth. The consumables that we use to analyse water samples should be of the highest purity to obtain reliable results.«

Karen Thomsen
*Microbiology Central Laboratory,
Hamburg Water GmbH, Germany*



»Our DNA isolation protocols from both animal and plant material require grinding of tissue prior to and during the cell lysis process. As we work with high numbers of samples, breaking of consumables and subsequent sample loss can be critical. The excellent quality and stability of the 1 mL Deepwell Plates from Eppendorf convinced us as it improved the reliability of our process significantly.«

Dr. Paul Gooding
*Plant Genomics Centre, Australian
Genome Research Facility*

Premium Quality is Our Standard

Certified quality and purity

Continuous quality assurance throughout the entire production process—from the initial material to the finished product.

- > Manufactured from carefully selected, purest raw materials which comply with international purity criteria (FDA 21 CFR§177.1520 »Olefin Polymers«, FDA 21 CFR§178.2010 »Antioxidants and Stabilizers for Polymers«)
- > Fully automatic production in clean room conditions according to VDI 2083 class 6 and to U.S. Fed.Stand. 209D class 100.000—continuous bioburden and particle monitoring of the plant for impeccable production environment
- > SOP present for storage and every production step—elimination of human contamination risk
- > Frequent production tool checks plus regular in-process quality and functional checks for consistent quality guarantee and absolute reproducibility
- > Full traceability for each product—full control, each product to be traced back to material lot

Proof of compliance with requirements of standards, guidelines and regulations that apply to biological, diagnostic and industrial laboratories

- > Lot-specific purity certificates, documenting lot-specific testing by an independent lab for products with purity grades
 - > Eppendorf Biopur®
 - > PCR clean
 - > Sterile
 - > PCR clean and Sterile
- > By entering the respective lot number, lot-specific certificates can be downloaded
- > General quality certificates as e.g.
 - > General Certificates of Purity for PCR clean and Eppendorf Biopur products
 - > Certificates of Quality – Statement on BSE/TSE
- > Product specific certificates
 - > For values of trace metal
 - > For absence of surface active additives
 - > With special focus, e.g. filter efficiency for ep Dualfilter T.I.P.S.®
- > Certificates verifying Quality Management/ISO conformity



> All certificates may be downloaded on www.eppendorf.com/certificates

A System You Can Rely On

Tailored to your application needs

Eppendorf has set industry standards in consumable purity levels. Building on the famous Eppendorf Quality, four additional purity grades tailor Eppendorf Consumables to various applicational needs—Sterile, PCR clean and Biopur.

You have highest expectations on consistency and reproducibility? Make Eppendorf consumables your choice!

Paving the way for modern lab standards

All consumables are subject to internal process controls on function, tightness, precision, transparency, low wetting property and high chemical and thermal resistance. Additionally Sterile, PCR clean and Biopur products are tested by an independent, external analytical laboratory for compliance—lot by lot. The results of these testings are certified and can be downloaded from www.eppendorf.com/certificates



Eppendorf Quality™

Continuous quality controls

Example applications:
All routine lab applications




Sterile

Continuous, lot-specific quality controls

Certified:
> Pyrogen-free
> Sterile

Example applications:
Microbiology and cell culture applications

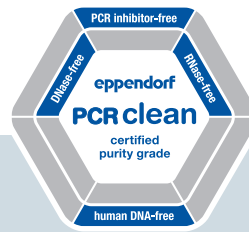


PCR clean

Continuous, lot-specific quality controls

Certified:
> Human DNA-free
> DNase-free
> RNase-free
> PCR inhibitor-free

Example applications:
The perfect choice for isolation and analytics (PCR/qPCR, microarrays) of RNA, DNA, nucleic acids

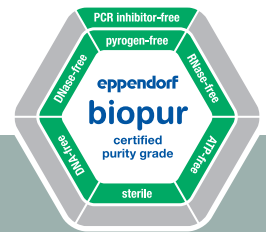

Biopur®

Continuous, lot-specific quality controls

Certified:
> Human DNA-free
> Bacterial DNA-free
> DNase-free
> RNase-free
> PCR inhibitor-free
> ATP-free
> Pyrogen-free
> Sterile

Example applications:
For highest purity demands in cell culture, nucleic acid analytics, or hygiene monitoring



| | Sterile* | PCR clean | PCR clean and Sterile* | Biopur* |
|--|----------|-----------|------------------------|---------|
| Lot testing (certified) for the following purity criteria | | | | |
| Human DNA-free | | • | • | • |
| DNA-free (Human- and bacterial DNA) | | | | • |
| DNase-free | | • | • | • |
| RNase-free | | • | • | • |
| PCR inhibitor-free | | • | • | • |
| ATP-free | | | | • |
| Pyrogen-free (endotoxin-free) | • | | • | • |
| Sterile (Ph. Eur./USP) | • | | • | • |
| Methods (examples) | | | | |
| Bacteria and yeast culture | ✓ | | ✓ | ✓ |
| Cell and tissue culture | ✓ | | ✓ | ✓✓ |
| Isolation and storage of DNA | | ✓✓ | ✓ | ✓ |
| Isolation and storage of RNA | | ✓ | ✓ | ✓✓ |
| DNA analysis (PCR, qPCR, restriction analysis, hybridization, microarrays, sequencing) | | ✓✓ | ✓ | ✓ |
| Mitochondrial DNA analysis | | | | ✓✓ |
| Bacterial DNA analysis | | | | ✓✓ |
| RNA analysis | | ✓ | ✓ | ✓✓ |
| Application Areas (Examples) | | | | |
| Molecular biology | | ✓✓ | ✓ | ✓ |
| Microbiology | ✓ | | ✓ | ✓ |
| Cell technology: | ✓ | | ✓ | ✓✓ |
| > Stem cell research | | | | |
| > Transgenic animals/plants | | | | |
| Research: | | | ✓ | ✓✓ |
| > Medical research | | | | |
| > Agriculture and aquaculture research | | | | |
| Quality control: | | | ✓ | ✓✓ |
| > Food and beverage | | | | |
| > Water supply | | | | |
| > Environmental monitoring | | | | |
| Forensic | | ✓ | ✓ | ✓✓ |

✓ Recommended, ✓✓ Highly recommended

* Increased safety due to individually packaged/single-blistered products.

The Perfect Fit—epT.I.P.S.®

Each of your valuable samples deserve best treatment. See for yourself how Eppendorf pipette tips will save time and reduce costs.

With respect to material, fit, design and operating forces our pipette tips set new standards. The close environment of each sample should be adapted to its specific quality and purity needs. This can involve a specific purity level or the absence of certain substances, but also stability, reliability or geometry. The epT.I.P.S. pipette tips from Eppendorf are designed to cover all of the specific needs of your samples.


Our epT.I.P.S.—Eppendorf Totally Integrated Pipetting System—have been developed to work in perfect combination with our pipettes.

The results are reduced-force tip attachment and ejection with complete sealing. Additionally, the universal nose cone design allows the tips to be used with pipettes from other manufacturers.

In the Eppendorf owned production facilities in the north of Germany, we can maintain the highest standards in the selection and processing of plastic materials. Ideal wetting properties, high transparency and special certified purity levels are the visible expression of this production philosophy.

Every tip is manufactured to match specifically to its respective Eppendorf pipettes, thus ensuring the maximum precision and reliability you have come to expect and to rely on.

How can you recognize original eppendorf pipette tips? The unmistakable feature of our tips is the »ring« of drop-shaped relief elements at the edge of the tip crown. Additionally, every tip has the »Eppendorf« name embossed on the edge.

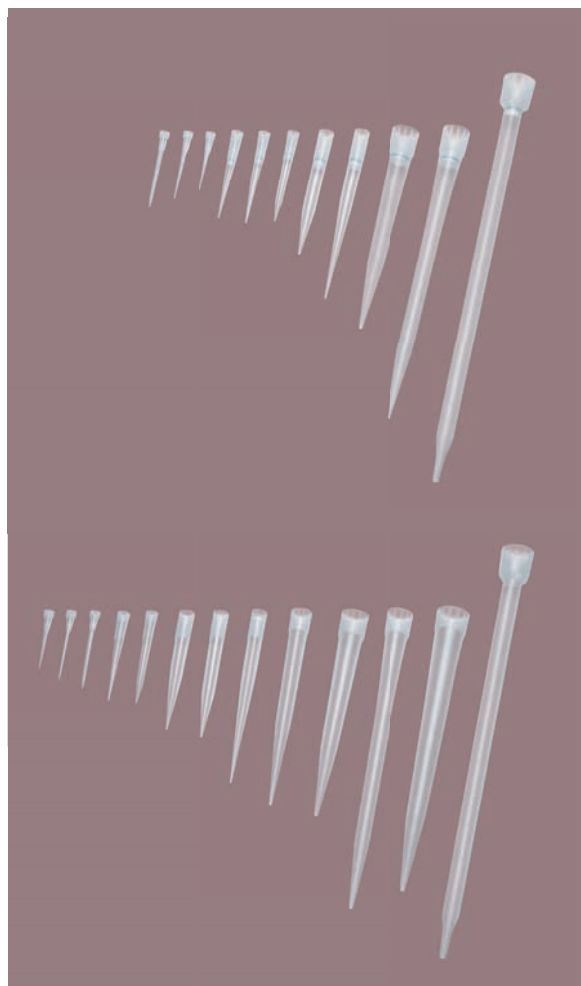

epT.I.P.S.[®] Racks

- > Eppendorf Biopur[®] pipette tips provide maximum biological purity. Guaranteed PCR clean, sterile, pyrogen-free, ATP-free, bacterial DNA-free, they meet the most stringent requirements of the medical, pharmaceutical and foodstuffs industry as well as those of molecular biology and cell technology
- > Continuous control of each batch by an independent laboratory—Batch-related certificates available on request, or go to www.eppendorf.com
- > Packaged in racks of 96, 48 or 24 tips



epT.I.P.S.[®] Box/ epT.I.P.S.[®] Set

- > Contamination-free transfer of trays to the working box
- > System optimized for use with multichannel pipettes
- > Color-coded trays for simple identification of tips and matching Eppendorf pipette
- > Tips can be attached to the pipette from the refill trays
- > Refill trays and epT.I.P.S. boxes are entirely autoclavable for future use

epT.I.P.S.[®] Reloads

- > All components are 100 % recyclable
- > Reduced waste as compared with disposable racks
- > The refill system, depending on tip size, is packaged as either dual-sided or in stack form
- > Reloads are available in two purity levels: Eppendorf Quality and PCR clean

epT.I.P.S.[®] Standard

- > Original, high-quality Eppendorf pipette tips packaged in resealable bags
- > Available in sizes from 10 μ L to 10 mL
- > 200 μ L, 300 μ L and 1,000 μ L tips are also available color-coded yellow and blue

epT.I.P.S.[®] Singles

- > Individually wrapped pipette tips in Eppendorf Biopur: guaranteed sterile, RNase-, DNA-, ATP- and pyrogen-free
- > Batch number and expiration date printed on each blisterpack
- > Continuous quality control of each batch by an independent laboratory—Batch-specific purity certificates available upon request, or go to www.eppendorf.com



Long-Distance Tips

Xtra long for best pipetting results!

With the epT.I.P.S. L pipette tips you will always receive best results while pipetting from and into conical tubes, high reagent bottles, narrow deep vessels, cell culture flasks or deepwell plates. The long and slim design of these elongated pipette tips give you free access to your sample with reduced risk of touching the sides of tubes or wells. Cross-contamination during pipetting can be avoided to a minimum.

Product features and benefits

- > Safe access to Eppendorf 5 mL tubes, 15 mL conical tubes, test tubes, cell culture flasks, deepwell plates and other deep vessels
- > Highest protection for pipette and sample with ep Dualfilter T.I.P.S. L
- > Available in the purity standards Eppendorf Quality, PCR clean, PCR clean/Sterile and Eppendorf Biopur

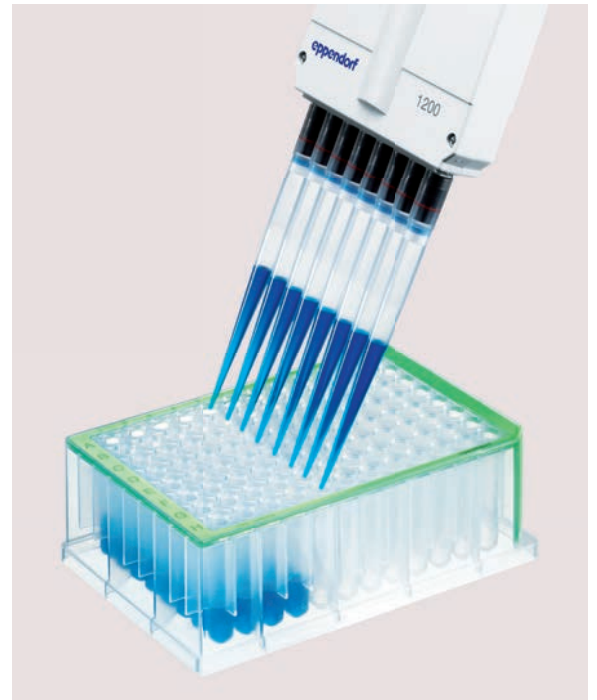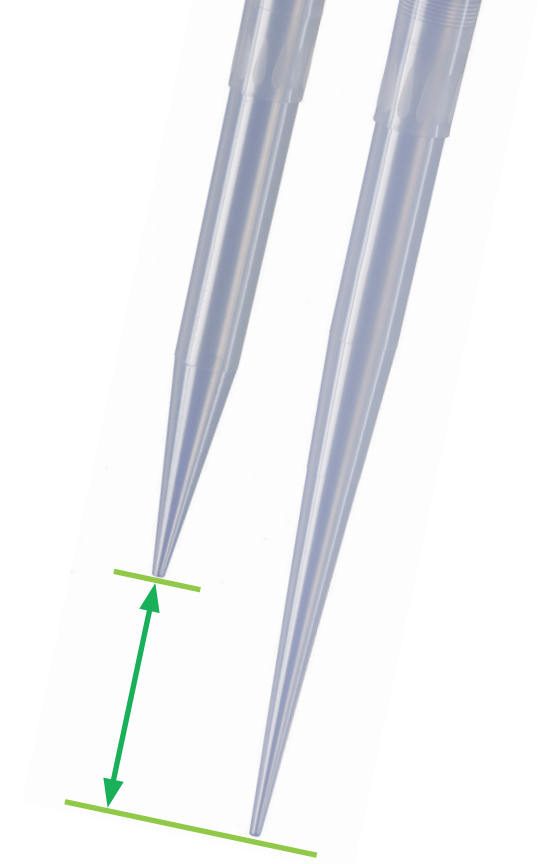

> epT.I.P.S. L are available in volumes of: 0.5–20 μ L L 46 mm long, 50–1,250 μ L L 103 mm long, 0.2–5 mL L 175 mm long* and 0.5–10 mL L 243 mm long

* Coming soon.

ep Dualfilter T.I.P.S.[®]

Two filter layers are better than one.

Eppendorf ep Dualfilter T.I.P.S. are the first filter tips with a two-phase filter for contamination protection.

The unmistakable blue and white filter layers are made of flexible, hydrophobic material to fit perfectly in the tip cone and retain practically 100% of all aerosols* and biomolecules. This unique filtering effect is achieved using various well-defined pore sizes in the two filter layers.

The white layer that faces the sample retains drops, splashes and aerosols. The blue layer forms a highly-effective double barrier that reliably binds biomolecules.

* An aerosol is a colloid of fine solid particles or liquid droplets, in air or another gas.

The air flow rate through the filter is the same as with single-layer filters—it simply makes your rapid pipetting tasks much safer.

The ep Dualfilter T.I.P.S. hydrophobic filter surface also enables easier and complete recovery of samples.

Product features und benefits

- > Two filter layers provide dual protection
- > Dual protection for pipette and sample
- > Dual protection from aerosols and biomolecules
- > Free of PCR inhibitor additives
- > Eppendorf PCR clean, Sterile and pyrogen-free
- > Batch certified

Dual protection against contamination right from the start.

ep Dualfilter T.I.P.S. are manufactured to the highest possible quality from pure, non-recycled materials under cleanroom conditions. They are characterized by defined flow dynamics, low wettability and high thermal stability.

ep Dualfilter T.I.P.S. are sterilized by electron beams and certified pyrogen-free and PCR clean (free from human DNA, DNase, RNase and PCR-inhibiting substances).

Labels: Pipette, Pipette cone, Aerosols, Pipette tip, Sample, Filter.

Without a filter **1**, the pipette is exposed to contamination by samples and aerosols. Conventional single-layer filters **2** do not fully block particles and molecules. Only ep Dualfilter T.I.P.S. **3** provide reliable protection even against the finest impurities.

> Individual batch certification can be downloaded from www.eppendorf.com/certificates

Scientifically Proven


Excellent protection with ep Dualfilter T.I.P.S.[®]

For the introduction of the ep Dualfilter T.I.P.S.[®] filter tips in 2006, filter tips from various manufacturers were tested in a blind trial at the Fraunhofer Institute for Toxicology and Experimental Medicine (ITEM) in Hanover, Germany. Testing focused on the prevention of tip cone contamination by aerosols. Filter effect with regard to salt aerosols and biomolecules (DNA fragments) was examined. Quantitative determination of DNA fragments was performed using real-time PCR.

The results indicate that ep Dualfilter T.I.P.S. filter tips are much better at protecting against aerosols than the other filter tips tested during the trial. With regard to particle permeability, ep Dualfilter T.I.P.S. are 55 to 677 times more effective. ep Dualfilter T.I.P.S. were 21 to 600 times more effective when protection against DNA fragments was examined for the first time.

Both filters, ep Dualfilter T.I.P.S. and ep Dualfilter T.I.P.S. SealMax, are made of the same Polyethylene (PE) and have comparable pore size structure. Thus a further study with ep Dualfilter T.I.P.S. SealMax was conducted based on the standard EN 1822. This standard deals with filtration performance testing of filters as used, for example, for applications in clean room technology or pharmaceutical industry. The results of the filter efficiency testing done by a certified institute showed a minimum particle collection efficiency of the ep Dualfilter T.I.P.S. SealMax of 99.5% with NaCl aerosol particle sizes of 0.05–0.5 µm.

Aerosol/particle penetration of various filter tips*

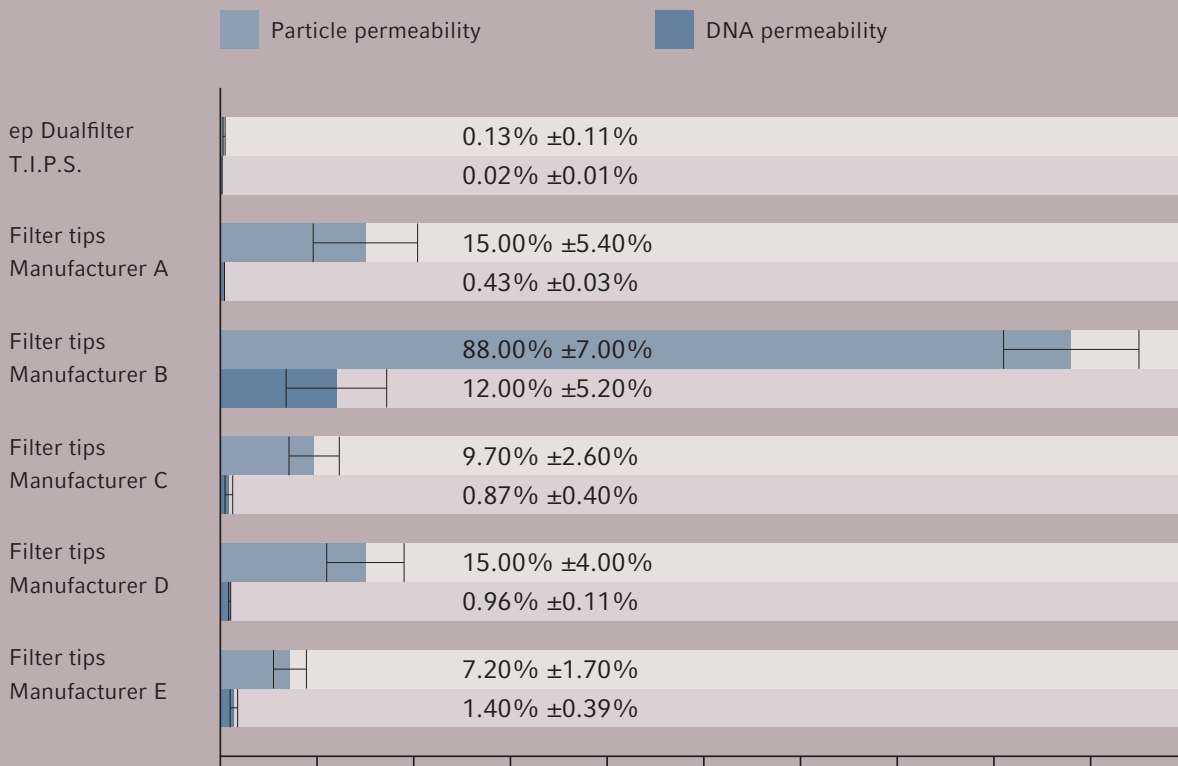

* Measurement report »Etablierung einer Methode zur Bestimmung des Penetrationsgrades von Pipettenspitzen mit eingebauten Filtern mittels Biomolekülen« (Establishment of a method to determine the degree of penetration of pipette tips with integrated filters by means of biomolecules), Roman Halter, Gerhard Pohlmann, Fraunhofer Institute for Toxicology and Experimental Medicine (ITEM), Hanover, Germany, July 2005


Don't Panic

Your pipette is safe with ep Dualfilter T.I.P.S. SealMax
Eppendorf's new ep Dualfilter T.I.P.S. SealMax filter tips provide comprehensive protection against contamination of your pipette and sample. When it comes to an accidental over-pipetting situation, the new violet layer of ep Dualfilter T.I.P.S. SealMax becomes a reliable barrier against sample liquid—no liquid will pass through the filter! Your pipette is safe at all times! Additionally, nearly 100% of aerosols* and biomolecules are retained while the hydrophobic white layer protects against splashes and droplets. Further information on specific sample recovery features, PCR inhibition and aerosol protection are available on your local Eppendorf website.

* An aerosol is a colloid of fine solid particles or liquid droplets, in air or another gas.

Product features and benefits

- > Reliable blocking of liquid for comprehensive pipette protection
- > Maximum protection against aerosols and biomolecules for pipette and sample
- > No PCR-inhibition in case of sample contact
- > Available in PCR clean, Sterile (sterile, pyrogen-free)
- > Lot specific purity certificates

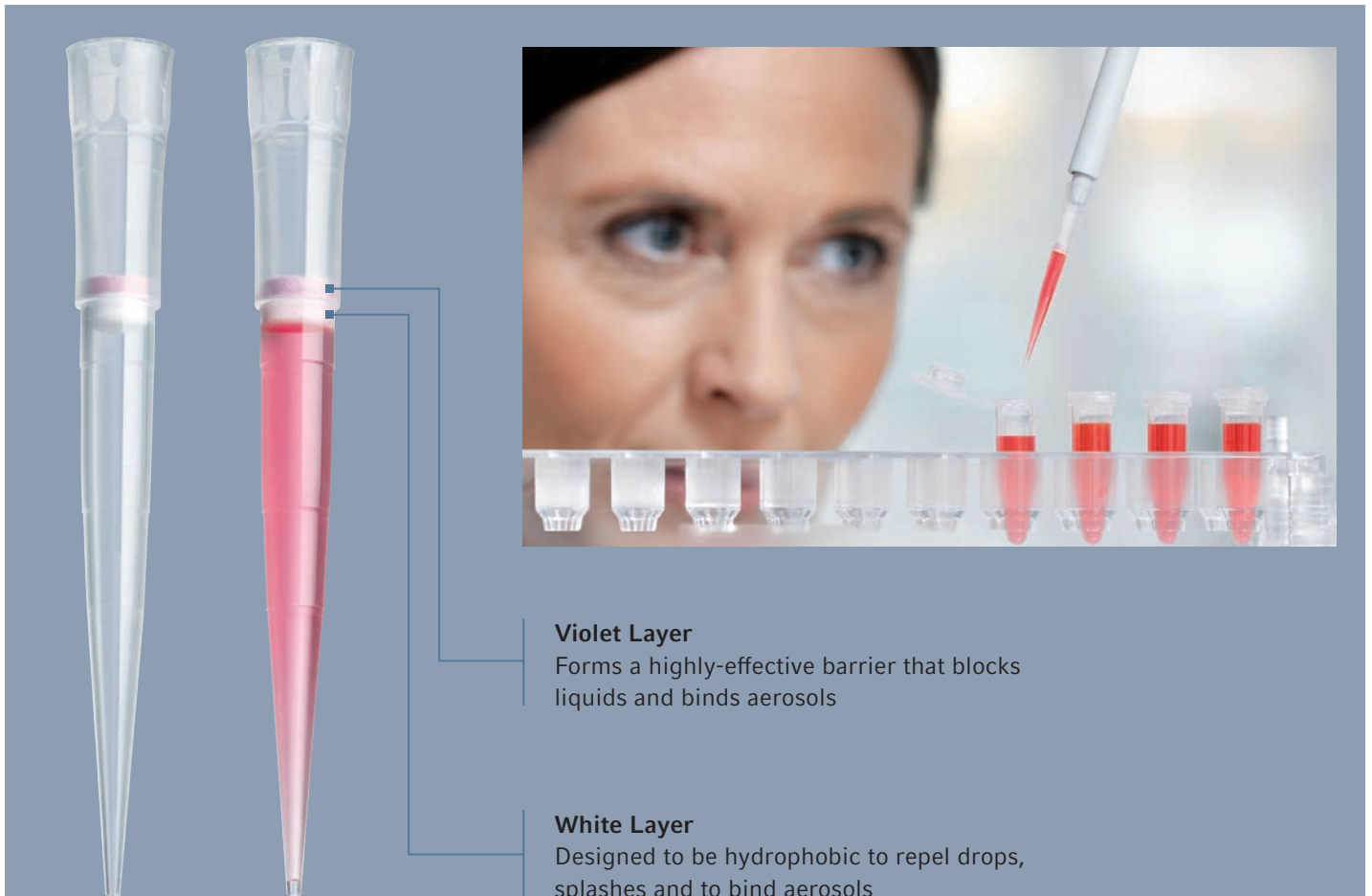

epT.I.P.S.[®] LoRetention

Increased sensitivity with maximum reproducibility

Do you depend on your precious samples? Then don't leave any valuable residues in your tips when pipetting solutions that contain detergents! Based on a unique and innovative method of material modification, the new Eppendorf epT.I.P.S. LoRetention pipette tips allow you to transfer almost 100 % of the liquid picked up («pearl effect») — for maximum recovery with unbeaten precision.


The pipetting of liquids containing detergents is ubiquitous in modern laboratory processes. The lower surface tension involved has a marked influence on the flow and dosing properties of samples, thus making it noticeably more difficult to perform reproducible work with these liquids.

The «Pearl Effect»

The ultrahydrophobic extremely homogenous surface of the epT.I.P.S. LoRetention pipette tips is achieved through an innovative treatment at the molecular level—the Pearl Effect Technology.

The tips are uncoated, free from additives, and do not bleed into the sample. Liquids roll off completely, so that only a tiny drop remains in the tip. In comparison, the adhesive layer of liquid in standard tips retains considerably more sample material.

For users of the new epT.I.P.S. LoRetention pipette tips, this unique material property guarantees maximum sample recovery, improved reproducibility and higher sensitivity. Available in purity grade of PCR clean and Eppendorf Quality, also as ep Dualfilter T.I.P.S..

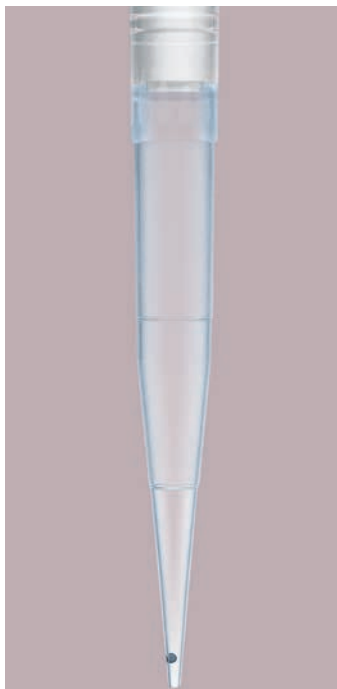

The «Pearl Effect» — Maximum recovery with epT.I.P.S. LoRetention



Simple differentiation of racks and reloads through a new, clear lid with a white seal. The reusable seal supplied with the reloads can

also be used to label existing epT.I.P.S. boxes. When open, the writing «LoRetention» can be seen from the inside.

The Facts Speak for Themselves

Reduce costs, increase reproducibility

In particular with sensitive PCR and real-time PCR applications, low-retention surfaces can increase reproducibility and notably reduce the loss of expensive reagents—depending on the reagent solution, this can be by up to 5%!

Versatile application


epT.I.P.S. LoRetention pipette tips are ideally suited for all applications with samples that contain detergents:

- > PCR and real-time PCR preparation of enzyme solutions and master mixes
- > Isolation, purification and denaturation of proteins
- > Restriction digestions and ligations
- > Pipetting DNA ladders during gel electrophoresis


epT.I.P.S. LoRetention perform considerably better in terms of precision and sample recovery than low-retention tips from other manufacturers.

Product features and benefits

- > Ultrahomogeneous surface for maximum reproducibility
- > No retention of liquid when detergent containing for reagent and cost saving
- > Significantly reduced build up of foam when pipetting
- > No coating thus no assay interference, no influencing of reaction/analysis by leachables
- > Chemically resistant and autoclavable

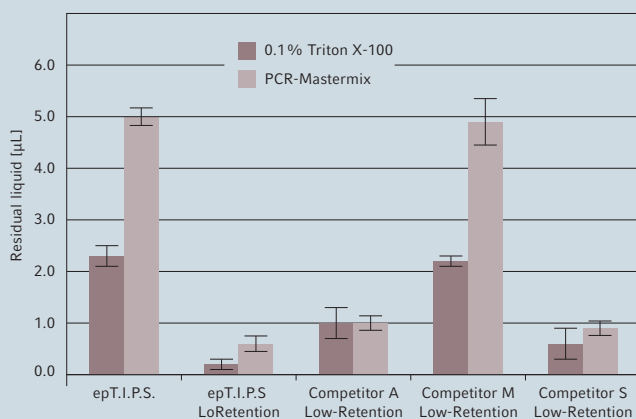

Comparison of residual liquid from epT.I.P.S. LoRetention with standard epT.I.P.S. pipette tips and low-retention tips of other manufacturers. 200 µL tips were used to pipette a commercial PCR master mix and the detergent 0.1 % Triton X-100. Once the liquid had been dispensed, the residual liquid in the tips was determined. The measurements were repeated several times and the standard deviation determined. epT.I.P.S. LoRetention pipette tips resulted in the lowest residual liquid, as well as the smallest standard deviation.

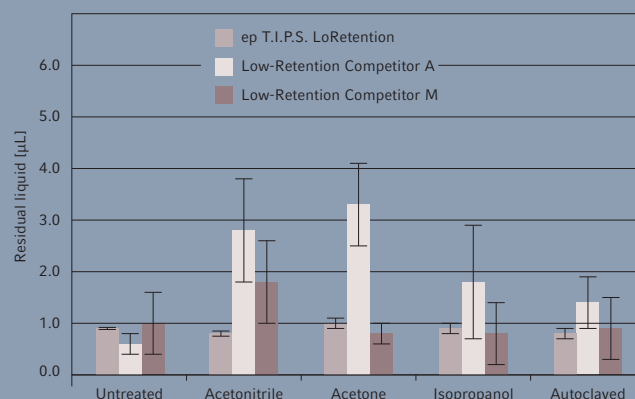

Comparison of resistance of the low-retention property when subjected to solvents or autoclaving. 200 µL epT.I.P.S. LoRetention pipette tips and low-retention tips of other manufacturers were autoclaved or treated with the solvent specified. An enzymatic buffer containing detergents was then pipetted and the residual moisture was determined. epT.I.P.S. LoRetention pipette tips resulted in reproducible low levels of residual liquid with minimum standard deviation. Low-retention tips from other manufacturers are less resistant to chemicals or autoclaving; this resulted in an increase in residual liquid and standard deviation.

* Further information can be found in our Application Notes 218 (Nov. 2009) and 192 (July 2008) at www.eppendorf.com/consumables.


Loading and Filling

Eppendorf Microloader

These unique tips are ideally suited to fill microcapillaries used for microinjection. The extremely long, fine and flexible tip provides also the ideal solution for all kind of applications in which additional reach is needed while pipetting smallest volumes.

Product features and benefits

- > Rack package can be sterilized
- > Ideal for recovering surplus solution from the capillary

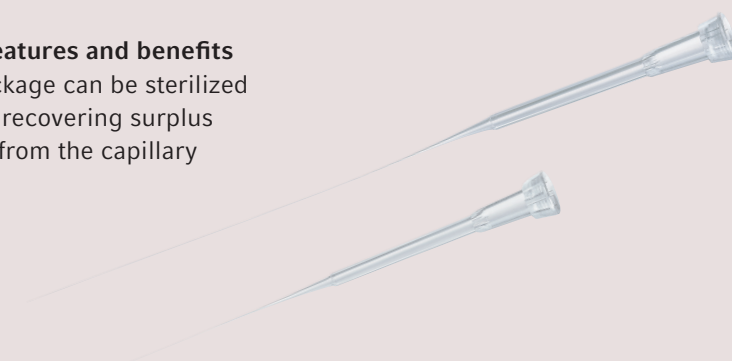


Eppendorf GELoader®

Eppendorf's GELoader Tips have been designed to simplify the loading of samples onto polyacrylamide gels. These flexible, long and narrow tips prevent the gels from being damaged whilst allowing optimal handling of smallest volumes.

Product features and benefits

- > Special tip for gel electrophoresis
- > GELoader Tips and rack autoclavable (121°C, 20 min.)
- > Highest precision and accuracy when used with Eppendorf pipettes for 0.5 to 10 µL (gray control button)

Rear filling of Microcapillary using Microloader and Research plus pipette

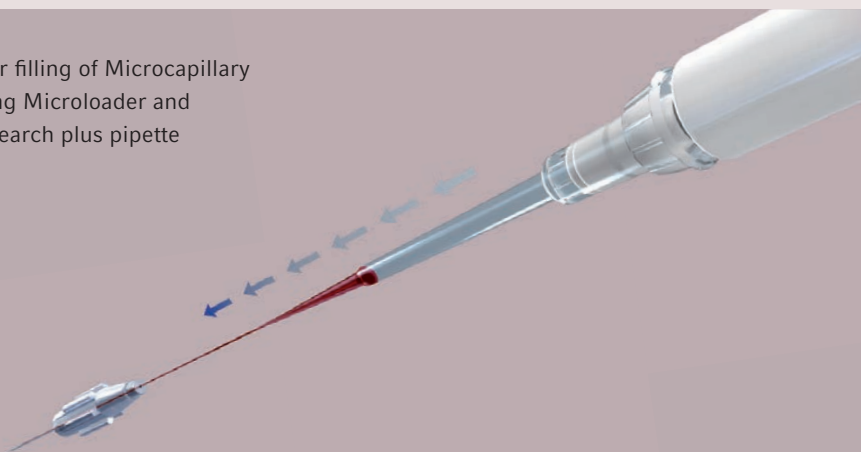

Unmistakably Eppendorf—The new Eppendorf Combitips advanced® also feature the characteristic Eppendorf 3D design on the upper rim of the Combitip, allowing you to easily and confidently recognize the Eppendorf Original!



Ergonomic design—Unique funnel geometry prevents damage to gloves and ensures comfortable handling



Elongated tips (for 2.5 mL, 5 mL and 10 mL)—Complete emptying of all common tubes prevents sample loss (see table on page 17)

Time for a New Original



Eppendorf Combitips advanced®

The invention of the Eppendorf hand dispenser Multipipette®/ Repeater®* and the Eppendorf Combitips® marks a milestone in the field of Liquid Handling!

Executing long pipetting respectively dispensing series is remarkably simpler and faster. With the innovative sensor technology for automatic Combitip recognition, this dispensing system has become an indispensable tool for every laboratory.

The increasingly high demands of modern laboratory work have also heightened the requirements for a high-precision dispensing tip. Our experts have thoroughly optimized the Combitips using sophisticated engineering.

The result: A revolutionary 360° evolution! Setting a completely new standard in dispensing tips.

* in USA

The »Combitips« principle

- > Positive displacement principle (comparable to a syringe)
- > High-precision-dispensing regardless of the density and viscosity of the liquid
- > Prevents aerosol contamination with hermetically sealed piston for secure dispensing and provides protection from radioactive and toxic substances
- > Quick dispensing of long series with precise, repeated dispensing of identical volumes (in combination with the Multipipette/Repeater®)
- > Individually color coded—Quick identification of the desired Combitips facilitates the workflow
- > Variety and selection—With 9 volume sizes (0.1 mL to 50 mL) and several purity levels you will always find the perfect Combitip for your application. The tips in Eppendorf Biopur are individually blister-wrapped and feature an access tab which makes them easier to open, even with gloves

Technical Specifications

Combitips advanced®

Technical specifications

| Combitips advanced | Min./max. volume | Increment/step size | Max. dispensing volume | Volume | Inaccuracy*1 | Imprecision*1 |
|---|------------------|---------------------|------------------------|-----------|--------------|---------------|
| for Multipipette M4 and Combitips advanced system | | | | | | |
| 0.1 mL | 1 µL | 1 µL | 20 µL | 2 µL | ±1.6 % | ±3.0 % |
| | 20 µL | | | 20 µL | ±1.0 % | ±2.0 % |
| 0.2 mL | 2 µL | 2 µL | 40 µL | 4 µL | ±1.3 % | ±2.0 % |
| | 40 µL | | | 40 µL | ±0.8 % | ±1.5 % |
| 0.5 mL | 5 µL | 5 µL | 100 µL | 10 µL | ±0.9 % | ±1.5 % |
| | 100 µL | | | 100 µL | ±0.8 % | ±0.6 % |
| 1 mL | 10 µL | 10 µL | 200 µL | 20 µL | ±0.9 % | ±0.9 % |
| | 200 µL | | | 200 µL | ±0.6 % | ±0.4 % |
| 2.5 mL | 25 µL | 25 µL | 500 µL | 50 µL | ±0.8 % | ±0.8 % |
| | 500 µL | | | 500 µL | ±0.5 % | ±0.3 % |
| 5 mL | 50 µL | 50 µL | 1,000 µL | 100 µL | ±0.6 % | ±0.6 % |
| | 1,000 µL | | | 1,000 µL | ±0.5 % | ±0.25 % |
| 10 mL | 100 µL | 100 µL | 2,000 µL | 200 µL | ±0.5 % | ±0.6 % |
| | 2,000 µL | | | 2,000 µL | ±0.5 % | ±0.25 % |
| 25 mL | 250 µL | 250 µL | 5,000 µL | 500 µL | ±0.4 % | ±0.6 % |
| | 5,000 µL | | | 5,000 µL | ±0.3 % | ±0.25 % |
| 50 mL | 500 µL | 500 µL | 10,000 µL | 1,000 µL | ±0.3 % | ±0.5 % |
| | 10,000 µL | | | 10,000 µL | ±0.3 % | ±0.3 % |
| for Multipipette stream/Xstream/E3/E3x and Combitips advanced system | | | | | | |
| 0.1 mL | 1 µL | 0.1 µL | 0.1 mL | 10 µL | ±1.6 % | ±2.5 % |
| | | | | 50 µL | ±1.0 % | ±1.5 % |
| | 100 µL | | | 100 µL | ±1.0 % | ±0.5 % |
| 0.2 mL | 2 µL | 0.2 µL | 0.2 mL | 20 µL | ±1.3 % | ±1.5 % |
| | | | | 100 µL | ±1.0 % | ±1.0 % |
| | 200 µL | | | 200 µL | ±1.0 % | ±0.5 % |
| 0.5 mL | 5 µL | 0.5 µL | 0.5 mL | 50 µL | ±0.9 % | ±0.8 % |
| | | | | 250 µL | ±0.9 % | ±0.5 % |
| | 500 µL | | | 500 µL | ±0.9 % | ±0.3 % |
| 1 mL | 10 µL | 1 µL | 1 mL | 100 µL | ±0.9 % | ±0.55 % |
| | | | | 500 µL | ±0.6 % | ±0.3 % |
| | 1 mL | | | 1,000 µL | ±0.6 % | ±0.2 % |
| 2.5 mL | 2.5 mL | 2.5 µL | 2.5 mL | 250 µL | ±0.8 % | ±0.45 % |
| | | | | 1,250 µL | ±0.5 % | ±0.3 % |
| | 2.5 mL | | | 2,500 µL | ±0.5 % | ±0.15 % |
| 5 mL | 50 µL | 5 µL | 5 mL | 500 µL | ±0.8 % | ±0.35 % |
| | | | | 2,500 µL | ±0.5 % | ±0.25 % |
| | 5 mL | | | 5,000 µL | ±0.5 % | ±0.15 % |
| 10 mL | 100 µL | 10 µL | 10 mL | 1,000 µL | ±0.5 % | ±0.25 % |
| | | | | 5,000 µL | ±0.4 % | ±0.25 % |
| | 10 mL | | | 10,000 µL | ±0.4 % | ±0.15 % |
| 25 mL | 250 µL | 25 µL | 25 mL | 2,500 µL | ±0.3 % | ±0.35 % |
| | | | | 12,500 µL | ±0.3 % | ±0.25 % |
| | 25 mL | | | 25,000 µL | ±0.3 % | ±0.15 % |
| 50 mL | 500 µL | 50 µL | 50 mL | 5,000 µL | ±0.3 % | ±0.5 % |
| | | | | 25,000 µL | ±0.3 % | ±0.2 % |
| | 50 mL | | | 50,000 µL | ±0.3 % | ±0.15 % |

*1 The data for imprecision (random error) and inaccuracy (systematic error) according to EN ISO 8655 only apply when using original Eppendorf Combitips advanced.

Compatibility of Combitips advanced with standard laboratory tubes

| Combitips advanced | Eppendorf Safe-Lock Tubes | | | Conical tubes | | Eppendorf Deepwell Plates | | |
|--------------------|---------------------------|--------|--------|---------------|-------|---------------------------|------------|------------|
| | 0.5 mL | 1.5 mL | 2.0 mL | 15 mL | 50 mL | 96/500 µL | 96/1000 µL | 96/2000 µL |
| 0.2 mL | + | + | + | – | – | + | + | + |
| 0.5 mL | + | + | + | – | – | + | + | + |
| 1 mL | + | + | + | – | – | + | + | + |
| 2.5 mL | ⊕ | ⊕ | ⊕ | ⊕ | ⊕ | + | ⊕ | ⊕ |
| 5 mL | + | + | + | ⊕ | ⊕ | + | + | + |
| 10 mL | + | + | + | – | ⊕ | + | + | + |
| 25 mL | – | + | + | – | + | + | + | + |
| 50 mL | – | + | + | – | – | + | + | + |

⊕ = improved compatibility compared to the Combitip plus generation



Intelligent packaging

With the redesigned dispenser box (for »Eppendorf Quality« and PCR clean) including an extricable chute, Combitips advanced are always within reach. Portioning in 4 bags (25 pcs. each) ensures dust-free storage.

Always easy to reach

The new Combitip Rack features trouble-free and single-hand operation during dispensing work. The rack holds up to 8 Combitips advanced in sizes ranging from 0.1 mL to 10 mL. The transparent tray allows for the sizes of the Combitips to be easily identified using the color coding.

A perfect team

Whether Multipipette M4, Multipipette stream/E3 or Xstream/E3x—the Combitips advanced are optimized for all Multipipette models to create the perfect connection!

Positive Displacement Tips



For Eppendorf Biomaster®

Enjoy contamination-free pipetting with the Positive Displacement Pipette, an adjustable pipette for the volume range from 1 μ L to 20 μ L. With the Positive Displacement Tip, the Eppendorf Mastertip®, the pipette functions according to the positive-displacement principle; thus, the formation of aerosols—a common cause of contamination—is eliminated.

Product features and benefits

- > Ready-for-use piston integrated in tip for contamination-free pipetting in positive displacement mode
- > Ideal for liquids that are viscous (e.g. glycerol) or have high vapour pressure (e.g. ethanol, acetone)
- > 96 positive displacement tips per rack


For Eppendorf Varipette®*

The pipette tips for the Varipette are tailored to different vessels. For example, the Eppendorf Varitips® P is designed for aspirating 1 mL to 10 mL from beakers, and it pipettes according to the positive displacement principle. The Eppendorf Varitips® S for 2.5 mL to 10 mL forms a system with the Varipette, and it is used for aspirating liquids from tall, narrow-neck vessels.

Product features and benefits

- > Varitips P (fig. 1) for aspirating 1 mL to 10 mL from beakers using a positive displacement technique
- > Varitips S (fig. 2): The dosing part, in combination with the Eppendorf Maxitip (fig. 3) forms a system for aspiration of between 2.5 mL and 10 mL liquid from high, narrow-mouth vessels—air displacement principle
- > No carryover as Maxitips are exchanged
- > Valve for Maxitip ensures drip-free dispensing of liquids with a high vapour pressure
- > Maxitip is graduated for dispensing accurate volumes of nonaqueous liquids

* US product names: Maxipettor®, Maxitip P, Maxitip S

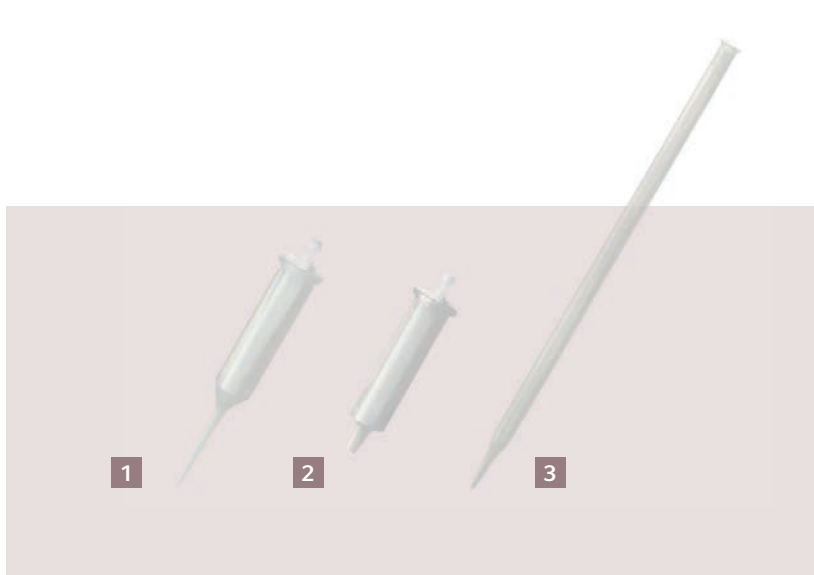

Reliable Unattended Automation



epT.I.P.S.[®] Motion—Eppendorf Totally Integrated Pipetting System for Automation

For routine pipetting with the best precision in the industry it is important to have a perfect system of tips and automate. The epT.I.P.S. Motion have been developed to work in perfect combination with our epMotion[®] platform.

The epT.I.P.S. Motion racks consist of two options for simple integration into the epMotion liquid handling workstation: A single use box with color coded trays for easy volume identification and a TipHolder adapter to use the trays as a Reload-System. With color coded trays you can quickly and easily see the volume range of the tip and assure an even higher level of sample safety. The new Reload trays also offer a more environmentally friendly option and will be delivered in a PET blister packaging with a sealed cap assuring the same well established quality as all Eppendorf tips.

The optional TipHolder adapter (autoclavable) replaces the normal tip box component and reduces waste by up to 40 %. A perfect match for the most demanding user.

Product features and benefits

More safety


- > Each tip is checked for straightness before packaging
- > Color-coded trays for direct volume identification
- > Dispensing tool design assures an optimal fit each time
- > Optical sensor automatically identifies tip type

Flexible handling

- > Three volume ranges available with and w/o filter in multiple purity grades
- > Easily exchangeable, no additional labware files required
- > Available as SafeRacks for tip reuse with individual tip compartments

Eco-friendly Reload System

- > Easy conversion with TipHolder adapter
- > TipHolder adapter is autoclavable
- > Reduces waste by up to 40 %

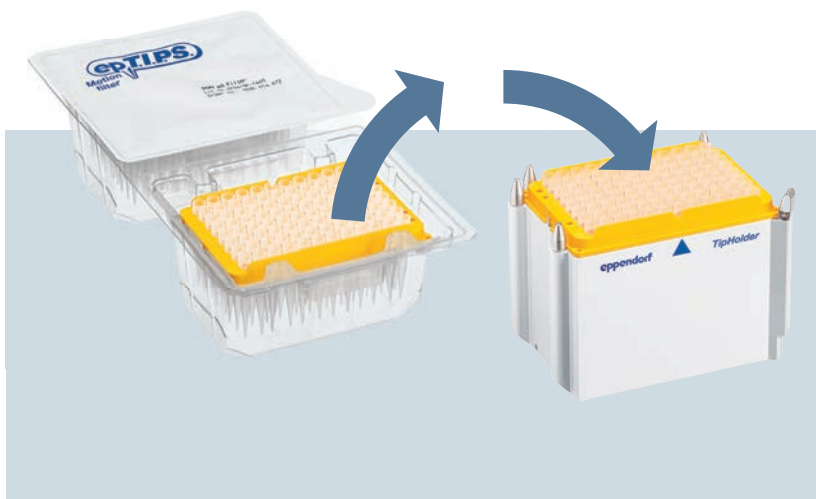

**Volume range/epT.I.P.S.[®]
epT.I.P.S.[®] LoRetention**
(all tips shown are actual size)

Standard

Reloads

Reloads



0.1–10 µL, 34 mm ■ dark gray



Eppendorf Quality™

2 bags of 500 tips = 1,000 tips

■ 0030 000.811

Eppendorf Quality™

10 trays of 96 tips = 960 tips

■ 0030 073.363
■ 0030 072.049 LoRetention

PCR clean

10 trays of 96 tips = 960 tips

■ 0030 073.746
■ 0030 072.006 LoRetention

0.1–20 µL, 40 mm ■ medium gray



2 bags of 500 tips = 1,000 tips

■ 0030 000.838

10 trays of 96 tips = 960 tips

■ 0030 073.380

10 trays of 96 tips = 960 tips

■ 0030 073.762

0.5–20 µL L, 46 mm □ light gray



2 bags of 500 tips = 1,000 tips

□ 0030 000.854

10 trays of 96 tips = 960 tips

□ 0030 073.401
□ 0030 072.057 LoRetention

10 trays of 96 tips = 960 tips

□ 0030 073.789
□ 0030 072.014 LoRetention

2–200 µL, 53 mm ■ yellow



2 bags of 500 tips = 1,000 tips

■ 0030 000.889
■ 0030 000.870 yellow

10 trays of 96 tips = 960 tips

■ 0030 073.428
■ 0030 072.065 LoRetention

10 trays of 96 tips = 960 tips

■ 0030 073.800
■ 0030 072.022 LoRetention

20–300 µL, 55 mm ■ orange



2 bags of 500 tips = 1,000 tips

■ 0030 000.900
■ 0030 000.897 yellow

10 trays of 96 tips = 960 tips

■ 0030 073.444

10 trays of 96 tips = 960 tips

■ 0030 073.827

50–1,000 µL, 71 mm ■ blue



2 bags of 500 tips = 1,000 tips

■ 0030 000.927
■ 0030 000.919 blue

10 trays of 96 tips = 960 tips

■ 0030 073.460
■ 0030 072.073 LoRetention

10 trays of 96 tips = 960 tips

■ 0030 073.843
■ 0030 072.030 LoRetention

50–1,250 µL, 76 mm ■ green



2 bags of 500 tips = 1,000 tips

■ 0030 000.935

10 trays of 96 tips = 960 tips

■ 0030 073.487

10 trays of 96 tips = 960 tips

■ 0030 073.860

Box

Set

Singles

Racks



Eppendorf Quality™

Eppendorf Quality™

Eppendorf Biopur® (sterile)

Eppendorf Biopur® (sterile)

1 reusable box incl. 96 tips

■ 0030 073.002

1 reusable box
incl. 5 trays of 96 tips

■ 0030 073.207
■ 0030 072.251 LoRetention

1 reusable box incl. 96 tips

■ 0030 073.029

1 reusable box
incl. 5 trays of 96 tips

■ 0030 073.223

1 reusable box incl. 96 tips

□ 0030 073.045

1 reusable box
incl. 5 trays of 96 tips

□ 0030 073.240
□ 0030 072.260 LoRetention

1 reusable box incl. 96 tips

■ 0030 073.061

1 reusable box
incl. 5 trays of 96 tips

■ 0030 073.266
■ 0030 072.278 LoRetention

1 reusable box incl. 96 tips

■ 0030 073.088

1 reusable box
incl. 5 trays of 96 tips

■ 0030 073.282

1 reusable box incl. 96 tips

■ 0030 073.100

1 reusable box
incl. 5 trays of 96 tips

■ 0030 073.304
■ 0030 072.286 LoRetention

1 reusable box incl. 96 tips

■ 0030 073.126

1 reusable box
incl. 5 trays of 96 tips

■ 0030 073.320

100 tips, individually wrapped

■ 0030 010.019

100 tips, individually wrapped

■ 0030 010.035

100 tips, individually wrapped

■ 0030 010.051

5 racks of 96 tips = 480 tips

■ 0030 075.005

5 racks of 96 tips = 480 tips

■ 0030 075.021

5 racks of 96 tips = 480 tips

■ 0030 075.048

5 racks of 96 tips = 480 tips

■ 0030 075.064

5 racks of 96 tips = 480 tips

■ 0030 075.080

Volume range/epT.I.P.S.®
(all tips shown are actual size)

Standard

Reloads

Reloads



Eppendorf Quality™

Eppendorf Quality™

PCR clean

50–1,250 µL L, 103 mm ■ dark green

2 bags of 500 tips = 1,000 tips
■ 0030 000.730

10 trays of 96 tips = 960 tips
■ 0030 073.606

10 trays of 96 tips = 960 tips
■ 0030 073.614




0.25–2.5 mL, 115 mm ■ red

5 bags of 100 tips = 500 tips
■ 0030 000.951

10 trays of 48 tips = 480 tips
■ 0030 073.509

10 trays of 48 tips = 480 tips
■ 0030 073.886

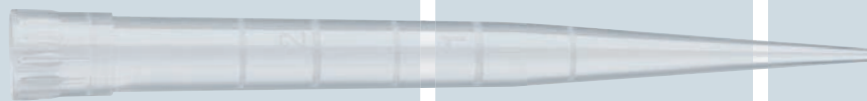

0.1–5 mL, 120 mm ■ violet

5 bags of 100 tips = 500 tips
■ 0030 000.978



0.2–5 mL L, 175 mm ■ violet

3 bags of 100 tips = 300 tips
■ 0030 000.650*



0.5–10 mL, 165 mm ■ turquoise

2 bags of 100 tips = 200 tips
■ 0030 000.765







0.5–10 mL L, 243 mm ■ turquoise
(Image reduced in size)

2 bags of 100 tips = 200 tips
■ 0030 000.781

Fits in deep and narrow test tubes, cell culture bottles and measuring pistons



*Coming soon.
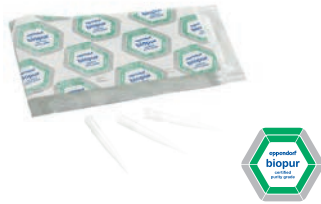
| Box | Set | Singles | Racks |
|---|---|--|---|
|  |  |  |  |
| <p>Eppendorf Quality™</p> | <p>Eppendorf Quality™</p> | <p>Eppendorf Biopur® (sterile)</p> | <p>Eppendorf Biopur® (sterile)</p> |
| <p>1 reusable box incl. 96 tips ■ 0030 073.622</p> | | | <p>5 racks of 96 tips = 480 tips ■ 0030 075.129</p> |
| <p>1 reusable box incl. 48 tips ■ 0030 073.142</p> | <p>1 reusable box incl. 5 trays of 48 tips ■ 0030 073.347</p> | | <p>5 racks of 48 tips = 240 tips ■ 0030 075.102</p> |
| <p>1 reusable box incl. 24 tips ■ 0030 073.169</p> | | | <p>5 racks of 24 tips = 120 tips ■ 0030 075.137</p> |
| | | | <p>5 racks of 24 tips = 120 tips ■ 0030 075.188*</p> |
| | | | <p>5 racks of 24 tips = 120 tips ■ 0030 075.145</p> |
| | | | |



Volume range/ ep Dualfilter T.I.P.S.[®],
ep Dualfilter T.I.P.S.[®] SealMax,
ep Dualfilter T.I.P.S.[®] LoRetention
(all tips shown are actual size)



0.1–10 μL S, 34 mm ■ dark gray



0.1–10 μL M, 40 mm ■ medium gray



0.5–20 μL L, 46 mm □ light gray



2–20 μL , 53 mm ■ yellow



2–100 μL , 53 mm ■ yellow



2–200 μL , 55 mm ■ yellow



20–300 μL , 55 mm ■ orange



50–1,000 μL , 76 mm ■ blue



50–1,250 μL L, 103 mm ■ dark green



Racks



PCR clean, sterile incl. pyrogen-free

10 racks of 96 tips = 960 tips

■ 0030 077.504
■ 0030 077.610 LoRetention
■ 0030 077.806 SealMax

10 racks of 96 tips = 960 tips

■ 0030 077.512

10 racks of 96 tips = 960 tips

□ 0030 077.520
□ 0030 077.628 LoRetention
□ 0030 077.814 SealMax

10 racks of 96 tips = 960 tips

■ 0030 077.539

10 racks of 96 tips = 960 tips

■ 0030 077.547
■ 0030 077.644 LoRetention
■ 0030 077.822 SealMax

10 racks of 96 tips = 960 tips

■ 0030 077.555
■ 0030 077.830 SealMax

10 racks of 96 tips = 960 tips

■ 0030 077.563
■ 0030 077.636 LoRetention
■ 0030 077.849 SealMax

10 racks of 96 tips = 960 tips

■ 0030 077.571
■ 0030 077.652 LoRetention
■ 0030 077.857 SealMax

5 Racks of 96 Tips = 480 Tips

■ 0030 077.750

Volume range/ep Dualfilter T.I.P.S.[®]
ep Dualfilter T.I.P.S.[®] LoRetention
(all tips shown are actual size)

Racks




PCR clean, sterile incl. pyrogen-free

0.1–5 mL, 120 mm ■ violet




5 racks of 24 tips = 120 tips
■ 0030 077.580

0.2–5 mL L, 175 mm ■ violet

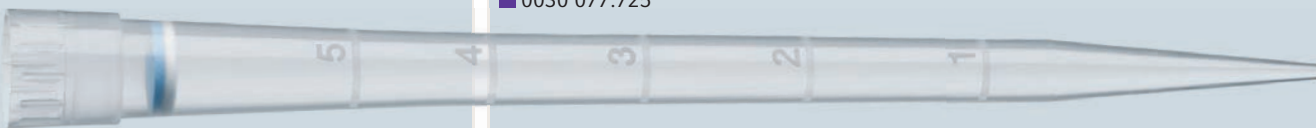

5 racks of 24 tips = 120 tips
■ 0030 077.725*

0.5–10 mL L, 243 mm ■ turquoise
(Image reduced in size)

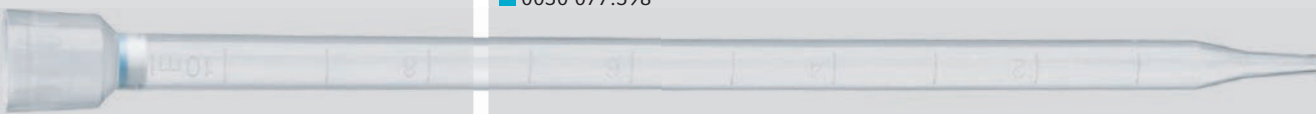


100 tips, individually wrapped
■ 0030 077.598

* Coming soon.

Special tips

Volume range/pipette tip

0.5–20 µL, 62 mm □ light gray

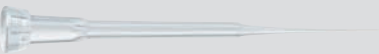

GEloader[®]

2 racks of 96 tips = 192 tips
□ 0030 001.222



0.5–20 µL, 100 mm ■ medium gray



Microloader

2 racks of 96 tips = 192 tips
■ 5242 956.003



20 µL, 52 mm □ light gray



Mastertip[®]

5 racks of 96 tips = 480 tips
□ 0030 001.320



Combitips advanced[®] and Accessories

| Combitips advanced [®] | Color coding | Eppendorf Quality [™] 100 pcs. (4 bags × 25 pcs.) | PCR clean 100 pcs. (4 reclosable bags × 25 pcs.) | Eppendorf Biopur [®] 100 pcs. (individually wrapped) |
|---------------------------------|--------------|---|--|---|
| 0.1 mL | □ white | 0030 089.405 | 0030 089.766 | 0030 089.618 |
| 0.2 mL | ■ light blue | 0030 089.413 | 0030 089.774 | 0030 089.626 |
| 0.5 mL | ■ violet | 0030 089.421 | 0030 089.782 | 0030 089.634 |
| 1 mL | ■ yellow | 0030 089.430 | 0030 089.790 | 0030 089.642 |
| 2.5 mL | ■ green | 0030 089.448 | 0030 089.804 | 0030 089.650 |
| 5 mL | ■ blue | 0030 089.456 | 0030 089.812 | 0030 089.669 |
| 10 mL | ■ orange | 0030 089.464 | 0030 089.820 | 0030 089.677 |
| 25 mL* | ■ red | 0030 089.472 | 0030 089.839 | 0030 089.685 |
| 50 mL* | ■ light gray | 0030 089.480 | 0030 089.847 | 0030 089.693 |

* 4 boxes of 25 pcs. each; each box containing one adapter.

| Accessories | Color coding | Eppendorf Quality [™] | Eppendorf Biopur [®] |
|---|--------------|--------------------------------|-------------------------------|
| 25 mL adapter (1 pc.) | ■ red | 0030 089.715 | |
| 50 mL adapter (1 pc.) | ■ light gray | 0030 089.723 | |
| 25 mL adapter (7 pcs.) | ■ red | | 0030 089.731 |
| 50 mL adapter (7 pcs.) | ■ light gray | | 0030 089.740 |
| Combitip Rack (for 8 Combitips advanced [®] , 0.1 mL–10 mL) | | 0030 089.758 | |
| Combitips advanced [®] | | 0030 089.936 | |
| Assortment pack (1 Combitip of each size, incl. adapters) | | | |

Eppendorf Varitips[®] P and S for Varipette[®]

| Description | | Order no. |
|---|---|--------------|
| Eppendorf Varitips[®] P | | |
| Starter Kit (100 Maxitips, 10 dispensing parts, 10 valves) | for aspirating from smaller vessels | 0030 050.525 |
| Eppendorf Varitips[®] P (set of 100) | | 0030 048.130 |
| Eppendorf Varitips[®] S—consisting of: | | |
| Barrels with pistons (for Eppendorf Varitips [®] S, set of 30) | for aspirating from Narrow-neck vessels and volumetric flasks | 0030 050.533 |
| Graduated Maxitip (for Eppendorf Varitips [®] S, set of 200) | | 0030 050.568 |
| Valves (for Eppendorf Varitips [®] S, set of 100) | | 0030 050.541 |

epMotion® Automated Pipetting System

| Description | Order no. |
|---|--------------|
| epT.I.P.S.® Motion pipette tips | |
| Automatic pipette tips in individual racks for use on the epMotion®. The tip type and size is automatically recognized on the device. 96 epT.I.P.S.®/ rack, 10 racks per set. The refill racks from the reload products can be placed in an autoclavable TipHolder adapter. Three purity levels: Eppendorf Quality™, PCR clean and sterile. Pipette tips available with or without a filter. Testing of production batches (certified). | |
| Without filter | |
| 50 µL, sterile, free of pyrogens, volume range 1–50 µL, 10 × 96 tips in racks | 0030 015.207 |
| 300 µL, sterile, free of pyrogens, volume range 20–300 µL, 10 × 96 tips in racks | 0030 015.223 |
| 1,000 µL, sterile, free of pyrogens, volume range 40–1,000 µL, 10 × 96 tips in racks | 0030 015.240 |
| 50 µL, Eppendorf Quality™, volume range 1–50 µL, 10 × 96 tips in racks | 0030 014.405 |
| 300 µL, Eppendorf Quality™, volume range 20–300 µL, 10 × 96 tips in racks | 0030 014.448 |
| 1,000 µL, Eppendorf Quality™, volume range 40–1,000 µL, 10 × 96 tips in racks | 0030 014.480 |
| 50 µL, Eppendorf Quality™, Reloads, volume range 1–50 µL, 24 × 96 tips | 0030 014.421 |
| 300 µL, Eppendorf Quality™, Reloads, volume range 20–300 µL, 24 × 96 tips | 0030 014.464 |
| 1,000 µL, Eppendorf Quality™, Reloads, volume range 40–1,000 µL, 24 × 96 tips | 0030 014.502 |
| SafeRacks, 50 µL, Eppendorf Quality™, volume range 1–50 µL, 10 × 96 tips | 0030 014.600 |
| SafeRacks, 300 µL, Eppendorf Quality™, volume range 20–300 µL, 10 × 96 tips | 0030 014.626 |
| SafeRacks, 1,000 µL, Eppendorf Quality™, volume range 40–1,000 µL, 10 × 96 tips | 0030 014.642 |
| With filter | |
| 50 µL, PCR clean, sterile, free of pyrogens, volume range 1–50 µL, 10 × 96 tips in racks | 0030 015.215 |
| 300 µL, PCR clean, sterile, free of pyrogens, volume range 20–300 µL, 10 × 96 tips in racks | 0030 015.231 |
| 1,000 µL, PCR clean, sterile, free of pyrogens, volume range 40–1,000 µL, 10 × 96 tips in racks | 0030 015.258 |
| 50 µL, PCR clean, volume range 1–50 µL, 10 × 96 tips in Racks | 0030 014.413 |
| 300 µL, PCR clean, volume range 20–300 µL, 10 × 96 tips in Racks | 0030 014.456 |
| 1,000 µL, PCR clean, volume range 40–1,000 µL, 10 × 96 tips in Racks | 0030 014.499 |
| 50 µL, PCR clean, Reloads, volume range 1–50 µL, 24 × 96 tips | 0030 014.430 |
| 50 µL, PCR clean, sterile, free of pyrogens, Reloads, volume range 1–50 µL, 24 × 96 tips | 0030 014.537 |
| 300 µL, PCR clean, Reloads, volume range 20–300 µL, 24 × 96 tips | 0030 014.472 |
| 300 µL, PCR clean, sterile, free of pyrogens, Reloads, volume range 20–300 µL, 24 × 96 tips | 0030 014.529 |
| 1,000 µL, PCR clean, Reloads, volume range 40–1,000 µL, 24 × 96 tips | 0030 014.510 |
| SafeRacks, 50 µL, PCR clean, volume range 1–50 µL, 10 × 96 tips | 0030 014.618 |
| SafeRacks, 300 µL, PCR clean, volume range 20–300 µL, 10 × 96 tips | 0030 014.634 |
| SafeRacks, 1,000 µL, PCR clean, volume range 40–1,000 µL, 10 × 96 tips | 0030 014.650 |
| Reservoir rack modules | |
| Are inserted in a Reservoir rack. They can be temperature controlled with a thermal module for heating and cooling. | |
| PCR 0.2 mL, for 8 × 0.2 mL PCR tubes | 5075 799.049 |
| Eppendorf Tubes® 5.0 mL, for 4 × 5 mL tubes | 5075 799.340 |
| Eppendorf Safe-Lock, for 4 × 0.5/1.5/2 mL tubes | 5075 799.081 |
| Tubes Ø 12 mm, for 4 × Ø 12 mm tubes | 5075 799.103 |
| Tubes Ø 16 mm, for 4 × Ø 16 mm tubes | 5075 799.120 |
| 15 mL conical tubes, for 4 × Ø 17 mm tubes | 5075 799.162 |
| 50 mL conical tubes, for 2 × Ø 29 mm tubes | 5075 799.189 |
| 30 mL reservoir, for use with reservoir rack, 5 × 10 large volume reservoirs, PCR clean | 0030 126.505 |
| 100 mL reservoir, for use with reservoir rack, 5 × 10 large volume reservoirs, PCR clean | 0030 126.513 |
| 400 mL reservoir, also for use in the epMotion® VAC, 10 pcs./set, made of PP | 5075 751.364 |

Eppendorf Pipette/epT.I.P.S.® Combinations

| Eppendorf Research® plus | | 0.1 µL– 10 µL (S) | 0.1 µL– 10 µL (M) | 0.1 µL– 20 µL | 0.5 µL– 20 µL L | 2 µL– 20 µL | 2 µL– 100 µL | 2 µL– 200 µL | 20 µL– 300 µL |
|---|---------------|----------------------|----------------------|------------------|--------------------|----------------|-----------------|-----------------|------------------|
| | | ■ dark gray | ■ medium gray | ■ light gray | ■ yellow | ■ yellow | ■ yellow | ■ yellow | ■ orange |
| ep T.I.P.S.® | | ✓ | | ✓ | ✓ | ✓ | ✓ | ✓ | ✓ |
| ep T.I.P.S.® LoRetention | | ✓ | | | ✓ | | | ✓ | |
| ep Dualfilter T.I.P.S.® | | ✓ | ✓ | | ✓ | ✓ | ✓ | ✓ | ✓ |
| ep Dualfilter T.I.P.S.® LoRetention | | ✓ | | | ✓ | | ✓ | | ✓ |
| ep Dualfilter T.I.P.S.® SealMax | | ✓ | | ✓ | | | ✓ | ✓ | ✓ |
| Fixed volume | | | | | | | | | |
| 10 µL | ■ medium gray | ✓ | ✓ | ✓ | ✓ | | | | |
| 20 µL | ■ light gray | | | | ✓ | | | | |
| 10 µL | ■ yellow | | | | | ✓ | ✓ | ✓ | |
| 20 µL | ■ yellow | | | | | ✓ | ✓ | ✓ | ✓ |
| 25 µL, 50 µL, 100 µL | ■ yellow | | | | | | ✓ | ✓ | ✓ |
| 200 µL | ■ yellow | | | | | | | ✓ | ✓ |
| 200 µL, 250 µL, 500 µL, 1,000 µL | ■ blue | | | | | | | | |
| Adjustable volume | | | | | | | | | |
| 0.1 µL–2.5 µL | ■ dark gray | ✓ | ✓ | ✓ | | | | | |
| 0.5 µL–10 µL | ■ medium gray | ✓ | ✓ | ✓ | ✓ | | | | |
| 2 µL–20 µL | ■ light gray | ✓ | ✓ | ✓ | ✓ | | | | |
| 2 µ–20 µL | ■ yellow | | | | | ✓ | ✓ | ✓ | ✓ |
| 10 µL–100 µL | ■ yellow | | | | | ✓ | ✓ | ✓ | ✓ |
| 20 µL–200 µL | ■ yellow | | | | | ✓ | ✓ | ✓ | ✓ |
| 30 µL–300 µL | ■ orange | | | | | ✓ | ✓ | ✓ | ✓ |
| 100 µL–1,000 µL | ■ blue | | | | | | | | |
| 0.5 mL–5 mL | ■ violet | | | | | | | | |
| 1 mL–10 mL | ■ turquoise | | | | | | | | |
| | | | | | | | | | |
| Eppendorf Xplorer® Eppendorf Xplorer® plus | | 0.1 µL– 10 µL (S) | 0.1 µL– 10 µL (M) | 0.1 µL– 20 µL | 0.5 µL– 20 µL L | 2 µL– 20 µL | 2 µL– 100 µL | 2 µL– 200 µL | 20 µL– 300 µL |
| | | ■ dark gray | ■ medium gray | ■ light gray | ■ yellow | ■ yellow | ■ yellow | ■ yellow | ■ orange |
| ep T.I.P.S.® | | ✓ | | ✓ | ✓ | ✓ | ✓ | ✓ | ✓ |
| ep T.I.P.S.® LoRetention | | ✓ | | | ✓ | | | ✓ | |
| ep Dualfilter T.I.P.S.® | | ✓ | ✓ | | ✓ | ✓ | ✓ | ✓ | ✓ |
| ep Dualfilter T.I.P.S.® LoRetention | | ✓ | | | ✓ | | ✓ | | ✓ |
| ep Dualfilter T.I.P.S.® SealMax | | ✓ | | ✓ | | | ✓ | ✓ | ✓ |
| 0.5 µL–10 µL | ■ medium gray | ✓ | ✓ | ✓ | ✓ | | | | |
| 5 µL–100 µL | ■ yellow | | | | | ✓ | ✓ | ✓ | ✓ |
| 15 µL–300 µL | ■ orange | | | | | ✓ | ✓ | ✓ | ✓ |
| 50 µL–1000 µL | ■ blue | | | | | | | | |
| 50 µL–1200 µL | ■ green | | | | | | | | |
| 0.25 mL–5 mL | ■ violet | | | | | | | | |
| 0.5 mL–10 mL | ■ turquoise | | | | | | | | |

✓: Compatible, ✓: Limited volume

| 50 µL– 1,000 µL ■ blue | 50 µL– 1,250 µL ■ green | 50 µL– 1,250 µL L ■ dark green | 0.25 mL– 2.5 mL ■ red | 0.1 mL– 5 mL ■ violet | 0.2 mL– 5 mL L ■ violet | 0.5 mL– 10 mL ■ turquoise | 0.5 mL– 10 mL L ■ turquoise | GELoader Microloader ■ light gray |
|------------------------------|-------------------------------|--------------------------------------|-----------------------------|-----------------------------|-------------------------------|---------------------------------|-----------------------------------|---|
| ✓ | ✓ | ✓ | ✓ | ✓ | ✓ | ✓ | ✓ | ✓ |
| ✓ | | | | | | | | |
| ✓ | | ✓ | | ✓ | ✓ | ✓ | ✓ | |
| ✓ | | | | | | | | |
| ✓ | | | | | | | | |
| | | | | | | | | |
| | | | | | | | | ✓ |
| | | | | | | | | ✓ |
| | | | | | | | | |
| | | | | | | | | |
| | | | | | | | | |
| ✓ | ✓ | ✓ | | | | | | |
| | | | | | | | | |
| | | | | | | | | ✓ |
| | | | | | | | | ✓ |
| | | | | | | | | |
| | | | | | | | | |
| ✓ | ✓ | ✓ | | | | | | |
| | | | | ✓ | ✓ | | | |
| | | | | | | ✓ | ✓ | |
| | | | | | | | | |
| | | | | | | | | |
| | | | | | | | | |
| | | | | | | | | |
| | | | | | | | | |
| | | | | | | | | |
| ✓ | ✓ | ✓ | | | | | | ✓ |
| | ✓ | ✓ | | | | | | |
| | | | | ✓ | ✓ | | | |
| | | | | | | ✓ | ✓ | |

Eppendorf Pipette/epT.I.P.S.® Combinations

| Eppendorf Reference® 2 | 0.1 µL– 10 µL (S) | 0.1 µL– 10 µL (M) | 0.1 µL– 20 µL | 0.5 µL– 20 µL L | 2 µL– 20 µL | 2 µL– 100 µL | 2 µL– 200 µL | 20 µL– 300 µL |
|---|----------------------|----------------------|------------------|--------------------|----------------|-----------------|-----------------|------------------|
| | ■ dark gray | ■ medium gray | ■ light gray | ■ yellow | ■ yellow | ■ yellow | ■ orange | |
| ep T.I.P.S.® | ✓ | | ✓ | ✓ | ✓ | ✓ | ✓ | ✓ |
| ep T.I.P.S.® LoRetention | ✓ | | | ✓ | | | ✓ | |
| ep Dualfilter T.I.P.S.® | ✓ | ✓ | | ✓ | ✓ | ✓ | ✓ | ✓ |
| ep Dualfilter T.I.P.S.® LoRetention | ✓ | | | ✓ | | ✓ | | ✓ |
| ep Dualfilter T.I.P.S.® SealMax | ✓ | | ✓ | | | ✓ | ✓ | ✓ |
| Fixed volume | | | | | | | | |
| 1 µL, 2 µL ■ dark gray | ✓ | ✓ | ✓ | | | | | |
| 5 µL, 10 µL ■ medium gray | ✓ | ✓ | ✓ | ✓ | | | | |
| 20 µL ■ light gray | | | ✓ | ✓ | | | | |
| 10 µL ■ yellow | | | | | ✓ | ✓ | ✓ | |
| 20 µL ■ yellow | | | | | ✓ | ✓ | ✓ | ✓ |
| 25 µL, 50 µL, 100 µL ■ yellow | | | | | | ✓ | ✓ | ✓ |
| 200 µL ■ yellow | | | | | | | ✓ | ✓ |
| 200 µL, 250 µL, 500 µL, 1,000 µL ■ blue | | | | | | | | |
| 2 mL, 2.5 mL ■ red | | | | | | | | |
| Adjustable volume | | | | | | | | |
| 0.1 µL–2.5 µL ■ dark gray | ✓ | ✓ | ✓ | | | | | |
| 0.5 µL–10 µL ■ medium gray | ✓ | ✓ | ✓ | ✓ | | | | |
| 2 µL–20 µL ■ light gray | ✓ | ✓ | ✓ | ✓ | | | | |
| 2 µL–20 µL ■ yellow | | | | | ✓ | ✓ | ✓ | ✓ |
| 10 µL–100 µL ■ yellow | | | | | ✓ | ✓ | ✓ | ✓ |
| 20 µL–200 µL ■ yellow | | | | | ✓ | ✓ | ✓ | ✓ |
| 30 µL–300 µL ■ orange | | | | | ✓ | ✓ | ✓ | ✓ |
| 100 µL–1,000 µL ■ blue | | | | | | | | |
| 0.25 mL–2.5 mL ■ red | | | | | | | | |
| 0.5 mL–5 mL ■ violet | | | | | | | | |
| 1 mL–10 mL ■ turquoise | | | | | | | | |

✓ : Compatible, ✓ : Limited volume

| 50 µL– 1,000 µL ■ blue | 50 µL– 1,250 µL ■ green | 50 µL– 1,250 µL L ■ dark green | 0.25 mL– 2.5 mL ■ red | 0.1 mL– 5 mL ■ violet | 0.2 mL– 5 mL L ■ violet | 0.5 mL– 10 mL ■ turquoise | 0.5 mL– 10 mL L ■ turquoise | GEloader Microloader ■ light gray |
|------------------------------|-------------------------------|--------------------------------------|-----------------------------|-----------------------------|-------------------------------|---------------------------------|-----------------------------------|---|
| ✓ | ✓ | ✓ | ✓ | ✓ | ✓ | ✓ | ✓ | ✓ |
| ✓ | | | | | | | | |
| ✓ | | ✓ | | ✓ | ✓ | ✓ | ✓ | |
| ✓ | | | | | | | | |
| ✓ | | | | | | | | |
| | | | | | | | | |
| | | | | | | | | ✓ |
| | | | | | | | | ✓ |
| | | | | | | | | |
| | | | | | | | | |
| | | | | | | | | |
| ✓ | ✓ | ✓ | | | | | | |
| | | | ✓ | | | | | |
| | | | | | | | | ✓ |
| | | | | | | | | ✓ |
| | | | | | | | | |
| | | | | | | | | |
| | | | | | | | | |
| ✓ | ✓ | ✓ | | | | | | |
| | | | ✓ | | | | | |
| | | | | ✓ | ✓ | | | |
| | | | | | | ✓ | ✓ | |



Make Your Lab a Better Place

The Eppendorf Liquid Handling Instruments Portfolio

As the inventor of the microliter system, we at Eppendorf have more than fifty years of experience in precise manual and automatic pipetting and transferring of the smallest quantities of liquids. The standards defined by Eppendorf quality assurance (for example, precision and accuracy guidelines) clearly surpass the standards required to obtain these results, and guarantee the reproducibility that our customers have come to depend on.

Our »PhysioCare Concept« incorporates growing ergonomic demands into all that we do. The use of our liquid handling products has been proven to reduce physical strain at the workstation to a minimum. This has been perfected by the ergonomic design of our products and flawless fit of device to consumable (e.g., measurably reduced ejecting force when pipetting).

For more information visit www.eppendorf.com



Your local distributor: www.eppendorf.com/contact

Eppendorf AG · 22331 Hamburg · Germany
eppendorf@eppendorf.com · www.eppendorf.com

www.eppendorf.com