



# **Total DNA, RNA, and protein isolation**

## **User manual**






NucleoSpin® TriPrep

February 2012/Rev.04

# Total DNA, RNA, and protein isolation

## Protocol-at-a-glance (Rev.04)

### NucleoSpin® TriPrep

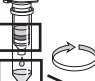




|  |   |  |
|--|---|--|
| <b>1 Homogenize sample</b>                     |  | Up to 30 mg  |
| <b>2 Lyse sample</b>                           |  | 350 µL RP1<br>3.5 µL β-mercaptoethanol<br>(or comparable reducing agent) |
| <b>3 Filtrate lysate</b>                       |  | 1 min, 11,000 x g  |
| <b>4 Adjust DNA and RNA binding conditions</b> |  | 350 µL ethanol (70 %)  |
| <b>5 Bind DNA and RNA</b>                      |  | Load sample<br>30 s, 11,000 x g  |


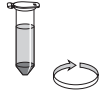
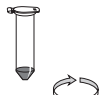
#### DNA and RNA Purification

(both bound to the silica membrane)

#### Protein Purification

(protein in the column flow-through)

|  |   |  |
|--|---|--|
| <b>6 Wash silica membrane</b>          |    | 1 <sup>st</sup> and 2 <sup>nd</sup> wash each:<br>500 µL <i>DNA Wash</i><br>1 min, 11,000 x g  |
| <b>7 Dry membrane</b>                  |   | RT,<br>3 min (with open lid)   |
| <b>8 Elute DNA</b>                     |   | 100 µL <i>DNA Elute</i><br>Incubate 1 min<br>1 min, 11,000 x g   |
| <b>9 Digest residual DNA</b>           |  | 95 µL DNase reaction mixture<br>RT, 15 min   |
| <b>10 Wash and dry silica membrane</b> |  | 1 <sup>st</sup> wash 200 µL RA2<br>2 <sup>nd</sup> wash 600 µL RA3<br>30 s, 11,000 x g<br>3 <sup>rd</sup> wash 250 µL RA3<br>2 min, 11,000 x g |
| <b>11 Elute highly pure RNA</b>        |  | 60 µL H <sub>2</sub> O (RNase-free)<br>1 min, 11,000 x g   |

|                                  |   |   |
|----------------------------------|---|---|
| <b>12 Precipitate protein</b>    |    | 10–700 µL flow-through<br>1 vol PP<br>RT, 10 min<br>5 min, 11,000 x g |
| <b>13 Wash protein pellet</b>    |   | 500 µL ethanol (50 %)<br>1 min, 11,000 x g                            |
| <b>14 Dry protein pellet</b>     |   | RT, 5–10 min  |
| <b>15 Prepare protein sample</b> |  | 20–100 µL PSB-TCEP<br>3 min, 95–98 °C<br>1 min, 11,000 x g            |

## Table of contents

|     |  |    |
|-----|--|----|
| 1   | Components   | 4  |
| 1.1 | Kit contents   | 4  |
| 1.2 | Reagents, consumables, and equipment to be supplied by user  | 5  |
| 1.3 | About this user manual   | 5  |
| 2   | Product description  | 6  |
| 2.1 | The basic principle  | 6  |
| 2.2 | Kit specifications   | 7  |
| 2.3 | Handling, preparation, and storage of starting materials   | 12 |
| 2.4 | Guideline for appropriate sample amount, precipitation volume, and resolubilization volume for protein isolation | 14 |
| 2.5 | Elution procedures for DNA   | 15 |
| 2.6 | Elution procedures for RNA   | 15 |
| 3   | Storage conditions and preparation of working solutions  | 16 |
| 4   | Safety instructions – risk and safety phrases  | 19 |
| 4.1 | Risk and safety phrases  | 19 |
| 4.2 | GHS classification   | 20 |
| 5   | Protocols  | 23 |
| 5.1 | DNA, RNA, and protein purification from cultured cells and tissue  | 23 |
| 5.2 | Total RNA preparation from RNA <i>later</i> <sup>®</sup> treated samples   | 30 |
| 5.3 | rDNase digestion in solution   | 31 |
| 6   | Appendix   | 33 |
| 6.1 | Protein quantification   | 33 |
| 6.2 | Troubleshooting  | 40 |
| 6.3 | References   | 45 |
| 6.4 | Ordering information   | 46 |
| 6.5 | Product use restriction/warranty   | 47 |

# 1 Components

## 1.1 Kit contents

| REF   | NucleoSpin® TriPrep   |                       |                         |
|---|-----------------------|-----------------------|-------------------------|
|   | 10 preps<br>740966.10 | 50 preps<br>740966.50 | 250 preps<br>740966.250 |
| Lysis Buffer RP1  | 10 mL                 | 25 mL                 | 125 mL                  |
| Buffer <i>DNA Wash</i><br>(Concentrate)*                                    | 4 mL                  | 12 mL                 | 3 x 22.5 mL             |
| Buffer <i>DNA Elute</i>   | 1.2 mL                | 6 mL                  | 3 x 12 mL               |
| Wash Buffer RA2   | 15 mL                 | 15 mL                 | 80 mL                   |
| Wash Buffer RA3<br>(Concentrate)*   | 5 mL                  | 12.5 mL               | 3 x 25 mL               |
| RNase-free H <sub>2</sub> O   | 5 mL                  | 15 mL                 | 65 mL                   |
| Protein Precipitator PP   | 9 mL                  | 45 mL                 | 225 mL                  |
| Protein Solving Buffer PSB<br>(without reducing agent)                      | 2 x 1 mL              | 7.5 mL                | 5 x 7.5 mL              |
| Reducing Agent TCEP   | 2 x 14 mg             | 107 mg                | 5 x 107 mg              |
| Reaction Buffer for rDNase  | 3 mL                  | 7 mL                  | 30 mL                   |
| rDNase, RNase-free<br>(lyophilized)*  | 1 vial (size C)       | 1 vial (size D)       | 5 vials (size D)        |
| NucleoSpin® Filters<br>(violet rings)                                       | 10                    | 50                    | 250                     |
| NucleoSpin® TriPrep<br>Columns (light blue rings,<br>plus Collection Tubes) | 10                    | 50                    | 250                     |
| Collection Tubes (2 mL)   | 30                    | 150                   | 750                     |
| Collection Tubes (1.5 mL)   | 20                    | 100                   | 500                     |
| User manual   | 1                     | 1                     | 1                       |

\* For preparation of working solutions and storage conditions see section 3.

## 1.2 Reagents, consumables, and equipment to be supplied by user

### Reagents

- 50 % ethanol (to prepare Buffer *DNA Wash*)
- 70 % ethanol (to adjust RNA binding conditions)
- 96–100 % ethanol (to prepare Wash Buffer RA3)
- Reducing agent ( $\beta$ -mercaptoethanol, DTT (dithiothreitol), or TCEP (BisTris (Bis-(2-hydroxyethyl)-imino-tris(hydroxymethyl)-methane)) to supplement lysis buffer

### Consumables

- 1.5 mL microcentrifuge tubes (for sample lysis and DNA elution)
- Disposable RNase-free pipette tips

### Equipment

- Manual pipettors
- Centrifuge for microcentrifuge tubes
- Vortex mixer
- Thermal heating block
- Equipment for sample disruption and homogenization
- Personal protection equipment (lab coat, gloves, goggles)

Additional material is furthermore needed for protein quantification, see section 6.1.

## 1.3 About this user manual

It is strongly recommended reading the detailed protocol sections of this user manual if the **NucleoSpin® TriPrep** kit is used for the first time. Experienced users, however, may refer to the Protocol-at-a-glance instead. The Protocol-at-a-glance is designed to be used only as a supplemental tool for quick referencing while performing the purification procedure.

All technical literature is available on the internet at **[www.mn-net.com](http://www.mn-net.com)**.

Please contact Technical Service regarding information about changes of the current user manual compared to previous revisions.

## 2 Product description

### 2.1 The basic principle

#### Introduction

Studies of gene expression at the level of transcription and translation by quantification of RNA and protein combined with verification of genomic sequence (e.g., transgene integration site) are often hampered by the small sample size and the necessity of different – often incompatible – techniques for DNA, RNA, and protein isolation. Samples may comprise biopsies, tumors, tissues, transgene organisms, and others. The **NucleoSpin® TriPrep** kit enables isolation of DNA, RNA, and protein from diverse sample types. DNA, RNA, and protein are isolated without splitting the sample prior to extraction. Thus, DNA, RNA, and protein are obtained from one and the same sample and not from three similar portions of one sample. This is especially valuable for unique, small, and precious samples. DNA and RNA are eluted separately from the NucleoSpin® TriPrep Column, with a low salt buffer and water, respectively. Isolated DNA and RNA are suitable for all common downstream applications. Isolated protein is suitable for SDS-PAGE, Western Blot analysis, and quantification.

#### DNA, RNA, and protein isolation

One of the most important aspects in the isolation of DNA, RNA, and protein is to prevent their degradation during the isolation procedure. With the **NucleoSpin® TriPrep** method, cells are lysed by incubation in a solution containing large amounts of chaotropic ions. This lysis buffer immediately inactivates virtually all enzymes (e.g., DNases, RNases, proteases, and phosphatases) which are present in almost all biological materials. The buffer dissolves even hardly soluble protein, creates appropriate binding conditions which favor adsorption of DNA and RNA to the silica membrane, and enables protein to pass the specially treated NucleoSpin® TriPrep Column virtually quantitatively. Expensive and harmful proteinase inhibitors or inhibitor cocktails are not necessary due to the denaturing properties of the lysis buffer. After two special washing steps, DNA is eluted with a low salt buffer (*DNA Elute*) which selectively elutes DNA and keeps RNA quantitatively on the column. Eluted DNA is immediately ready for downstream applications without further purification. DNA elution prior to RNA elution does neither compromise RNA quantity nor quality. RNA isolated with the **NucleoSpin® TriPrep** kit is of identical quality as RNA isolated with the well proven NucleoSpin® RNA II kit. Residual DNA still bound to the silica membrane is removed by an rDNase solution which is directly applied onto the silica membrane during the preparation (RNase-free rDNase is supplied with the kit). Simple washing steps with two different buffers remove salts, metabolites, and macromolecular cellular components. Pure RNA is finally eluted under low ionic strength conditions with RNase-free water (supplied).

Protein is isolated from the column flow-through. Protein is precipitated in denatured form with a special buffer (Protein Precipitator PP) which effectively precipitates protein. After a washing step the protein pellet is dissolved in Protein Solving Buffer (PSB) containing the odourless reducing agent TCEP. The protein can thus readily be applied to SDS-PAGE analysis. **The kit is not recommended for isolation of native proteins.**

DNA, RNA, and protein preparation using **NucleoSpin® TriPrep** kits can be performed at room temperature. The DNA and RNA eluates, however, should be treated with care because RNA is very sensitive to trace contaminations of RNases, often found on general lab ware, fingerprints, and dust. DNA can be stored at 4 °C for short term and at -20 °C for long term storage. To ensure RNA stability keep RNA frozen at -20 °C for short-term or -70 °C for long-term storage. Recovered protein dissolved in Protein Solving Buffer is unproblematic concerning stability.

## 2.2 Kit specifications

- **NucleoSpin® TriPrep** kits are recommended for the isolation of total DNA, RNA, and protein from cultured cells, tissue, and other biological samples.

**Table 1: Kit specifications at a glance**

| Parameter                             | NucleoSpin® TriPrep   |                               |                      |
|---------------------------------------|---|-------------------------------|----------------------|
| Format                                | Mini spin column  |                               |                      |
| Sample material                       | < 5 x 10 <sup>6</sup> cultured cells*,<br>< 30 mg human / animal tissue*,<br>< 100 mg plant tissue* |                               |                      |
|                                       | <b>Total RNA</b>  | <b>Total DNA</b>              | <b>Total protein</b> |
| Fragment size                         | > 200 b   | < 30 kbp                      | 15–300 kDa           |
| Typical yield                         | < 70 µg   | < 6 µg                        | < 1200 µg            |
| A <sub>260</sub> /A <sub>280</sub>    | 1.9–2.1   | 1.7–1.9                       | –                    |
| Typical RIN<br>(RNA integrity number) | > 9   | –                             | –                    |
| Elution volume<br>(RNA and DNA)       | 40–120 µL   | 100 µL                        | 10–100 µL            |
| Resolubilization volume<br>(protein)  |   |                               |                      |
| Preparation time                      | 30 min/6 preps  | 45 min/6 preps<br>(RNA + DNA) | 35 min/6 preps       |
| Binding capacity                      | 200 µg  | 10 µg**                       | –                    |

\* For samples larger than approx. one million cells or 5 mg tissue, DNA yield does not increase linearly with the sample amount. DNA yield may decrease utilizing five million cells, 30 mg tissue, or more. Typically however, DNA yield is still sufficient for PCR analysis.

\*\*Binding capacity of DNA ≤ 10 µg, strongly depending on RNA amount bound to the membrane

### DNA characteristics

- Isolated DNA is of high molecular weight and typically exceeds 20 kb.
- DNA is commonly stable, even at 37°C for 2 h with or without addition of a typical restriction enzyme buffer, showing the absence of DNases.
- DNA is digestible with restriction enzymes and suitable for PCR.

### RNA characteristics


- The **NucleoSpin® TriPrep** kit allows purification of pure RNA with an  $A_{260}/A_{280}$  ratio generally exceeding 1.9 (measured in TE buffer, pH 7.5).
- The isolated RNA is ready to use for applications like reverse transcriptase PCR (RT-PCR), primer extension, or RNase protection assays.
- RNA of high integrity can be isolated with NucleoSpin® TriPrep kit. RIN (RNA Integrity Number) of RNA isolated from fresh high quality sample material (e.g., eukaryotic cells or fresh mouse liver) generally exceeds 9.0. However, RNA integrity strongly depends on the sample quality. RNA integrity was examined using the Agilent 2100 Bioanalyzer in conjunction with the RNA 6000 Nano or Pico assay.
- RNA prepared with **NucleoSpin® TriPrep** is generally free of residual DNA. However, minute traces of DNA may remain, if large amounts of material rich in nucleic acids are used. If the isolated RNA will be used as template in a RT-PCR-reaction, we recommend using lower quantities of sample material, depending on cell or tissue type (in the range of  $1 \times 10^6$  cells or 10 mg of tissue resulting in about 20 µg of RNA).

### Protein characteristics

- Small (17 kDa) to large (250 kDa) proteins, as well as glycoproteins, membrane proteins, lipoproteins, phosphorylated proteins, and structural proteins have been analyzed successfully.
- The isolated protein is ready to use for SDS-PAGE, Western Blot analysis and protein quantification with the Protein Quantification Assay (see ordering information).



## Typical yields of DNA, RNA, and protein

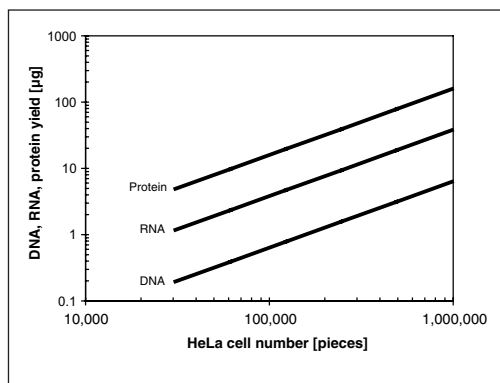


**Figure 1: DNA, RNA, and protein yield from different amounts of HeLa cells**

DNA, RNA, and protein were isolated from different amounts of HeLa cells. For cell numbers exceeding 250,000 only a fraction of the protein solution was used for precipitation and quantification; total protein yield was calculated.

DNA, RNA, and protein were isolated as described in Figure 1. Obtained correlation coefficients between HeLa sample amount and DNA, RNA, and protein yield are shown in Table 2.

**Table 2: Correlation between sample amount, nucleic acid and protein yield**

|         | $3 \times 10^4$ – $5 \times 10^5$ cells | $3 \times 10^4$ – $1 \times 10^6$ cells |
|---------|---|---|
| DNA     | > 0.99                                  | > 0.95                                  |
| RNA     | > 0.98                                  | > 0.98                                  |
| Protein | > 0.99                                  | > 0.99                                  |

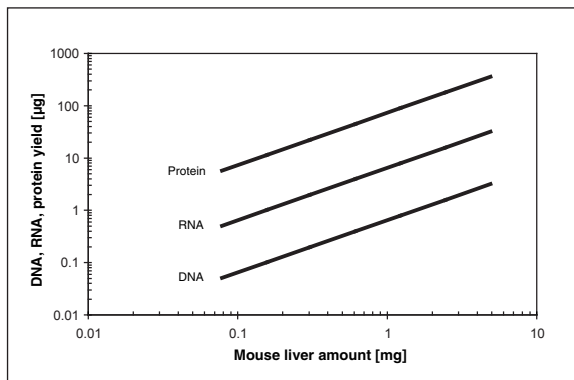

**Figure 2: DNA, RNA, and protein yield from different amounts of mouse liver tissue**

DNA, RNA, and protein were isolated from different amounts of mouse liver tissue. For tissue amounts exceeding 2.5 mg only a fraction of the protein solution was used for precipitation and quantification; total protein yield was calculated.

DNA, RNA, and protein were isolated as described in Figure 2. Obtained correlation coefficients between sample amount and DNA, RNA, and protein yield are shown in Table 3.

**Table 3: Correlation between sample amount, nucleic acid and protein yield**

|         | 0.08–1.25 mg mouse liver | 0.08–2.5 mg mouse liver | 0.08–5 mg mouse liver |
|---------|--------------------------|-------------------------|-----------------------|
| DNA     | > 0.99                   | > 0.95                  | > 0.67                |
| RNA     | > 0.98                   | > 0.98                  | > 0.98                |
| Protein | > 0.99                   | > 0.99                  | > 0.99                |

**Protein yield**

Protein yield depends on sample type, amount and quality as well as on homogenization efficiency. Further the utilized quantification method influences determined protein yield. The following values were determined with MACHEREY-NAGEL's Protein Quantification Assay (see ordering information) and shall serve as a guideline for expected protein yield. It is assumed that the complete sample amount is processed, that is the complete lysed sample – after ethanol addition – is loaded onto the column and the complete 700 µL flow through is subjected to protein precipitation. Note that in many cases precipitation of only a portion of the column flow through (e.g., 100 µL) is recommended and will yield enough protein in terms of absolute amount and concentration for SDS-PAGE and Western Blot analysis.

As a guideline for appropriate precipitation volumes see section 2.4.

**Table 4: Typical protein yield**

| <b>Sample type and amount</b>                           | <b>Protein yield</b> |
|---|----------------------|
| Cultured human cells (e.g., HeLa, approx. $10^6$ cells) | ~ 50–150 µg          |
| Plants (e.g., garden cress, approx. 100 mg)             | ~ 150–350 µg         |
| Animal tissue (e.g., pig liver, approx. 30 mg)          | ~ 500–1200 µg        |

## 2.3 Handling, preparation, and storage of starting materials

RNA is not protected against digestion until the sample material is flash frozen or disrupted in the presence of RNase inhibiting or denaturing agents. Therefore it is important that samples are flash frozen in liquid N<sub>2</sub> immediately and stored at -70 °C or processed as soon as possible. Samples can be stored in Lysis Buffer RP1 after disruption at -70 °C for up to one year, at +4 °C for up to 24 hours, or up to several hours at room temperature. Frozen samples are stable up to 6 months. Frozen samples in Buffer RP1 should be thawed slowly before starting with the isolation of total RNA.

*Wear gloves at all times during the preparation. Change gloves frequently.*

**Cultured animal cells** are collected by centrifugation and directly lysed by adding Buffer RP1 according to step 2 of the standard protocol (see section 5.1).

### Cell lysis of adherent growing cells in a culture dish:

Completely aspirate cell-culture medium and continue immediately with the addition of Lysis Buffer RP1 to the cell-culture dish. Avoid incomplete removal of the cell culture medium in order to allow full lysis activity of the lysis buffer.

### To trypsinize adherent growing cells:

Aspirate cell culture medium and add an equal amount of PBS in order to wash the cells. Aspirate PBS. Add 0.1–0.3 % trypsin in PBS and incubate for an appropriate time to detach the cells from the dish surface. After cell detachment, add medium, transfer cells to an appropriate tube (not supplied), and pellet by centrifugation for 5 min at 300 x g. Remove supernatant and continue with the addition of lysis buffer to the cell pellet.

**Human / animal and plant tissues** are often solid and must therefore be broken up mechanically as well as lysed. Depending on the disruption method, the viscosity of the lysed sample has to be reduced further for optimal results. It is essential for efficient RNA preparation that all RNA contained in the sample is released from the cells by disruption and that the viscosity of the sample is reduced by homogenization.

The most commonly used technique for disruption of animal tissues is grinding with a **pestle and mortar**. Grind the sample to a fine powder in presence of liquid N<sub>2</sub>. Take care that the sample does not thaw during or after grinding or weighing and add the frozen powder to an appropriate aliquot of Buffer RP1 containing β-mercaptoethanol and mix immediately. The broken-up tissue must then be homogenized with a **NucleoSpin® Filter** or by passing ≥5 times through a 0.9 mm syringe needle.

Thawing of undisrupted animal tissue should exclusively be done in presence of Buffer RP1 during simultaneous mechanical disruption, for example with a **rotor-stator homogenizer**. This ensures that RNA is not degraded by RNases before the preparation starts. The spinning rotor disrupts and simultaneously homogenizes the sample by mechanical shearing within seconds up to minutes (homogenization time

depends on sample). Take care to keep the rotor tip submerged in order to avoid excess foaming. To degenerate evolved foam, centrifuge 1 min at 400 x *g*. Select a suitably sized homogenizer (5–7 mm diameter rotors can be used for homogenization in microcentrifuge tubes).

**Bacteria and yeast** have to be incubated in lysozyme or lyticase/zymolase solutions, respectively. By this treatment, the robust cell walls of these organisms are digested or at least weakened, which is essential for effective cell lysis by Buffer RP1. For microorganisms with extremely resistant cell walls – like some Gram-positive bacterial strains – it may be necessary to optimize the conditions of the treatment with lytic enzymes or the cultivation conditions. After lysis, homogenization is achieved by using a **NucleoSpin® Filter** or the syringe-needle method.

## 2.4 Guideline for appropriate sample amount, precipitation volume, and resolubilization volume for protein isolation

The following Table 5 shall serve as a first guide for choosing appropriate amounts of sample material, precipitation volume, and resolubilization volume. Depending on sample type and downstream application (e.g., Coomassie or silver stain, sensitivity of antibody, detection system) appropriate volumes might deviate from the table below and have to be determined experimentally.

**Table 5: Guideline for appropriate sample amount**

| Amount of   | Cultivated cells<br>(e.g., HeLa) |                 |                 | Animal tissue<br>(e.g., liver) |        |        | Plant tissue<br>(e.g., garden<br>cress leaf) |        |        |
|---|----------------------------------|-----------------|-----------------|--------------------------------|--------|--------|--|--------|--------|
|   | 10 <sup>6</sup>                  | 10 <sup>5</sup> | 10 <sup>4</sup> | 30 mg                          | 3 mg   | 0.3 mg | 100 mg                                       | 10 mg  | 1 mg   |
| Lysis Buffer RP1<br>incl. reducing<br>agent                             | 350 µL                           |                 |                 |                                |        |        |  |        |        |
| Ethanol (70 %)  | 350 µL                           |                 |                 |                                |        |        |  |        |        |
| Column flow-<br>through to be<br>precipitated*                          | 35 µL                            | 350 µL          | 700 µL          | 35 µL                          | 350 µL | 700 µL | 35 µL  | 350 µL | 700 µL |
| Volume of Protein<br>Precipitator PP                                    | 35 µL                            | 350 µL          | 700 µL          | 35 µL                          | 350 µL | 700 µL | 35 µL  | 350 µL | 700 µL |
| Buffer PSB used<br>for protein pellet<br>solubilization                 | 100 µL                           | 100 µL          | 20 µL           | 100 µL                         | 100 µL | 20 µL  | 100 µL                                       | 100 µL | 20 µL  |
| Protein sample to<br>be analyzed by<br>SDS-PAGE with<br>Coomassie stain | 10 µL                            |                 |                 |                                |        |        |  |        |        |
| Protein sample to<br>be analyzed by<br>SDS-PAGE with<br>silver stain    | 1 µL                             |                 |                 |                                |        |        |  |        |        |
| Protein sample to<br>be analyzed by<br>Western Blot                     | 1–10 µL                          |                 |                 |                                |        |        |  |        |        |

\* Protein pellets with a diameter of up to approximately 1–2 mm in size are ideally suited for subsequent solubilization. Protein pellets exceeding volumes of approximately 10 µL should be avoided as large protein pellets are harder to dissolve than small pellets. To obtain small protein pellets, adapt the volume of column flow through in respect to the amount of sample material. Commonly small and even invisible protein pellets yield sufficient protein for SDS PAGE and Western Blot analysis.

## Solubilization of protein pellets and reduction of protein disulfide bonds

The NucleoSpin® TriPrep kit provides a protein sample buffer (Protein Solving Buffer PSB) and the Reducing Agent TCEP.

The Protein Solving Buffer PSB is similar in composition and function to the buffer commonly known as “Laemmli” buffer. For most applications, PSB may be substituted by “Laemmli” buffer. However, for applications with large protein pellets ( approx. 1 mm) PSB is recommended.

TCEP is a powerful, multi-purpose, and odourless reducing agent. It is non-volatile and unlike commonly used reducing agents like DTT and  $\beta$ -mercaptoethanol resistant to air oxidation. TCEP reduces disulfide bonds as effectively as dithiothreitol (DTT). TCEP reduces even most stable water-soluble alkyl disulfides selectively and completely over a wide pH range.

Solubilization of TCEP in PSB according to the instruction, results in a PSB-TCEP solution with a concentration of 50 mM TCEP (see section 6.1 for composition). This provides sufficient molar excess to reduce peptide and protein disulfide bonds effectively within a few minutes (in a range up to a protein concentration of approximately 1  $\mu$ g/ $\mu$ L).

## 2.5 Elution procedures for DNA

Elution of DNA is carried out under selective conditions to make sure that only DNA is eluted while RNA is still bound to the membrane. The DNA washing solution *DNA Wash* and the DNA elution buffer *DNA Elute* are finely tuned to achieve this. Therefore, the DNA elution volume should only be altered moderately, in the range of 60–150  $\mu$ L. Furthermore, the temperature of the *DNA Elute* solution shall not exceed 30 °C, otherwise RNA will partly elute with the *DNA Elute* solution. *DNA Elute* solution may stay for 1 min up to 15 min on the column, a 1–5 min incubation time is recommended. Eluted DNA is immediately ready for downstream applications without further purification.

## 2.6 Elution procedures for RNA

It is possible to adapt elution method and elution volume of water to the subsequent application of interest. In addition to the standard method described in the individual protocols (recovery rate about 70–90 %) there are several modifications possible.

- **High yield:** Perform two elution steps with the volume indicated in the individual protocol. About 90–100 % of bound nucleic acid will be eluted.
- **High yield and high concentration:** Elute with the standard elution volume and apply the eluate once more onto the column for reelution.

Eluted RNA should immediately be put on ice and always kept on ice for optimal stability and to avoid RNA degradation by almost omnipresent RNases (general lab ware, fingerprints, dust). For short-term storage freeze at -20 °C, for long-term storage freeze at -70 °C.

### 3 Storage conditions and preparation of working solutions

**Attention:** Buffers RP1 and RA2 contain chaotropic salt. Wear gloves and goggles!

**CAUTION:** Buffers RP1 and RA2 contain guanidinium thiocyanate which can form highly reactive compounds when combined with bleach (sodium hypochlorite). DO NOT add bleach or acidic solutions directly to the sample-preparation waste.

Storage conditions:

- Store lyophilized **rDNase (RNase-free)** at 4 °C on arrival (stable up to 1 year).
- Store **Reducing Agent TCEP** at 4 °C on arrival.
- All other kit components should be stored at room temperature (18–25 °C) and are stable up to one year. Storage at lower temperatures may cause precipitation of salts.

Before starting any **NucleoSpin® TriPrep** protocol prepare and consider the following:

- Check that 70% ethanol is available as additional solution to adjust binding conditions in the RP1-lysate.
- Check that 50% ethanol is available as additional solution to wash the protein pellet and to prepart Buffer *DNA Wash*.
- Check that 96–100% ethanol is available to prepare Wash Buffer RA3.
- The ***DNA Wash*** solution is delivered as a concentrate. To prepare the final *DNA Wash* solution add four volumes of 50% ethanol to the *DNA Wash* Concentrate as indicated in the table below. Mark the label of the bottle to indicate that the ethanol is added. Keep the bottle tightly close to avoid evaporation of ethanol.
- Due to its composition the ***DNA Elute*** solution (DNA elution buffer) does not inhibit DNases, that is *DNA Elute* does not contain substances (e.g., EDTA) to complex divalent cations. Therefore, be aware not to contaminate *DNA Elute* with DNases! Further, due to its composition, *DNA Elute* solution does not inhibit microbial growth. Therefore, make sure not to contaminate *DNA Elute* with any source of microbial contamination.
- **rDNase, RNase-free:** Add indicated volume of RNase-free H<sub>2</sub>O (see table below) to the rDNase vial and incubate for 1 min at room temperature. Gently swirl the vials to completely dissolve the rDNase. Be careful not to mix rDNase vigorously as rDNase is sensitive to mechanical agitation. Dispense into aliquots and store at -18 °C. The frozen working solution is stable for 6 months. Do not freeze/thaw the aliquots more than three times.



- **Wash Buffer RA3:** Add the indicated volume of 96–100 % ethanol (see table below) to Buffer RA3 Concentrate. Mark the label of the bottle to indicate that the ethanol is added. Store Wash Buffer RA3 at room temperature (18–25 °C) for up to one year. Keep the bottle tightly close to avoid evaporation of ethanol.
- **Protein Solving Buffer PSB and Reducing Agent TCEP:** For SDS-PAGE under reducing conditions (most common type of SDS-PAGE) transfer PSB from one vial to one vial of the Reducing Agent TCEP. Mix gently to avoid excessive foaming until the reducing agent is dissolved completely (this process will require several minutes)\*. Protein Solving Buffer containing Reducing Agent TCEP (PSB-TCEP) is stable for several days at room temperature (18–25 °C) and several months at 4 °C. For long term storage of PSB-TCEP keep at -20 °C.
- If SDS-PAGE under non-reducing conditions is intended consider the following:  
A: Omit addition of the Reducing Agent TCEP to Buffer PSB.  
B: Omit addition of  $\beta$ -mercaptoethanol (or other reducing agent) to Lysis Buffer RP1.
- If other reducing agents than TCEP are preferred (e.g., DTT,  $\beta$ -mercaptoethanol), appropriate amounts should be added to PSB. Please consider limited stability of DTT compared to TCEP.
- If PSB-TCEP is turbid, warm up PSB-TCEP to > 25 °C before use until solution is completely clear (i.e., all precipitate is dissolved completely). PSB-TCEP has a half-life of approximately 5 months if stored at 4 °C and approximately 7 months if stored at -20 °C.

---

\* For 50 and 250 prep kits: For better handling, PSB-TCEP may be transferred into the original PSB vial (with screw cap).

| <b>NucleoSpin® TriPrep</b>                 |  |  |   |
|--|--|--|---|
| <b>REF</b>                                 | <b>10 preps<br/>740966.10</b>                                | <b>50 preps<br/>740966.50</b>                                | <b>250 preps<br/>740966.250</b>   |
| Buffer<br><i>DNA Wash</i><br>(Concentrate) | 4 mL<br>Add 16 mL<br>50% ethanol                             | 12 mL<br>Add 48 mL<br>50% ethanol                            | 3 x 22.5 mL<br>Add 90 mL 50 %<br>ethanol to each vial                         |
| Wash<br>Buffer RA3<br>(Concentrate)        | 5 mL<br>Add 20 mL<br>96–100 % ethanol                        | 12.5 mL<br>Add 50 mL<br>96–100 % ethanol                     | 3 x 25 mL<br>Add 100 mL<br>96–100 % ethanol<br>to each vial                   |
| rDNase,<br>RNase-free<br>(lyophilized)     | 1 vial (size C)<br>Add 230 µL<br>RNase-free H <sub>2</sub> O | 1 vial (size D)<br>Add 540 µL<br>RNase-free H <sub>2</sub> O | 5 vials (size D)<br>Add 540 µL<br>RNase-free H <sub>2</sub> O<br>to each vial |
| Reducing<br>Agent TCEP*                    | 2 x 14 mg<br>Add 1 mL PSB each                               | 107 mg<br>Add 7.5 mL PSB                                     | 2 x 107 mg<br>Add 7.5 mL PSB each   |

## 4 Safety instructions – risk and safety phrases

The following components of the **NucleoSpin® TriPrep** kits contain hazardous contents. *Wear gloves and goggles and follow the safety instructions given in this section.*

### 4.1 Risk and safety phrases

| Component          | Hazard contents  | Hazard symbol     | Risk phrases           | Safety phrases |
|--------------------|--|-------------------|------------------------|----------------|
| Inhalt             | Gefahrstoff  | Gefahrstoffsymbol | R-Sätze                | S-Sätze        |
| RP1                | Guanidinium thiocyanate 30–60 %<br><i>Guanidiniumthiocyanat 30–60 %</i>  | ✘ Xn*             | R 20/21/22-32-52/53    | S 13-61        |
| RA2                | Guanidinium thiocyanate 30–60 % + ethanol 20-35 %<br><i>Guanidiniumthiocyanat 30–60 % + Ethanol 20–35 %</i>        | ✘ Xn*             | R 10-20/21/22-32-52/53 | S 13-16-61     |
| TCEP               | Tris(2-carboxyethyl) phosphine hydrochloride 70–100 %<br><i>Tris(2carboxy-ethyl)phoshine Hydrochlorid 70–100 %</i> | ✘ Xi**            | R 36/38                | S 26-37/39     |
| rDNase, RNase-free | rDNase, lyophilized<br><i>rDNase, lyophilisiert</i>  | ✘ Xn              | R 42/43                | S 22-24        |

#### Risk phrases

|            |   |
|------------|---|
| R 10       | Flammable.<br><i>Entzündlich.</i>   |
| R 20/21/22 | Harmful by inhalation, in contact with skin and if swallowed.<br><i>Gesundheitsschädlich beim Einatmen, Verschlucken und Berührung mit der Haut.</i>  |
| R 32       | Contact with acids liberates very toxic gas.<br><i>Entwickelt bei Berührung mit Säure sehr giftige Gase.</i>  |
| R 36/38    | Irritating to eyes and skin<br><i>Reizt die Augen und die Haut.</i>   |
| R 42/43    | May cause sensitization by inhalation and skin contact<br><i>Sensibilisierung durch Einatmen und Hautkontakt möglich.</i>   |
| R 52/53    | Harmful to aquatic organisms, may cause long-term adverse effects in the aquatic environment.<br><i>Schädlich für Wasserorganismen, kann in Gewässern längerfristig schädliche Wirkungen haben.</i> |

\* Hazard labeling not necessary if quantity per bottle below 125 g or mL (certificate of exemption according to 67/548/EEC Art. 25, 1999/45/EC Art. 12 and German GefStoffV § 20 (3) and TRGS 200 7.1). For further information see Material Safety Data Sheet.

\*\*Hazard labeling not necessary if quantity per bottle below 25 g or mL (certificate of exemption according to 67/548/EEC Art. 25, 1999/45/EC Art. 12 and German GefStoffV § 20 (3) and TRGS 200 7.1). For further information see Material Safety Data Sheet.





### Safety phrases

- S 13 Keep away from food, drink and animal foodstuffs.  
*Von Nahrungsmitteln, Getränken und Futtermitteln fernhalten.*
- S 16 Keep away from sources of ignition – No smoking.  
*Von Zündquellen fernhalten – Nicht rauchen!*
- S 22 Do not breathe dust.  
*Staub nicht einatmen.*
- S 24 Avoid contact with the skin.  
*Berührung mit der Haut vermeiden.*
- S 26 In case of contact with eyes, rinse immediately with plenty of water and seek medical advice.  
*Bei Berührung mit den Augen gründlich mit Wasser abspülen und Arzt konsultieren.*
- S 37/39 Wear suitable protective clothing and gloves.  
*Bei der Arbeit geeignete Schutzhandschuhe und Schutzkleidung tragen.*
- S 61 Avoid release to the environment. Refer to special instructions / safety data sheet.  
*Freisetzung in die Umwelt vermeiden. Besondere Anweisungen einholen / Sicherheitsdatenblätter zu Rate ziehen.*

## 4.2 GHS classification

Only harmful features do not need to be labeled with H and P phrases until 125 mL or 125 g.

*Mindergefährliche Eigenschaften müssen bis 125 mL oder 125 g nicht mit H- und P-Sätzen gekennzeichnet werden.*

| Component             | Hazard contents   | GHS symbol  | Hazard phrases            | Precaution phrases   |
|-----------------------|---|---|---------------------------|--|
| <i>Inhalt</i>         | <i>Gefahrstoff</i>  | <i>GHS Symbol</i>   | <i>H-Sätze</i>            | <i>P-Sätze</i>   |
| RP1                   | Guanidinium thiocyanate<br>30–60 %<br><i>Guanidiniumthiocyanat</i><br>30–60 %   |    | Warning<br><i>Achtung</i> | 302, 412,<br>EUH031<br><br>260, 273,<br>301+312, 330                               |
| RA2                   | Guanidinium thiocyanate<br>30–60 % + ethanol 20-<br>35 %<br><i>Guanidiniumthiocyanat</i><br>30–60 % + Ethanol 20–35 %               |  | Warning<br><i>Achtung</i> | 226, 302,<br>412,<br>EUH031<br><br>210, 233, 260,<br>273, 301+312,<br>330, 403+235 |
| TCEP                  | Tris(2-carboxyethyl)<br>phosphine hydrochloride<br>70–100 %<br><i>Tris(2carboxy-ethyl)phosphine</i><br><i>Hydrochlorid</i> 70–100 % |  | Warning<br><i>Achtung</i> | 315, 319<br><br>280, 302+352,<br>305+351+338,<br>332+313,<br>337+313               |
| rDNase,<br>RNase-free | rDNase, lyophilized<br><i>rDNase, lyophilisiert</i>   |  | Danger<br><i>Gefahr</i>   | 317, 334<br><br>261, 280,<br>302+352,<br>304+341,<br>333+313,<br>342+311, 363      |

### Hazard phrases

- H 226 Flammable liquid and vapour.  
*Flüssigkeit und Dampf entzündbar.*
- H 302 Harmful if swallowed.  
*Gesundheitsschädlich bei Verschlucken.*
- H 315 Causes skin irritation.  
*Verursacht Hautreizungen.*
- H 317 May cause an allergic skin reaction.  
*Kann allergische Hautreaktionen verursachen.*
- H 319 Causes serious eye irritation.  
*Verursacht schwere Augenreizung.*
- H 334 May cause allergy or asthma symptoms or breathing difficulties if inhaled.  
*Kann bei Einatmen Allergie, asthmaartige Symptome oder Atembeschwerden verursachen.*
- H 412 Harmful to aquatic life with long lasting effects.  
*Schädlich für Wasserorganismen, mit langfristiger Wirkung.*
- EUH031 Contact with acids liberates toxic gas.  
*Entwickelt bei Berührung mit Säure giftige Gase.*

### Precaution phrases

- P 210 Keep away from heat / sparks / open flames / hot surfaces – No smoking.  
*Von Hitze / Funken / offener Flamme / heißen Oberflächen fernhalten.*
- P 233 Keep container tightly closed  
*Behälter dicht verschlossen halten.*
- P 260 Do not breathe vapours.  
*Dampf nicht einatmen.*
- P 261 Avoid breathing dust.  
*Einatmen von Staub vermeiden.*
- P 273 Avoid release to the environment.  
*Freisetzung in die Umwelt vermeiden.*
- P 280 Wear protective gloves / eye protection.  
*Schutzhandschuhe / Augenschutz tragen.*
- P 301+312 IF SWALLOWED: Call a POISON CENTER or doctor / physician if you feel unwell.  
*BEI VERSCHLUCKEN: Bei Unwohlsein GIFTINFORMATIONSZENTRUM oder Arzt anrufen.*
- P 302+352 IF ON SKIN: Wash with plenty of soap and water.  
*BEI KONTAKT MIT DER HAUT: Mit viel Wasser und Seife waschen.*
- P 304+341 IF INHALED: If breathing is difficult, remove to fresh air and keep at rest in a position comfortable for breathing.  
*BEI EINATMEN: Bei Atembeschwerden an die frische Luft bringen und in einer Position ruhigstellen, die das Atmen erleichtert.*
- P 305+351+338 IF IN EYES: Rinse continuously with water for several minutes. Remove contact lenses if present and easy to do – continue rinsing.  
*BEI KONTAKT MIT DEN AUGEN: Einige Minuten lang behutsam mit Wasser spülen. Vorhandene Kontaktlinsen nach Möglichkeit entfernen. Weiter spülen.*
- P 330 Rinse mouth.  
*Mund ausspülen.*
- P 332+313 If skin irritation occurs: Get medical advice / attention.  
*Bei Hautreizung: Ärztlichen Rat einholen / ärztliche Hilfe hinzuziehen.*
- P 333+313 IF skin irritation or a rash occurs: Get medical advice / attention.  
*Bei Hautreizung oder -ausschlag: Ärztlichen Rat einholen / ärztliche Hilfe hinzuziehen.*

### Precaution phrases

- P 337+313 Get medical advice / attention.  
*Bei anhaltender Hautreizung: Ärztlichen Rat einholen / ärztliche Hilfe hinzuziehen.*
- P 342+311 If experiencing respiratory symptoms: Call a POISON CENTER or doctor / physician.  
*Bei Symptomen der Atemwege: GIFTINFORMATIONSZENTRUM oder Arzt anrufen.*
- P 403+235 Store in a well ventilated place. Keep cool.  
*Kühl an einem gut belüfteten Ort aufbewahren.*
- P 363 Wash contaminated clothing before reuse.  
*Kontaminierte Kleidung vor erneutem Tragen waschen.*

For further information please see Material Safety Data Sheets ([www.mn-net.com](http://www.mn-net.com)).  
*Weiterführende Informationen finden Sie in den Sicherheitsdatenblättern ([www.mn-net.com](http://www.mn-net.com)).*

## 5 Protocols

### 5.1 DNA, RNA, and protein purification from cultured cells and tissue

*Joint protocol steps for DNA, RNA, and protein purification.*

*DNA purification: steps 1–8*

*RNA purification: steps 1–11*

*Protein purification: steps 1–5 and 12–15*

#### Before starting the preparation:

- Check if Buffer *DNA Wash*, Wash Buffer RA3, rDNase, and Reducing Agent TCEP were prepared according to section 3.

#### 1 Homogenize sample

Disrupt up to **30 mg of human / animal tissue** or up to **100 mg of plant tissue** (for homogenization methods see section 2.3).



**Disrupt  
sample**

Up to **5 x 10<sup>6</sup> eukaryotic cultured cells** are collected by centrifugation and lysed by addition of Buffer RP1 directly.

#### 2 Lyse sample

Add **350 µL Buffer RP1** and **3.5 µL β-mercaptoethanol** (β-ME) to the cell pellet or to ground tissue and vortex vigorously.



**+ 350 µL RP1**  
**+ 3.5 µL  
β-ME**

*Note: As alternative to β-ME the reducing agent DTT or TCEP may be used. Use a final concentration of 10–20 mM DTT or TCEP within the Lysis Buffer RP1 (e.g. ad 7–14 µL of a 500 mM DTT or TCEP solution).*

### 3 Filtrate lysate

Reduce viscosity and clear the lysate by filtration through NucleoSpin® Filter: Place **NucleoSpin® Filter (violet ring)** in a Collection Tube (2 mL), apply the mixture, and centrifuge for **1 min** at **11,000 x g**.

*The lysate may be passed alternatively  $\geq 5$  times through a 0.9 mm needle (20 gauge) fitted to a syringe.*

In case of visible pellet formation (depending on sample amount and nature), transfer supernatant without any formed pellet to a new 2 mL centrifuge tube (not included).

***Important:** To process higher amounts of cells ( $> 1 \times 10^6$ ) or tissue ( $> 10$  mg), the lysate should first be homogenized using the 0.9 mm needle (20 gauge), followed by filtration through NucleoSpin® Filter.*



**1 min**  
**11,000 x g**



### 4 Adjust DNA and RNA binding conditions

Discard the NucleoSpin® Filter and add **350  $\mu$ L ethanol (70 %)** to the homogenized lysate and mix by pipetting up and down (approx. 5 times).

Alternatively, transfer flow-through into a new 1.5 mL microcentrifuge tube (not provided), add 350  $\mu$ L ethanol (70 %), and mix by vortexing (2 x 5 s).

*After addition of ethanol a stringy precipitate may become visible which will not affect the RNA isolation. Make sure to disaggregate any precipitate by mixing and load all of the disaggregated precipitate on the column as described in step 5. Do not centrifuge at this stage in order to avoid sedimentation of any precipitate.*



**+ 350  $\mu$ L**  
**70 % EtOH**

**Mix**



## 5 Bind DNA and RNA

For each preparation, take one **NucleoSpin® TriPrep Column (light blue ring)** placed in a Collection Tube and **load the lysate**. Centrifuge for **30 s** at **11,000 x g**. Place the NucleoSpin® TriPrep Column in a new Collection Tube (2 mL).

**RNA and DNA are bound to the column membrane. Protein is in the flow-through.**

*Maximal loading capacity of NucleoSpin® TriPrep Columns is 750 µL. Repeat the procedure if larger volumes are to be processed.*

For **DNA** and **RNA** isolation continue with step 6. *It is recommended to continue the DNA and RNA isolation protocol first and to perform the protein purification subsequently.*

For **protein** isolation recover flow-through and continue with step 12.

*The protein containing flow-through is stable for several hours at 4–8 °C.*



**Load sample**



**30 s**  
**11,000 x g**

## 6 Wash silica membrane

### 1<sup>st</sup> wash

Add **500 µL DNA Wash** to the NucleoSpin® TriPrep Column. Centrifuge for **1 min** at **11,000 x g**. Discard flow-through and reuse the Collection Tube.

### 2<sup>nd</sup> wash

Add again **500 µL DNA Wash** to the NucleoSpin® TriPrep Column. Centrifuge for **1 min** at **11,000 x g**. Discard Collection Tube with flow-through.

*Chaotropic salt is removed by these washing steps. DNA and RNA are still bound to the membrane. The membrane is prepared for subsequent DNA elution.*



**+ 500 µL**  
**DNA Wash**

**1 min**  
**11,000 x g**



**+ 500 µL**  
**DNA Wash**

**1 min**  
**11,000 x g**

## 7 Dry membrane

Insert the NucleoSpin® TriPrep Column into a 1.5 mL microcentrifuge tube (not provided). Open the lid of the column and let it stand for 3 min.

*This step ensures removal of residual ethanol.*

RT  
3 min  
(with open lid)

## 8 Elute DNA

Add **100 µL DNA Elute** directly onto the membrane and **incubate** for **1 min**. Elute DNA by centrifugation for **1 min** at **11,000 x g**.

*The temperature of the DNA Elute solution shall not exceed 30°C, otherwise RNA will partly elute with the DNA Elute solution. DNA Elute solution may stay for 1 min up to 15 min on the column before DNA is eluted. A 1–5 min incubation time is recommended. Eluted DNA is immediately ready for downstream applications without further purification.*



+ 100 µL  
DNA Elute

Incubate  
1 min



1 min  
11,000 x g

### Further steps for RNA purification (steps 9–11)

## 9 Digest residual DNA on-column

**Prepare rDNase reaction mixture** in a sterile microcentrifuge tube (not provided): For each isolation, add **10 µL** reconstituted **rDNase** (also see section 3) to **90 µL Reaction Buffer for rDNase**. Mix by flicking the tube.

Apply **95 µL rDNase reaction mixture** directly onto the center of the silica membrane of the column. Incubate at **room temperature** for **15 min**.

*Do not centrifuge, directly proceed with step 10. Depending on sample type and amount approximately 50–90 % of the DNA is eluted in the DNA elution step. Residual DNA on the column is digested on-column with rDNase.*



+ 95 µL  
rDNase  
reaction  
mixture

RT  
15 min

**10 Wash and dry silica membrane****1<sup>st</sup> wash**

Add **200 µL Buffer RA2** to the NucleoSpin® TriPrep Column. Centrifuge for **30 s** at **11,000 x g**. Place the NucleoSpin® TriPrep Column into a new Collection Tube (2 mL).

*Buffer RA2 will inactivate the rDNase.*

**+ 200 µL RA2**

**30 s**  
**11,000 x g**

**2<sup>nd</sup> wash**

Add **600 µL Buffer RA3** to the NucleoSpin® TriPrep Column. Centrifuge for **30 s** at **11,000 x g**. Discard flow-through and place the NucleoSpin® TriPrep Column back into the Collection Tube.



**+ 600 µL RA3**

**30 s**  
**11,000 x g**

**3<sup>rd</sup> wash**

Add **250 µL Buffer RA3** to the NucleoSpin® TriPrep Column. Centrifuge for **2 min** at **11,000 x g** to dry the membrane completely. Place the NucleoSpin® TriPrep Column into an RNase-free Collection Tube (1.5 mL, supplied).

*If for any reason, the liquid level in the Collection Tube has reached the NucleoSpin® TriPrep Column after centrifugation, discard flow-through and centrifuge again.*

**+ 250 µL RA3**

**2 min**  
**11,000 x g**

**11 Elute highly pure RNA**

Elute the RNA in 60 µL RNase-free H<sub>2</sub>O (supplied) and centrifuge at **11,000 x g** for **1 min**.

*If higher RNA concentrations are desired, elution can be done with 40 µL. Overall yield, however, will decrease when using smaller volumes.*

*For further alternative elution procedures see section 2.6.*



**+ 60 µL**  
**RNase-free**  
**H<sub>2</sub>O**

**1 min**  
**11,000 x g**

**Further steps for protein purification (steps 12–15)**

Use the NucleoSpin® TriPrep Column flow-through from step 5 (i.e., the ethanolic lysate which has been passed through the NucleoSpin® TriPrep Column) as starting point for protein precipitation.

## 12 Precipitate protein

**Transfer** an appropriate amount (10–700  $\mu\text{L}$ ) of flow-through into a fresh Collection Tube (1.5 mL, supplied).

*See section 2.4 as guideline for choosing an appropriate amount.*

Add **one volume PP** (Protein Precipitator). Mix vigorously.

**Incubate** mixture at room temperature for approximately 10 min.

*Note: For samples of moderate to high protein content (e.g., 100 mg young plant leaf, 30 mg liver) this incubation step may be omitted. For samples of low to medium protein content (e.g., 15 mg young plant leaf) the 10 min incubation increases protein yield relative to no incubation significantly. An incubation of longer than one hour does not further increase protein yield.*

Centrifuge for **5 min** at **11,000 x g**.

---

## 13 Wash protein

**Remove supernatant** by pipetting or decanting as complete as possible.

Add **500  $\mu\text{L}$**  of **50 % ethanol** to the pellet (mixing or incubation at this step is not necessary).

Centrifuge **1 min** at **11,000 x g**.

**Remove** supernatant by pipetting or decanting as completely as possible.

*Note: Protein precipitate at this stage is quite different in appearance depending on kind and amount of starting material. The consistence of precipitate might appear as no visible pellet (e.g., for 10,000 cells, 0.3 mg liver and 1 mg leaf samples); a greenish tube wall coating on one side of the tube (for e.g. leaf material); green or white pellet at the bottom of the tube (e.g., for leaf and liver samples, respectively); green or white crumbs at one side of the inner wall of the centrifuge tube (e.g., for leaf and liver samples, respectively). If no precipitate is visible, mark the side of the tube where a precipitate is expected in order to avoid touching this side of the inner tube wall with the pipette tip during the washing step. See also section 2.4 how to avoid very large protein pellets.*

---

## 14 Dry protein pellet

Dry precipitate for **5–10 min** at **room temperature**; keep lid open.

*Note: Large pellets (e.g., complete precipitation of 700  $\mu\text{L}$  column flow-through from a 30 mg liver sample) need longer time for drying. Samples which are dried incomplete may cause problems when loading the sample onto the gel due to residual ethanol content. No problems with over-drying have been observed with small-sized pellets. See also section 2.4 how to avoid very large protein pellets.*

---

## 15 Prepare protein sample

Add **20–100 µL PSB-TCEP** (Protein Solving Buffer, containing reducing agent).

*Assure that PSB-TCEP is clear (not turbid). If necessary, warm PSB-TCEP to > 25 °C to dissolve turbidity.*

*See section 2.4 as guideline for choosing an appropriate amount of PSB-TCEP for dissolving of protein pellets.*

*Disaggregate large and visible pellets with a pipette tip to facilitate subsequent protein dissolution; this is not necessary for small and invisible pellets.*

**Incubate** for **3 min** at **95–98 °C** for complete protein dissolving and denaturation.

Let sample cool down to room temperature.

---

**Centrifuge** for **1 min** at **11,000 x g** to pellet residual insoluble material.

*Note: Depending on sample amount and nature there might be no visible pellet of insoluble material up to large pellets of different size and structure. Do not disturb residual precipitates at this stage. Protein will be in the supernatant. Do not centrifuge samples at temperatures < 18 °C. SDS may precipitate at this temperature.*

Recover supernatant for further analysis. See section 6.1 for suitable protein quantification methods.

*Note: At this stage samples can be stored at -20 °C for several months or at 4 °C for several days. After storage, equilibrate sample to room temperature, mix, and then centrifuge briefly before withdrawal of sample aliquots. Repeated sample denaturing for 3 min at 95–98 °C is not necessary. Repetitive withdrawal, freezing, and thawing for at least three times has shown constant sample quality.*

*Due to the strong denaturing purification method protein is precipitated in denatured form with reduced solubility in water. Therefore resolubilization of the protein pellet in PSB-TCEP or in traditional Laemmli buffer is recommended. The use of Protein Solving Buffer PSB is not mandatory for dissolving protein. Alternatively to PSB, PSB-TCEP, or Laemmli buffer, precipitated protein can be dissolved in 1% SDS or 8 M urea. Further, the protein pellet can be dissolved in urea/thiourea/CHAPS buffers as used for 2-D electrophoresis. However, depending on the target protein, the overall yield of solubilized protein may be reduced compared to PSB or PSB-TCEP used as dissolving agent.*

---

## 5.2 Total RNA preparation from RNA/ater® treated samples

### Before starting the preparation:

- Check if Buffer *DNA Wash*, Wash Buffer RA3, rDNase, and Reducing Agent TCEP were prepared according to section 3.
- 

#### 1 Prepare sample

Remove RNA/ater® solution. Cut an appropriate amount of tissue.

---

#### 2 Lyse sample

Add **350 µL Buffer RP1** and **3.5 µL β-mercaptoethanol** (β-ME) to the sample.

Disrupt the sample material by using for example rotor-stated homogenizers (for homogenization methods see section 2.3).

*Note: As alternative to β-ME the reducing agent DTT or TCEP may be used. Use a final concentration of 10–20 mM DTT or TCEP within the Lysis Buffer RP1.*

---

Proceed with step 3 of the NucleoSpin® TriPrep standard protocol (section 5.1).

### 5.3 rDNase digestion in solution

The on-column rDNase digestion in the standard protocol is very efficient and results in minimal residual DNA content of the purified RNA. This DNA will not be detectable in most downstream applications, but there are still certain applications which require even lower contents of residual DNA. However, removal of DNA to a completely undetectable level is challenging and the efficiency of an on-column DNA digestion is sometimes not sufficient for downstream applications requiring lowest residual content of DNA.

A typical example for such a demanding application is an RT-PCR reaction in which the primer molecules do not differentiate between cDNA (derived from RNA) and contaminating genomic DNA. Especially, if

- high copy number targets are analyzed (e.g., multi gene family, mitochondrial, plasmid or plasmid targets (from transfections))
- the target gene is of a very low expression level
- the amplicon is relatively small (< 200 bp).

DNA digestion in solution can efficiently destroy contaminating DNA. However, stringent RNase control and subsequent repurification of the RNA (in order to remove buffer, salts, DNase and digested DNA) are usually required.

The high quality, recombinant RNase-free DNase (rDNase) in the NucleoSpin® TriPrep kits facilitates such a digestion in solution in order to remove even traces of contaminating DNA.

---

#### **A Digest DNA (reaction setup)**

Add **6 µL Reaction Buffer for rDNase** and **0.6 µL rDNase** to **60 µL eluted RNA**.

*Alternatively, premix 100 µL Reaction Buffer for rDNase and 10 µL rDNase and add 1/10 volume to one volume of RNA eluate.*

---

#### **B Incubate sample**

Incubate for **10 min** at **37 °C**.

---

## **C Repurify RNA**

Repurify RNA with a suitable RNA clean-up procedure, for example by use of the NucleoSpin® RNA Clean-up or NucleoSpin® RNA Clean-up XS kit (see ordering information) or by ethanol precipitation.

### **Ethanol precipitation, exemplary**

Add 0.1 volume of 3 M sodium acetate, pH 5.2 and 2.5 volumes of 96–100% ethanol to one volume of sample. Mix thoroughly.

Incubate several minutes to several hours at -20 °C or 4 °C.

*Note: Choose long incubation times if the sample contains low RNA concentration. Short incubation times are sufficient if the sample contains high RNA concentration.*

Centrifuge for 10 min at max. speed.

Wash RNA pellet with 70% ethanol.

Dry RNA pellet and resuspend RNA in RNase-free H<sub>2</sub>O.

---



## 6 Appendix

### 6.1 Protein quantification

Quantification of protein dissolved in sample buffer such as PSB, PSB-TCEP or traditional Laemmli buffer is occasionally required prior to SDS-PAGE and Western Blot analysis. However, major protein quantification assays are influenced by or are incompatible with SDS and/or reducing agents commonly present in protein sample buffers used for SDS-PAGE. A protein quantification procedure has to be chosen carefully to ensure appropriate compatibility of the method with the protein dissolution solution.

The NucleoSpin® RNA/Protein procedure allows several protein quantification methods at different steps of the procedure:

- Quantification of protein dissolved in PSB, PSB-TCEP or Laemmli buffer (recommended)
- Quantification of protein dissolved in alternative protein dissolution buffers (1 % SDS or 8 M urea)
- Quantification of protein within the column flow-through (i.e., prior to protein precipitation, at step 5 of the standard protocol)

**For most reliable results and convenience we recommend the MACHEREY-NAGEL Protein Quantification Assay** (for ordering information see section 6.4) to quantify protein dissolved in PSB, PSB-TCEP, or Laemmli buffer.

#### **Quantification of protein dissolved in PSB or PSB-TCEP**

The concentration of protein dissolved in PSB, PSB-TCEP or Laemmli buffer can be determined with several methods.

Below you find a choice of quantification methods, which are compatible with PSB, PSB-TCEP and Laemmli buffer, but not all of the methods show the same sensitivity. The following list compares the different sensitivities and gives assistance to find out the most suitable protein quantification assay.

#### **1: Protein Quantification Assay (highly recommended method!)**

This is the most sensitive and convenient method for protein quantification in PSB or PSB-TCEB. Highly recommended due to sensitivity and high compatibility. For ordering information see section 6.4.

#### **2: Method adapted from the publication Karlsson *et al.* 1994**

See below for details.

#### **3: Pierce BCA Protein Assay Kit (reducing agent compatible)**

Dilute the protein sample 1:5 with water to enable compatibility.

#### 4: Bio-Rad DC Protein Assay

Dilute the protein sample 1:10 with water to enable compatibility.

#### 5: Bio-Rad RC DC Protein Assay

Dilute the protein sample 1:5 with water to enable compatibility.

#### 6: Serva ProtQuant Assay Kit

According to manufacturer's instructions, this assay should be compatible with PSB, PSB-TCEP, and Laemmli buffer samples.

#### 7: G-Biosciences SPN™-Protein Assay

According to manufacturer's instructions, this assay should be compatible with PSB, PSB-TCEP, and Laemmli buffer samples.

#### 8: Bio-Rad Protein Assay (Bradford)

This method has a very low tolerance towards SDS (0.1 % SDS for the Standard Assay Procedure). Therefore, PSB, PSB-TCEP Laemmli buffer samples have to be diluted considerably with water to reduce interference. After dilution of the sample 1:20 with water protein can be quantified with the standard assay procedure. The microassay procedure however, is not compatible with such samples, even after 1:50 dilution of the sample with water.

### Compatibility of protein quantification methods with PSB and PSB-TCEP samples

| # | Protein quantification assay                              | Input sample volume (e.g., PSB-TCEP sample)                             | Acceptable protein amount per assay | Protein concentration quantification range (undiluted PSB-TCEP sample) |
|---|---|---|-------------------------------------|--|
| 1 | Protein Quantification Assay* (highly recommended)        | 20 µL standard (1–60 µL optional)                                       | 0.6–20 µg                           | 0.03–1 µg/µL (standard)<br>0.01–20 µg/µL (optional)                    |
| 2 | Karlsson protocol (recommended)*                          | 20 µL standard (1–60 µL optional)                                       | 0.6–20 µg                           | 0.03–1 µg/µL (standard)<br>0.01–20 µg/µL (optional)                    |
| 3 | Pierce BCA Protein Assay Kit (reducing Agent compatible)* | 25 µL of a 1:5 prediluted sample, corresponding to 5 µL original sample | 3.125–50 µg                         | 0.625–10 µg/µL   |

\* Method tested in MN laboratory for compatibility with PSB, PSB-TCEP and Laemmli buffer protein samples.

| #  | Protein quantification assay                                 | Input sample volume (e.g., PSB-TCEP sample)   | Acceptable protein amount per assay                   | Protein concentration quantification range (undiluted PSB-TCEP sample) |
|----|--|---|---|--|
| 4  | Bio-Rad DC Protein Assay*                                    | 100 $\mu$ L of a 1:10 prediluted sample, corresponding to 10 $\mu$ L original sample (standard assay) or 5 $\mu$ L of a 1:10 prediluted sample, corresponding to 0.5 $\mu$ L original sample (micro testtube assay) | 20–150 $\mu$ g (standard)<br>1.0–7.5 $\mu$ g (micro)  | 2–15 $\mu$ g/ $\mu$ L  |
| 5  | Bio-Rad RC DC Protein Assay*                                 | 100 $\mu$ L of a 1:5 prediluted sample, corresponding to 20 $\mu$ L original sample or (standard assay)<br>25 $\mu$ L of a 1:5 prediluted sample, corresponding to 5 $\mu$ L original sample (micro testtube assay) | 20–150 $\mu$ g (standard)<br>5.0–37.5 $\mu$ g (micro) | 1–7.5 $\mu$ g/ $\mu$ L   |
| 6  | Serva ProtQuant Assay Kit**                                  | 20 $\mu$ L  | 5–35 $\mu$ g  | 0.25–1.75 $\mu$ g/ $\mu$ L   |
| 7  | G-Bioscience SPN™-Protein Assay**                            | 1–10 $\mu$ L  | 0.5–10 $\mu$ g  | 0.05–10 $\mu$ g/ $\mu$ L   |
| 8a | Bio-Rad Protein Assay – Standard Assay Procedure* (Bradford) | 100 $\mu$ L of 1:20 prediluted sample, corresponding to 5 $\mu$ L original sample   | 20–140 $\mu$ g  | 4–28 $\mu$ g/ $\mu$ L  |
| 8b | Bio-Rad Protein Assay – Microassay Procedure* (Bradford)     | Not recommended   | –   | –  |

\* Method tested in MN laboratory for compatibility with PSB, PSB-TCEP, and Laemmli buffer protein samples.

\*\* Method compatible with PSB-TCEP protein samples referring to manufacturer's product information. Not tested in MN laboratories.

### Quantification of protein within the column flow-through

Alternative to quantification of protein dissolved in PSB-TCEP, protein can be quantified within the ethanolic lysate column flow-through. Knowledge of protein concentration in the column flow-through helps to choose an appropriate volume for subsequent precipitation with Protein Precipitator PP. The following methods are suitable to quantify protein in the column flow-through:

#### 1: Bio-Rad Protein Assay (Bradford)

The standard assay procedure is compatible with the column flow-through, however, protein standards have to be prepared in a ethanolic lysate buffer (mix Buffer RP1 and ethanol (70 %) in a ratio of 1:1). The microassay procedure is not compatible the column-flow-through!

#### 2: Pierce BCA Protein Assay Kit (reducing agent compatible)

This method is compatible with the column flow-through, however, protein standards have to be prepared in a ethanolic lysate buffer (mix Buffer RP1 and ethanol (70 %) in a ratio of 1:1).

#### 3: Bio-Rad DC Protein Assay

This method is compatible with the column flow-through, however, protein standards have to be prepared in a ethanolic lysate buffer (mix Buffer RP1 and ethanol (70 %) in a ratio of 1:1).

#### 4: Bio-Rad RC DC Protein Assay

This method is compatible with the column flow-through, however, protein standards have to be prepared in a ethanolic lysate buffer (mix Buffer RP1 and ethanol (70 %) in a ratio of 1:1).

#### 5: Roti-Nanoquant Assay

This method is compatible with the column flow-through, however, protein standards have to be prepared in a ethanolic lysate buffer (mix Buffer RP1 and ethanol (70 %) in a ratio of 1:1).

### Protein quantification in alternative protein pellet dissolution buffers

The use of the PSB or PSB-TCEP buffer is not mandatory for solving proteins.

Precipitated protein (protein pellets) may be dissolved in alternative solutions, such as **1 % SDS**, or **8 M urea**, or in **urea / thiourea / CHAPS buffers** as used for 2-D electrophoresis. However, depending on the target protein, the yield of solubilized protein may be reduced compared to PSB or PSB-TCEP. Check manufacturers product information to ensure compatibility of your protein quantification assay with your alternative protein dissolving solution.

### Protocol for protein quantification

For most reliable results and convenience we recommend the MACHEREY-NAGEL Protein Quantification Assay to quantify protein dissolved in PSB, PSB-TCEP, or Laemmli buffer (for ordering information see section 6.4).

The procedure presented below (based on the publication of Karlsson *et al.* 1994), is also suitable for quantification of protein dissolved in Protein Solving Buffer PSB-TCEP and may be followed alternatively.

Nucleic acids disturb protein quantification as described by Karlsson *et al.* 1994. Protein samples obtained with the NucleoSpin® RNA/Protein kit are virtually free of nucleic acids, thus, protein quantification is not affected.

Upon addition of TCA (Trichloroacetic acid) to the sample, protein precipitates and causes turbidity. The degree of turbidity is used for quantification relative to a sample with known protein concentration. This test enables determination of protein concentration in the range 5 ng/ $\mu\text{L}$ –20  $\mu\text{g}/\mu\text{L}$  by using variable sample volumes of 1–60  $\mu\text{L}$ .

| Recommended sample volume<br>(protein dissolved in PSB-TCEP) | For protein<br>concentration in the<br>range of | Protein amount<br>per well |
|--|---|----------------------------|
| 60 $\mu\text{L}$   | 0.01–0.33 $\mu\text{g}/\mu\text{L}$             | 0.6–20 $\mu\text{g}$       |
| 20 $\mu\text{L}$   | 0.03–1.00 $\mu\text{g}/\mu\text{L}$             | 0.6–20 $\mu\text{g}$       |
| 1 $\mu\text{L}$  | 0.6–20 $\mu\text{g}/\mu\text{L}$                | 0.6–20 $\mu\text{g}$       |

For the primary determination of protein concentration of the sample measure different amounts, for example 2  $\mu\text{L}$ , 10  $\mu\text{L}$ , and 50  $\mu\text{L}$ . This assures that at least one value of the three amounts tested falls within the range of the calibration curve. Further, for a first rough estimation of expected protein yield also consider table 4 in section 2.2 and table 5 in section 2.4.

### Materials

- TCA 60 % (Trichloroacetic acid, not supplied)
- Protein Solving Buffer with reducing agent (PSB-TCEP, see note below)
- BSA (Bovine Serum Albumin, not supplied)
- Multititer plate (not supplied)

*Note: The volume of PSB-TCEP, provided with the kit, might not be sufficient to quantify all isolated protein samples. Additional PSB-TCEP can either be ordered separately (see ordering information) or easily be prepared (see composition of PSB-TCEP below).*

### Composition of PSB-TCEP

- 125 mM BisTris (Bis-(2-hydroxyethyl)-imino-tris(hydroxymethyl)-methane)
- 2% SDS (sodium dodecyl sulphate)
- 50 mM TCEP (Tris(2-carboxyethyl)phosphine Hydrochloride)
- 20% glycerol
- 0.01% bromophenol blue
- pH 6.8

*Note: The composition of the previously used Protein Loading Buffer (PLB) has been improved and is now called Protein Solving Buffer PSB (reduced concentration of SDS and bromophenol blue).*

*The change in composition has increased the compatibility with protein quantification methods (see above). For details on the composition of previous Protein Loading Buffer (PLB) contact our technical service.*

### Method

Prepare a BSA stock solution with 40 mg/mL BSA in H<sub>2</sub>O.


Prepare a BSA dilution series:

| Tube | Add PSB to tube | Add BSA solution to tube                  | Resulting BSA concentration | Resulting BSA in 20 $\mu$ L |
|------|-----------------|---|-----------------------------|-----------------------------|
| #1   | 97.5 $\mu$ L    | 2.5 $\mu$ L BSA stock solution (40 mg/mL) | 1 $\mu$ g/ $\mu$ L          | 20 $\mu$ g                  |
| #2   | 50 $\mu$ L      | 50 $\mu$ L from tube #1                   | 0.5 $\mu$ g/ $\mu$ L        | 10 $\mu$ g                  |
| #3   | 50 $\mu$ L      | 50 $\mu$ L from tube #2                   | 0.25 $\mu$ g/ $\mu$ L       | 5 $\mu$ g                   |
| #4   | 50 $\mu$ L      | 50 $\mu$ L from tube #3                   | 0.125 $\mu$ g/ $\mu$ L      | 2.5 $\mu$ g                 |
| #5   | 50 $\mu$ L      | 50 $\mu$ L from tube #4                   | 0.063 $\mu$ g/ $\mu$ L      | 1.25 $\mu$ g                |
| #6   | 50 $\mu$ L      | 50 $\mu$ L from tube #5                   | 0.031 $\mu$ g/ $\mu$ L      | 0.625 $\mu$ g               |
| #7   | 50 $\mu$ L      | –   | 0 $\mu$ g/ $\mu$ L          | 0 $\mu$ g                   |

Make sure that the protein concentration of your sample lies within the range of the largest (#1) and smallest (#6) value of the calibration curve in order to obtain valid measurements. Do not extrapolate beyond this range!

The prepared BSA dilution series is sufficient for subsequent determination of two calibration curves.

- 1 Add 20  $\mu\text{L}$  of each dilution series sample # 1–7 in microtiter plate wells.
- 2 Add 20  $\mu\text{L}$  of samples (protein dissolved in PSB-TCEP) with unknown protein concentration to further wells (alternatively 1–60  $\mu\text{L}$ ).
- 3 Add 40  $\mu\text{L}$  PSB-TCEP to each well. Final volume: 60  $\mu\text{L}$  (alternatively add 0–59  $\mu\text{L}$  if other volumes than 20  $\mu\text{L}$  of sample are used in step 2).
- 4 Add 40  $\mu\text{L}$  TCA (60%) to each well.
- 5 Mix until complete colour change from blue to yellow.
- 6 Incubate for 30 min ( $\pm$  3 min) at room temperature.
- 7 Measure extinction at 570 nm\*.
- 8 Determine protein concentration of samples in relation to dilution series.

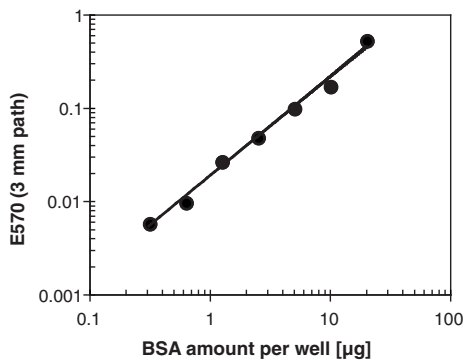

**Figure 1:** BSA standard curve for determination of protein in Protein Solving Buffer PSB.

\* Measurement of extinction in the range of 530–700 nm is suitable and will typically result in correlation coefficients of 0.99 (concentration of BSA dilution series vs. obtained absorption values).

## 6.2 Troubleshooting

| Problem   | Possible cause and suggestions  |
|---|---|
| DNA is contaminated with RNA                            | <p><i>Buffer temperature</i></p> <ul style="list-style-type: none"> <li>• DNA elution buffer <i>DNA Elute</i> exceeded 30 °C during application. Use <i>DNA Elute</i> with a temperature preferentially of 18–25 °C.</li> </ul>   |
| DNA yield lower than RNA yield                          | <p><i>Sample material</i></p> <ul style="list-style-type: none"> <li>• DNA and RNA yield depend very much on sample material. Ratio of RNA yield to DNA yield may vary from approximately 1–20.</li> </ul>  |
| DNA degrades upon storage                               | <p><i>DNase contamination</i></p> <ul style="list-style-type: none"> <li>• DNA elution buffer <i>DNA Elute</i> does not contain divalent cation complexing substances (e.g., EDTA). Therefore, DNA is not protected against DNases. Keep <i>DNA Elute</i> solution clean and avoid any contamination. As a precaution, keep DNA on ice for short term or at -20 °C for long term storage.</li> <li>• Some sample materials may contain DNase traces which are not sufficiently washed away by the standard procedure. Perform a wash step of the column with Buffer RA2 after loading the lysate onto the column and before starting the washing steps with <i>DNA Wash</i> solution: add 500 µL Buffer RA2 onto the column, centrifuge 1 min at 11,000 x g and continue with <i>DNA Wash</i> washing steps.</li> </ul> |
| Suboptimal performance of DNA in downstream application | <p><i>Divalent cations</i></p> <ul style="list-style-type: none"> <li>• Eluted DNA contains small amounts of divalent cations. If the downstream application comprises, for example 50% DNA eluate of the final reaction volume the divalent cations introduced into the reaction by the DNA eluate may alter the performance. Decrease the divalent cation concentration of the reaction by 1–5 mM for compensation.</li> </ul>  |
| Low DNA yield for large sample amounts                  | <p><i>Sample amount to large</i></p> <ul style="list-style-type: none"> <li>• Depending on the type of sample and its DNA content, DNA yield may not increase proportional to increased sample amount. Sample amounts larger than, for example 5 mg tissue or 10<sup>6</sup> cultured cells may yield less DNA than smaller sample amounts. Use smaller sample to ensure good correlation between sample amount and DNA yield.</li> </ul>   |



**Problem**

**Possible cause and suggestions**

---

RNA is degraded/  
no RNA  
obtained

*RNase contamination*

- Create an RNase-free working environment. Wear gloves during all steps of the procedure. Change gloves frequently. Use of sterile disposable polypropylene tubes is recommended. Keep tubes closed whenever possible during the preparation. Glassware should be oven-baked for at least 2 hours at 250 °C before use.
- 

Poor RNA  
quality or yield

*Reagents not applied or restored properly*

- Reagents not properly restored. Add the indicated volume of RNase-free water to rDNase vial and 96% ethanol to Buffer RA3 Concentrate and mix. Reconstitute and store lyophilized rDNase according to instructions given in section 3.
- Sample and reagents have not been mixed completely. Always vortex vigorously after each reagent has been added.
- No ethanol has been added after lysis. Binding of RNA to the silica membrane is only effective in the presence of ethanol.

*Kit storage*

- Reconstitute and store lyophilized rDNase according to instructions given in section 3.
- Store other kit components at room temperature. Storage at low temperatures may cause salt precipitation.
- Keep bottles tightly closed in order to prevent evaporation or contamination.

*Ionic strength and pH influence  $A_{260}$  absorption as well as ratio  $A_{260}/A_{280}$*

- For adsorption measurement, use 5 mM Tris pH 8.5 as diluent. Please see also:
    - Manchester, KL, 1995. Value of  $A_{260}/A_{280}$  ratios for measurement of purity of nucleic acids. *Biotechniques* 19, 208–209.
    - Wilfinger, WW, Mackey, K and Chomczynski, P, 1997. Effect of pH and ionic strength on the spectrophotometric assessment of nucleic acid purity. *Biotechniques* 22, 474-481.
-

| <b>Problem</b> | <b>Possible cause and suggestions</b> |
|----------------|---------------------------------------|
|----------------|---------------------------------------|

---

|   |   |
|---|---|
| Poor RNA quality or yield<br><i>(continued)</i> | <p><i>Sample material</i></p> <ul style="list-style-type: none"><li>• Sample material was not stored properly. Whenever possible, use fresh material. If this is not possible, flash freeze the samples in liquid N<sub>2</sub>. Samples should always be kept at -70 °C. Never allow tissues to thaw before addition of Buffer RP1. Perform disruption of samples in liquid N<sub>2</sub>.</li><li>• Insufficient disruption and/or homogenization of starting material: Ensure thorough sample disruption and use NucleoSpin® Filters for easy homogenization of disrupted starting material.</li></ul> |
|---|---|

---

|  |  |
|--|--|
| Clogged NucleoSpin® Column/<br>Poor RNA quality or yield | <p><i>Sample material</i></p> <ul style="list-style-type: none"><li>• Too much starting material used: Overloading may lead to decreased overall yield. Reduce amount of sample material or use larger volume of Buffer RP1.</li><li>• Insufficient disruption and/or homogenization of starting material: Ensure thorough sample disruption and use NucleoSpin® Filters for easy homogenization of disrupted starting material.</li></ul> |
|--|--|

---

|                                       |   |
|---------------------------------------|---|
| Contamination of RNA with genomic DNA | <p><i>rDNase not active</i></p> <ul style="list-style-type: none"><li>• Reconstitute and store lyophilized rDNase according to instructions given in section 3.</li></ul> <p><i>rDNase solution not properly applied</i></p> <ul style="list-style-type: none"><li>• Pipette rDNase solution directly onto the center of the silica membrane.</li></ul> |
|---------------------------------------|---|

---

|  |  |
|--|--|
|  | <p><i>Too much cell material used</i></p> <ul style="list-style-type: none"><li>• Reduce quantity of cells or tissue used.</li></ul> |
|--|--|

---

| Problem | Possible cause and suggestions |
|---------|--------------------------------|
|---------|--------------------------------|

|  |  |
|--|--|
| Contamination of RNA with genomic DNA<br>(continued) | <p data-bbox="325 202 689 233"><i>DNA detection system too sensitive</i></p> <ul data-bbox="325 239 984 691" style="list-style-type: none"> <li data-bbox="325 239 984 564">• The amount of DNA contamination is significantly reduced during the on-column rDNase digestion. Anyhow we can not guarantee that the purified RNA is 100% free of DNA. In very sensitive applications it might be possible to detect DNA. The NucleoSpin® TriPrep system is checked by the following procedure: One million HeLa cells are subjected to RNA isolation according to the protocol. RNA eluate is used as template for PCR detection of a 1 kb fragment in a 30 cycle reaction. Generally, no PCR fragment is obtained if the rDNase is applied, however, a strong PCR fragment is obtained if rDNase is omitted. The probability of DNA detection with PCR increases with:                             <ul data-bbox="364 571 984 691" style="list-style-type: none"> <li data-bbox="364 571 984 655">- the number of DNA copies per preparation: single copy target &lt; plasmid / mitochondrial target &lt; plasmid transfected into cells</li> <li data-bbox="364 662 698 691">- decreasing PCR amplicon size.</li> </ul> </li> </ul> |
|--|--|

|  |  |
|--|--|
| Contamination of RNA with genomic DNA<br>(continued) | <ul data-bbox="325 722 982 845" style="list-style-type: none"> <li data-bbox="325 722 982 778">• Use larger PCR targets (e.g., &gt; 500 bp) or intron spanning primers if possible.</li> <li data-bbox="325 785 982 845">• <b>Use support protocol 5.3 for subsequent rDNase digestion in solution.</b></li> </ul> |
|--|--|

|   |   |
|---|---|
| Suboptimal performance of RNA in downstream experiments | <p data-bbox="325 880 617 911"><i>Carry-over of ethanol or salt</i></p> <ul data-bbox="325 917 982 1121" style="list-style-type: none"> <li data-bbox="325 917 982 1026">• Do not let the flow-through touch the column outlet after the second RA3 wash. Be sure to centrifuge at the corresponding speed for the respective time in order to remove ethanolic Buffer RA3 completely.</li> <li data-bbox="325 1032 982 1121">• Check if Buffer RA3 has been equilibrated to room temperature before use. Washing at lower temperatures lowers efficiency of salt removal by Buffer RA3.</li> </ul> |
|   | <p data-bbox="325 1157 614 1187"><i>Store isolated RNA properly</i></p> <ul data-bbox="325 1193 982 1331" style="list-style-type: none"> <li data-bbox="325 1193 982 1331">• Eluted RNA should always be kept on ice for optimal stability since trace contaminations of omnipresent RNases (general lab ware, fingerprints, dust) will degrade the isolated RNA. For short term storage freeze at -20 °C, for long term storage freeze at -70 °C.</li> </ul>   |

| <b>Problem</b>   | <b>Possible cause and suggestions</b>   |
|--|---|
| Trouble with resolubilization of precipitated protein in PSB-TCEP                      | <p><i>Protein pellets exceeding several millimeters in size are hard to dissolve.</i></p> <ul style="list-style-type: none"> <li>• Use smaller volumes of column flow-through for protein precipitation in order to obtain small sized pellets. Even invisible protein pellets commonly yield enough protein for SDS-PAGE and Western Blot analysis.</li> </ul>   |
| Protein dissolved in PSB-TCEP flows out of SDS-PAGE gel slot immediately after loading | <p><i>Protein pellet has not been dried sufficiently and contains residual ethanol.</i></p> <ul style="list-style-type: none"> <li>• Increase drying time or decrease pellet size by precipitating a smaller volume of column flow-through.</li> </ul>  |
| Unclear results with commonly used protein quantification systems                      | <p><i>Most commonly used protein quantification systems are incompatible with concentrations of SDS and/or reducing agents present in Protein Solving Buffer.</i></p> <ul style="list-style-type: none"> <li>• Use a suitable quantification method as described in section 6.1.</li> <li>• If an other protein dissolution buffer than PSB or PSB-TCEP was used for dissolving the protein pellet, ensure compatibility of your buffer and quantification method of choice.</li> </ul> |
| No protein precipitate (pellet) visible  | <p><i>A small sample amount was used and/or a small volume of column flow-through was used for precipitation.</i></p> <ul style="list-style-type: none"> <li>• Formation of a visible protein pellet is not required for sufficient protein recovery. Even invisible protein pellets commonly yield enough protein for SDS-PAGE and Western Blot analysis.</li> </ul>   |
| PSB-TCEP turbid or partially solidified  | <p><i>PSB-TCEP may form a precipitate at temperatures below 18°C.</i></p> <ul style="list-style-type: none"> <li>• Warm up <math>\geq 25^{\circ}\text{C}</math> to dissolve turbidity completely.</li> </ul>  |
| No/low protein yield   | <p><i>Protein was resolubilized in water</i></p> <ul style="list-style-type: none"> <li>• Due to the strongly denatured form of the protein, the solubility in water is significantly reduced. Use PSB-TCEP for protein solubilization.</li> </ul>  |

## 6.3 References

**The following publications show the general usefulness of the parallel extraction of DNA, RNA, and protein from small and precious samples:**

**Coombs LM**, Pigott D, Proctor A, Eydmann M, Denner J, and Knowles MA (1990): Simultaneous isolation of DNA, RNA, and antigenic protein exhibiting kinase activity from small tumor samples using guanidine isothiocyanate. *Analytical Biochemistry* 188, 338-343.

**Banerjee S**, Smallwood A, Chambers AE, and Nicolaides K (2003) : Quantitative recovery of immunoreactive proteins from clinical samples following RNA and DNA isolation. *BioTechniques* 35 (3), 450-456.

**Hoemann CD**, Sun J, Chrzanowski V, and Buschmann MD (2002): A multivalent assay to detect glycosaminoglycan, protein, collagen, RNA, and DNA content in milligram samples of cartilage or hydrogel-based repair cartilage. *Analytical Biochemistry* 300, 1-10.

**The following publications cite the use of the NucleoSpin® kits for DNA, RNA, and protein isolation:**

**Rodríguez-Jiménez FJ**, Moreno-Manzano V, Lucas-Dominguez R, and Sánchez-Puelles JM (2008): Hypoxia Causes Down-Regulation of Mismatch Repair System and Genomic Instability in Stem Cells. *Stem Cells*, May 2008; 10.1634/stemcells.2007-1016.

**Bahn A**, Hagos Y, Reuter S, Balen D, Brzica H, Krick W, Burckhardt BC, Sabolic I, and Burckhardt G (2008): Identification of a new urate and high affinity nicotinate transporter - human organic anion transporter 10 (hOAT10, SLC22A13). *J. Biol. Chem.* published 14 April 2008, 10.1074/jbc.M800737200

**Weiske J**, Albring KF, and Huber O (2007): The tumor suppressor Fhit acts as a repressor of  $\beta$ -catenin transcriptional activity. *PNAS*, Dec 2007; 104: 20344 - 20349.

**The following publication describes the Reducing Agent TCEP:**

**Getz EB**, Xiao M, Chakrabarty T, Cooke R and Selvin PR (1999): A comparison between the sulfhydryl reductants Tris(2-carboxyethyl)phosphine and Dithiothreitol for use in protein biochemistry. *Analytical Biochemistry* 273, 73-80.

**The following publication describes a method for quantification of protein dissolved in sample buffer such as PSB:**

**Karlsson JO**, Ostwald K, Kabjörn C, and Andersson M (1994): A method for protein assay in Laemmli buffer. *Analytical Biochemistry* 219, 144-146.

## 6.4 Ordering information

| Product                                | REF                | Pack of   |
|--|--------------------|-----------|
| NucleoSpin® TriPrep*                   | 740966.10/.50/.250 | 10/50/250 |
| NucleoSpin® RNA/Protein                | 740933.10/.50/.250 | 10/50/250 |
| Protein Quantification Assay           | 740967.50/.250     | 50/250    |
| NucleoSpin® RNA II                     | 740955.10/.50/.250 | 10/50/250 |
| NucleoSpin® RNA XS                     | 740902.10/.50/.250 | 10/50/250 |
| NucleoSpin® RNA L                      | 740962.20          | 20        |
| NucleoSpin® RNA Clean-up               | 740948.10/.50/.250 | 10/50/250 |
| NucleoSpin® RNA Clean-up XS            | 740903.10/.50/.250 | 10/50/250 |
| NucleoSpin® RNA/DNA Buffer Set*        | 740944             | 100       |
| Buffer RP1                             | 740934.50          | 50 mL     |
| Buffer RP1                             | 740934.500         | 500 mL    |
| Protein Solving Buffer Set<br>PSB/TCEP | 740941             | 1 set     |
| rDNase Set                             | 740963             | 1 set     |
| NucleoSpin® Filters                    | 740606             | 50        |
| Collection Tubes (2 mL)                | 740600             | 1000      |
| Porablot                               | See price list     |           |
| Blotting Paper                         | See price list     |           |

Visit [www.mn-net.com](http://www.mn-net.com) for more detailed product information.

\* DISTRIBUTION AND USE OF NUCLEOSPIN® TRIPREP and NUCLEOSPIN® RNA/DNA BUFFER SET IN THE USA IS PROHIBITED FOR PATENT REASONS.

## 6.5 Product use restriction/warranty

**NucleoSpin® TriPrep** kit components are intended, developed, designed, and sold FOR RESEARCH PURPOSES ONLY, except, however, any other function of the product being expressly described in original MACHEREY-NAGEL product leaflets.

MACHEREY-NAGEL products are intended for GENERAL LABORATORY USE ONLY! MACHEREY-NAGEL products are suited for QUALIFIED PERSONNEL ONLY! MACHEREY-NAGEL products shall in any event only be used wearing adequate PROTECTIVE CLOTHING. For detailed information please refer to the respective Material Safety Data Sheet of the product! MACHEREY-NAGEL products shall exclusively be used in an ADEQUATE TEST ENVIRONMENT. MACHEREY-NAGEL does not assume any responsibility for damages due to improper application of our products in other fields of application. Application on the human body is STRICTLY FORBIDDEN. The respective user is liable for any and all damages resulting from such application.

DNA/RNA/PROTEIN purification products of MACHEREY-NAGEL are suitable for *IN-VITRO*-USES ONLY!

ONLY MACHEREY-NAGEL products specially labeled as IVD are also suitable for *IN-VITRO*-diagnostic use. Please pay attention to the package of the product. *IN-VITRO*-diagnostic products are expressly marked as IVD on the packaging.

IF THERE IS NO IVD SIGN, THE PRODUCT SHALL NOT BE SUITABLE FOR *IN-VITRO*-DIAGNOSTIC USE!

ALL OTHER PRODUCTS NOT LABELED AS IVD ARE NOT SUITED FOR ANY CLINICAL USE (INCLUDING, BUT NOT LIMITED TO DIAGNOSTIC, THERAPEUTIC AND/OR PROGNOSTIC USE).

No claim or representations is intended for its use to identify any specific organism or for clinical use (included, but not limited to diagnostic, prognostic, therapeutic, or blood banking). It is rather in the responsibility of the user or - in any case of resale of the products - in the responsibility of the reseller to inspect and assure the use of the DNA/RNA/protein purification products of MACHEREY-NAGEL for a well-defined and specific application.

MACHEREY-NAGEL shall only be responsible for the product specifications and the performance range of MN products according to the specifications of in-house quality control, product documentation and marketing material.

This MACHEREY-NAGEL product is shipped with documentation stating specifications and other technical information. MACHEREY-NAGEL warrants to meet the stated specifications. MACHEREY-NAGEL's sole obligation and the customer's sole remedy is limited to replacement of products free of charge in the event products fail to perform as warranted. Supplementary reference is made to the general business terms and conditions of MACHEREY-NAGEL, which are printed on the price list. Please contact us if you wish to get an extra copy.

There is no warranty for and MACHEREY-NAGEL is not liable for damages or defects arising in shipping and handling (transport insurance for customers excluded), or out of accident or improper or abnormal use of this product; defects in products or

components not manufactured by MACHEREY-NAGEL, or damages resulting from such non-MACHEREY-NAGEL components or products.

MACHEREY-NAGEL makes no other warranty of any kind whatsoever, and SPECIFICALLY DISCLAIMS AND EXCLUDES ALL OTHER WARRANTIES OF ANY KIND OR NATURE WHATSOEVER, DIRECTLY OR INDIRECTLY, EXPRESS OR IMPLIED, INCLUDING, WITHOUT LIMITATION, AS TO THE SUITABILITY, REPRODUCTIVITY, DURABILITY, FITNESS FOR A PARTICULAR PURPOSE OR USE, MERCHANTABILITY, CONDITION, OR ANY OTHER MATTER WITH RESPECT TO MACHEREY-NAGEL PRODUCTS.

In no event shall MACHEREY-NAGEL be liable for claims for any other damages, whether direct, indirect, incidental, compensatory, foreseeable, consequential, or special (including but not limited to loss of use, revenue or profit), whether based upon warranty, contract, tort (including negligence) or strict liability arising in connection with the sale or the failure of MACHEREY-NAGEL products to perform in accordance with the stated specifications. This warranty is exclusive and MACHEREY-NAGEL makes no other warranty expressed or implied.

The warranty provided herein and the data, specifications and descriptions of this MACHEREY-NAGEL product appearing in MACHEREY-NAGEL published catalogues and product literature are MACHEREY-NAGEL's sole representations concerning the product and warranty. No other statements or representations, written or oral, by MACHEREY-NAGEL's employees, agent or representatives, except written statements signed by a duly authorized officer of MACHEREY-NAGEL are authorized; they should not be relied upon by the customer and are not a part of the contract of sale or of this warranty.

Product claims are subject to change. Therefore please contact our Technical Service Team for the most up-to-date information on MACHEREY-NAGEL products. You may also contact your local distributor for general scientific information. Applications mentioned in MACHEREY-NAGEL literature are provided for informational purposes only. MACHEREY-NAGEL does not warrant that all applications have been tested in MACHEREY-NAGEL laboratories using MACHEREY-NAGEL products. MACHEREY-NAGEL does not warrant the correctness of any of those applications.

Last updated: 07/2010, Rev. 03

Please contact:  
MACHEREY-NAGEL GmbH & Co. KG  
Tel.: +49 24 21 969-270  
tech-bio@mn-net.com

---

**Trademarks:**

NucleoSpin® is a registered trademark of MACHEREY-NAGEL GmbH & Co KG

RNA<sup>later</sup>® is a registered trademark of AMBION, Inc.

SPN® is a registered trademark of G-Biosciences

All used names and denotations can be brands, trademarks, or registered labels of their respective owner – also if they are not special denotation. To mention products and brands is only a kind of information (i.e., it does not offend against trademarks and brands and can not be seen as a kind of recommendation or assessment). Regarding these products or services we can not grant any guarantees regarding selection, efficiency, or operation.