



FMS2, Mini Multi Shaker, Adjustable – 3D & 2D Motion

- BAE MECHANISM (patented, 5 step angle adjustment system) allows orbital motion and various twist movement
- Sealed construction to prevent heat transfer from motor to platform for use in incubator and cold chamber
- Variable speed control, from 5 to 90rpm
- Orbit of 18mm for best mixing result especially for gel and blotting
- 20cm x 20cm platform for efficient incubator operation
- Dual mat to accommodate each vessels in the proper mat.


Model	FMS2
Speed range	5 ~ 90rpm
Motor	DC Warm Geared motor
Temp. resist	Lower than 65°C (Hum 90%)
Loadable Weight	1kg
Weight	3kg
Dimensions (mm)	W200xD200xH127
Wattage	50W
Power	Free Voltage (AC100-240V)



The table for angle adjustment:

Motion Select	FMS2
Orbital Motion	orbit of 18mm 5°,10°, 15°, 20° tilt angle
3D Motion	

View of actual movements:

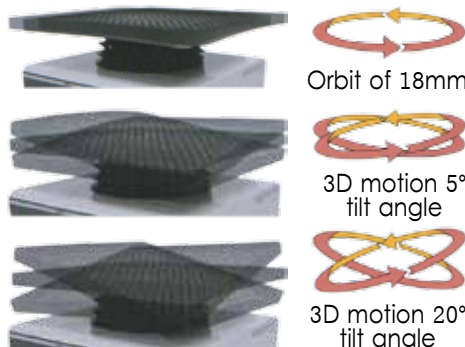

FMS3, Multi Shaker, Extremely quiet operation with various movement position, Adjustable – 3D & 2D Motion

- BAE MECHANISM (patented, 5 step angle adjustment system) allows orbital motion and various twist movement
- Variable speed control, from 30 to 300rpm
- Orbit of 18mm for best mixing result especially for gel and blotting
- High quality permanent DC brush less motor for gentle & powerful motion
- Continuous or timed operation
- Simple and easy of angle adjustment
- Dual mat to accommodate each vessels in the proper mat.



Model	FMS3
Speed range	30 ~ 300rpm
Motor	DC Warm Geared motor
Temp. resist	Lower than 55°C (Dry air 85%)
Loadable Weight	3kg
Timer	180min
Weight	12kg
Dimensions (mm)	W300xD360xH155
Wattage	30W
Power	Free Voltage (AC100-240V)

The table for angle adjustment:

Motion Select	FMS3
Orbital Motion	orbit of 18mm 5°,10°, 15°, 20° tilt angle
3D Motion	

Accessories:



Other Accessories:

Model	Description
CP50	50ml flask holder
CP100	100ml flask holder
CP250	250ml flask holder
CP500	500ml flask holder
FT302-50-13	Flask Rack for 50ml flask-Capacity: 13ea
FT302-100-13	Flask Rack for 100ml flask-Capacity: 13ea
FT302-250-8	Flask Rack for 250ml flask-Capacity: 8ea
FT302-500-5	Flask Rack for 500ml flask-Capacity: 5ea
FT302	Universal plate for flask holder