

DILUMAT S Options / Kit*




A.	TRACEABILITY	2
B.	DISPENSING / WEIGHING.....	3
C.	DISPENSING IN VIAL	4
C.1.	Height available with different nozzles	5
D.	Which tubing set to use ?	6
E.	Build the tubing set.....	6

*kits can be installed by the customer or an after sale service technician

**Options must be indicated on the offer to be installed by the manufacture

A. TRACEABILITY

Sample Identification

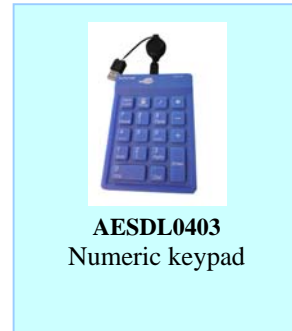

B. DISPENSING / WEIGHING



AESDL0300
Second pump kit

contains :
AESDL0394
Tubing set w/o priming nozzle (1u)



AESDL0420
Dilutube adapter kit

To use with sterile tubing DILUTUBE (AESDL0410)
or
dispensing nozzles DILUMAT 4

AESDI0011 : Dispensing nozzle
AESDI0011A : Dispensing angled nozzle 18°
AESDL0117 : Long dispensing nozzle



Max
2.5 KG




AESDL0398
Dispensing nozzle
Spécial Powder

To use with tubing AESDL0394

C. DISPENSING IN VIAL

LARGE VIAL



AESDI0043	AESDL0420
Bowl holding stand	Dilutube kit

With tubing set **AESDL0394**
and
AESDI0011A Dispensing angled nozzle 18°

LARGE VIAL



AESDI0043
Bowl holding stand

With tubing set **AESDL0394**

SMALL VIAL




AESDI0043	AESDL0420
Bowl holding stand	Dilutube kit

With **AESDL0113** Tubing set with long nozzle (15 cm)
(inner Ø pump tubing : 4.8 mm)


Low distribution volume
Sample weighing 1 to 10 g
(Final dilution: 100 ml)

C.1. Height available with different nozzles

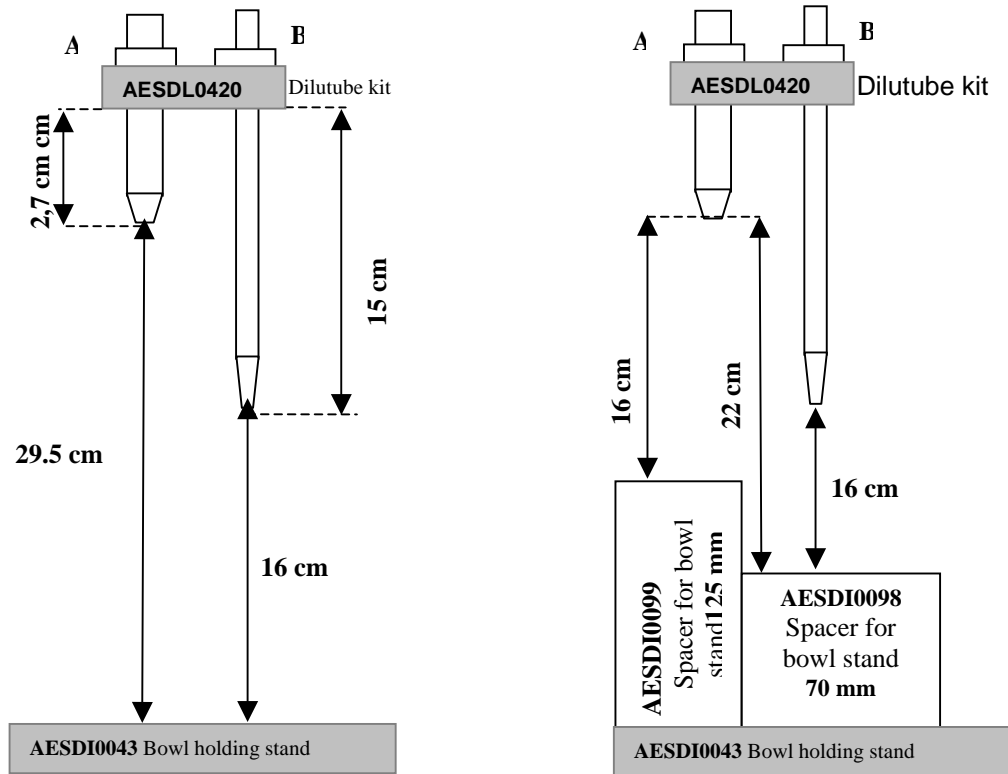


A = AESDI0011 Dispensing nozzle or AESDI0011A Dispensing angled nozzle 18°
With tubing AESDL0394

B = AESDL0117 Long dispensing nozzle (Tubing set AESDL0113)

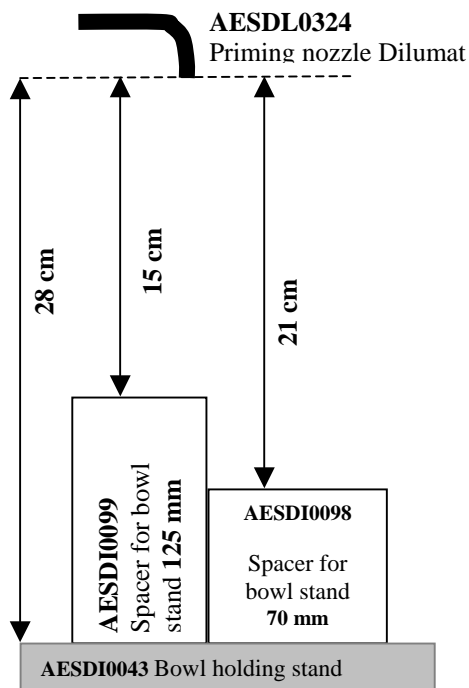

DILUMAT S Options / Kit*

D. Which tubing set to use ?

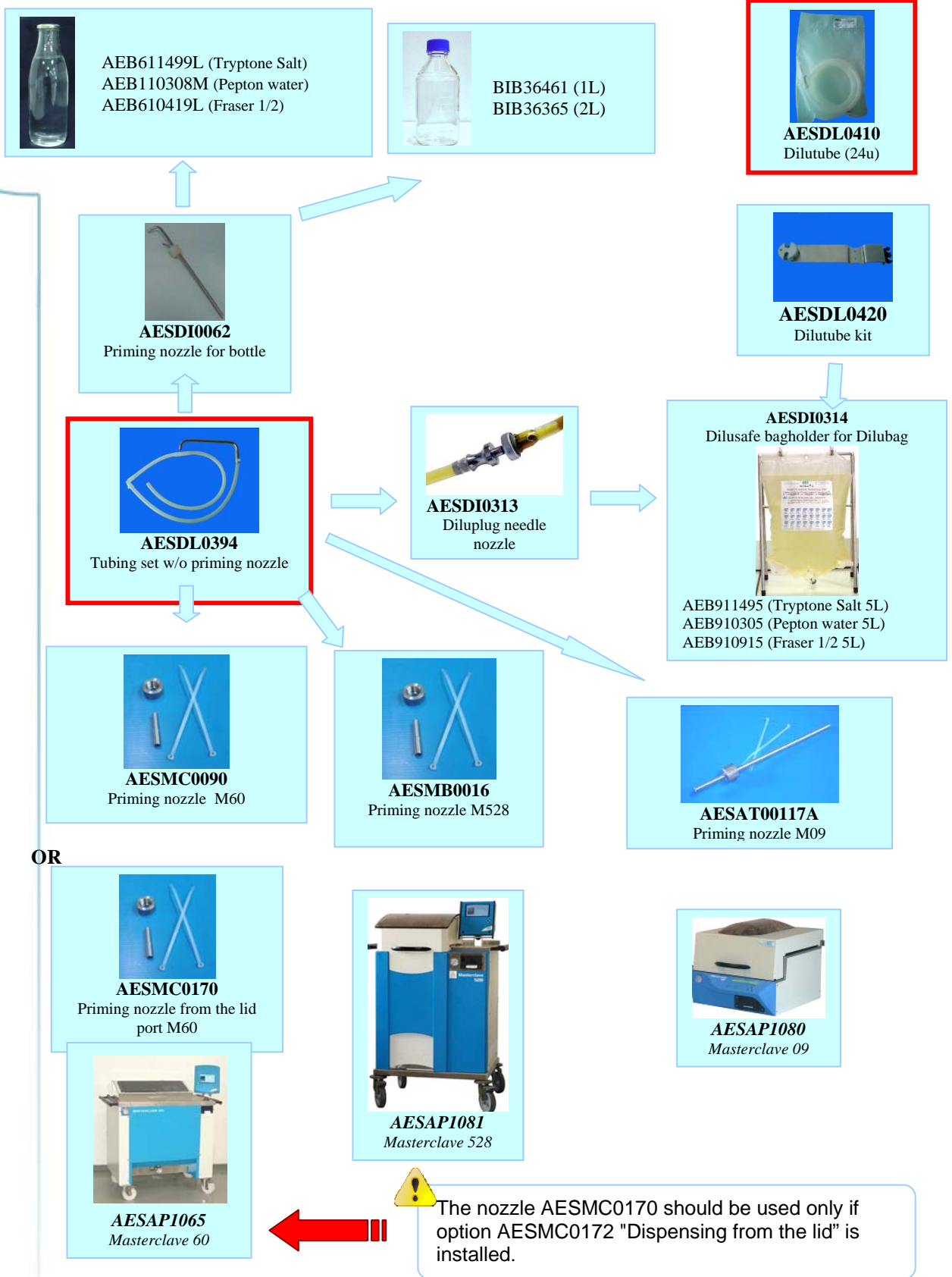


E. Build the tubing

DILUMAT S Options / Kit*

set

Complete tubing set Dilumat S for priming in bottle

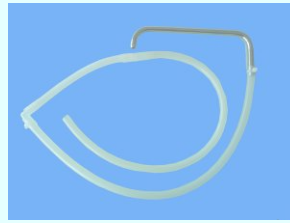

AESDL0394
Complete tubing set Dilumat S w/o
priming nozzle

+



AESDI0062
Priming nozzle for bottle
cap Ø280-340

AESAP0664
Tubing clamp

Accessories for priming nozzle AESDI0062 :



AESDI0308
SILICONE CAP W/
HOLE D36/44

AESDI0055
Silicone cap
Ø 280-340



AESDI0048
Filter 0,2µm (50u)

AESDL0394 Complete tubing set w/o priming nozzle



AESDL0371
Angled connector PE
(5u)

+



AESDL0038 *
Precut rotor tubing (5u)
or
AESDL0038A *
Rotor tubing (5m)

+



AESDL0372
Connector PE (5u)

+



AESDL0324
Dispensing nozzle
Dilumat S

Build with : - silicon tubing 6-10mm : **AESMN3928A** (5m) ou **AESMN3928** (25m)
- Clamp for tubing **AESAP0664**

* NB : The rotor tubing must be change each 100 hours

Tubing set Dilumat S SPECIAL POWDER



AESDL0371
Angled connector PE
(5u)

+



AESDL0038 *
Precut rotor tubing (5u)
or
AESDL0038A *
Rotor tubing (5m)

+



AESDL0372
connector PE (5u)

+



AESDL0398
Dispensing nozzle
Special powder

Build with : - silicon tubing 6-10mm : **AESMN3928A** (5m) ou **AESMN3928** (25m)
- Clamp for tubing **AESAP0664**

* NB : The rotor tubing must be change each 100 hours